


A11105 354519

Measurement Services:


REFERENCE

NIST
PUBLICATIONS

Radiance Temperature Calibrations

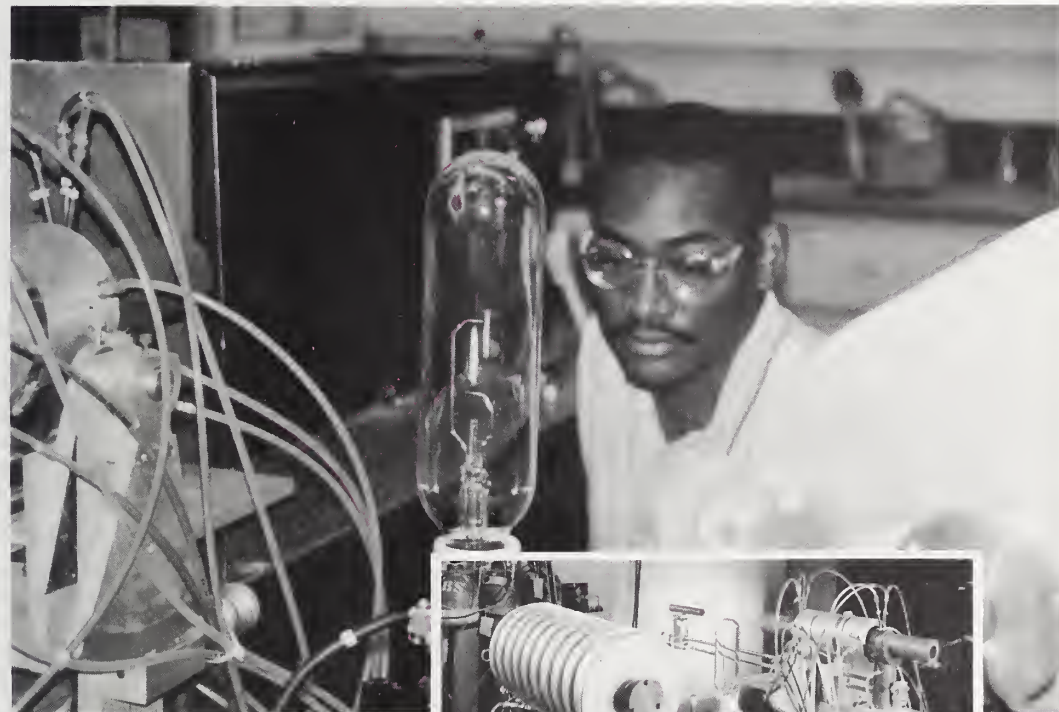

NIST
Special
Publication
250-43

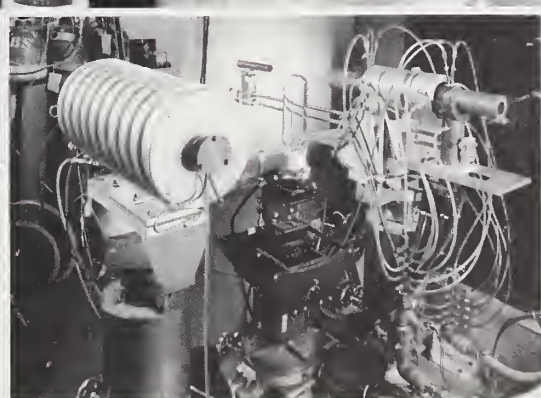

Charles E. Gibson, Benjamin K. Tsai, and Albert C. Parr

U.S. Department of Commerce
Technology Administration
National Institute of Standards and Technology

QC
100
U57
10.250-43
1998

The National Institute of Standards and Technology was established in 1988 by Congress to “assist industry in the development of technology . . . needed to improve product quality, to modernize manufacturing processes, to ensure product reliability . . . and to facilitate rapid commercialization . . . of products based on new scientific discoveries.”

NIST, originally founded as the National Bureau of Standards in 1901, works to strengthen U.S. industry’s competitiveness; advance science and engineering; and improve public health, safety, and the environment. One of the agency’s basic functions is to develop, maintain, and retain custody of the national standards of measurement, and provide the means and methods for comparing standards used in science, engineering, manufacturing, commerce, industry, and education with the standards adopted or recognized by the Federal Government.

As an agency of the U.S. Commerce Department’s Technology Administration, NIST conducts basic and applied research in the physical sciences and engineering, and develops measurement techniques, test methods, standards, and related services. The Institute does generic and precompetitive work on new and advanced technologies. NIST’s research facilities are located at Gaithersburg, MD 20899, and at Boulder, CO 80303. Major technical operating units and their principal activities are listed below. For more information contact the Publications and Program Inquiries Desk, 301-975-3058.

Office of the Director

- National Quality Program
- International and Academic Affairs

Technology Services

- Standards Services
- Technology Partnerships
- Measurement Services
- Technology Innovation
- Information Services

Advanced Technology Program

- Economic Assessment
- Information Technology and Applications
- Chemical and Biomedical Technology
- Materials and Manufacturing Technology
- Electronics and Photonics Technology

Manufacturing Extension Partnership Program

- Regional Programs
- National Programs
- Program Development

Electronics and Electrical Engineering Laboratory

- Microelectronics
- Law Enforcement Standards
- Electricity
- Semiconductor Electronics
- Electromagnetic Fields¹
- Electromagnetic Technology¹
- Optoelectronics¹

Chemical Science and Technology Laboratory

- Biotechnology
- Physical and Chemical Properties²
- Analytical Chemistry
- Process Measurements
- Surface and Microanalysis Science

Physics Laboratory

- Electron and Optical Physics
- Atomic Physics
- Optical Technology
- Ionizing Radiation
- Time and Frequency¹
- Quantum Physics¹

Materials Science and Engineering Laboratory

- Intelligent Processing of Materials
- Ceramics
- Materials Reliability¹
- Polymers
- Metallurgy
- NIST Center for Neutron Research

Manufacturing Engineering Laboratory

- Precision Engineering
- Automated Production Technology
- Intelligent Systems
- Fabrication Technology
- Manufacturing Systems Integration

Building and Fire Research Laboratory

- Structures
- Building Materials
- Building Environment
- Fire Safety Engineering
- Fire Science

Information Technology Laboratory

- Mathematical and Computational Sciences²
- Advanced Network Technologies
- Computer Security
- Information Access and User Interfaces
- High Performance Systems and Services
- Distributed Computing and Information Services
- Software Diagnostics and Conformance Testing

¹At Boulder, CO 80303.

²Some elements at Boulder, CO.

NIST Special Publication 250-43

NIST MEASUREMENT SERVICES: Radiance Temperature Calibrations

Charles E. Gibson, Benjamin K. Tsai, and Albert C. Parr

Optical Technology Division
Physics Laboratory
National Institute of Standards and Technology
Gaithersburg, MD 20899-0001

January 1998



U.S. DEPARTMENT OF COMMERCE

William M. Daley, Secretary

Technology Administration

Gary R. Bachula, Acting Under Secretary for Technology

National Institute of Standards and Technology

Raymond G. Kammer, Director

National Institute of Standards and Technology Special Publication 250-43
Natl. Inst. Stand. Technol. Spec. Publ. 250-43, 123 pages (Jan. 1998)
CODEN: NSPUE2

U.S. GOVERNMENT PRINTING OFFICE
WASHINGTON: 1998

For sale by the Superintendent of Documents, U.S. Government Printing Office, Washington, DC 20402-9325

PREFACE

The calibration and related measurement services of the National Institute of Standards and Technology are intended to assist the makers and users of precision measuring instruments in achieving the highest possible levels of accuracy, quality, and productivity. NIST offers over 300 different calibrations, special tests, and measurement assurance services. These services allow customers to directly link their measurement systems to measurement systems and standards maintained by NIST. These services are offered to the public and private organizations alike. They are described in NIST Special Publication (SP) 250, NIST Calibration Services Users Guide.

The Users Guide is supplemented by a number of Special Publications (designated as the “SP250 Series”) that provide detailed descriptions of the important features of specific NIST calibration services. These documents provide a description of the: (1) specifications for the services; (2) design philosophy and theory; (3) NIST measurement system; (4) NIST operational procedures; (5) assessment of the measurement uncertainty including random and systematic errors and an error budget; and (6) internal quality control procedures used by NIST. These documents will present more detail than can be given in NIST calibration reports, or than is generally allowed in articles in scientific journals. In the past, NIST has published such information in a variety of ways. This series will make this type of information more readily available to the user.

This document, SP250-43 (1997), NIST Measurement Services: Radiance Temperature Calibrations, is a revision of SP250-7 (1987). It covers the calibration of radiance temperature standards for disappearing filament optical pyrometers, ribbon filament lamps, and radiation thermometers (test numbers 35010C-35070S in SP250, NIST Calibration Services Users Guide). Inquiries concerning the technical content of this document or the specifications for these services should be directed to the author or to one of the technical contacts cited in SP250.

NIST welcomes suggestions on how publications such as this might be made more useful. Suggestions are also welcome concerning the need for new calibrations services, special tests, and measurement assurance programs.

Stanley D. Rasberry
Director
Measurement Services

Katharine B. Gebbie
Director
Physics Laboratory

ABSTRACT

This document describes the realization and dissemination of the International Temperature Scale of 1990 (ITS-90) above 700 °C at the National Institute of Standards and Technology (NIST). By using fundamental principles of blackbody physics, the ITS-90 scale is first fixed at the freezing point of gold ($T_{90} = 1064.18$ °C) and is then extended to temperatures between 700 °C and 2700 °C by determining the ratio of the spectral radiance of a tungsten ribbon filament lamp to that of a gold fixed-point blackbody at a wavelength of 655.3 nm. A description of the facilities in the NIST Radiance Temperature Calibration Laboratory is given, along with a discussion of the wavelength calibration, size of source, and linearity issues. The use and calibration of radiance temperature standards are explained. Values of stability and uncertainties in the scale are reported. The expanded uncertainties ($k = 2$) at 800 °C and 2300 °C are 0.6 °C and 1.4 °C for ribbon filament lamps, 4 °C and 7 °C for disappearing filament optical pyrometers, and 2 °C and 3 °C for infrared radiation thermometers.

KEYWORDS: gold-point blackbody, International Temperature Scale (ITS-90), pyrometers, pyrometry, radiance temperature calibration, radiation thermometry, tungsten ribbon filament lamps

TABLE OF CONTENTS

Abstract	iv
1. Introduction	1
2. Measurement theory	4
3. Realization of the U.S. national scale of radiance temperature	9
3.1 1990 NIST scale of radiance temperature	9
3.2 Temperature standards	9
3.2.1 Gold-point blackbody	9
3.2.2 Tungsten ribbon filament lamps	11
3.2.3 Working standard lamps	14
3.2.4 Variable-temperature blackbody	15
3.3 NIST photoelectric pyrometer	19
3.3.1 Measurement system	19
3.3.2 Optical system	21
3.3.3 Interference filters	22
3.3.4 Detectors	23
3.3.5 Linearity of response	24
3.3.6 Wavelength calibration	27
3.3.7 Size of source effect	29
3.4 Radiance temperature scale uncertainties	30
3.5 Quality control	32
4. Tungsten ribbon filament lamp calibrations	37
4.1 Lamp preparation	37
4.2 Lamp calibration	37
4.3 Lamp data analysis	38
4.4 Lamp calibration uncertainty	40
5. Disappearing filament optical pyrometer calibrations	44
5.1 Optical pyrometer preparation	44
5.2 Optical pyrometer calibration	44
5.3 Optical pyrometer data analysis	44
5.4 Optical pyrometer calibration uncertainty	45

TABLE OF CONTENTS (continued)

6. Radiation thermometer calibrations	49
6.1 Radiation thermometer preparation	49
6.2 Radiation thermometer calibration	49
6.3 Radiation thermometer data analysis	49
6.4 Radiation thermometer calibration uncertainty	49
7. Modifications to the NIST photoelectric pyrometer	53
Acknowledgments	54
References	55
Appendix A: Ribbon filament lamp sample calibration report	A1
Appendix B: Optical pyrometer sample calibration report	B1
Appendix C: Radiation thermometer sample calibration report	C1
Appendix D: Radiance temperature measurement services laboratory procedure manual	D1
Appendix E: How to request radiance temperature calibration services	E1

LIST OF FIGURES

Figure 1	Schematic of optics for the NIST photoelectric pyrometer.....	4
Figure 2	Blackbody spectral distribution. As the temperature increases, the peak moves toward shorter wavelengths, and the slope increases at each wavelength.....	5
Figure 3	Schematic of gold-point blackbody.....	10
Figure 4	Gold-point blackbody liquid-to-solid phase transition.....	11
Figure 5	Gold-point blackbody control system	12
Figure 6	Schematic of lamp current monitoring system.....	13
Figure 7	Typical lamp currents for the General Electric 30A/T24/6 ribbon filament lamp	14
Figure 8	Ribbon filament lamp.....	16
Figure 9	Calibration history of the working standard (no. SL20)	17
Figure 10	Schematic of the variable temperature blackbody.....	17
Figure 11	Variable temperature blackbody control system	18
Figure 12	Spatial scan of the variable temperature blackbody at 2000 °C. The temperature differences are from the center position (0,0).	18
Figure 13	Variable temperature blackbody PID control. After decreasing the set point by 100 °C, the blackbody stabilizes in less than 10 min.....	19
Figure 14	NIST photoelectric pyrometer.....	20
Figure 15	NIST radiance temperature laboratory measurement system.....	21
Figure 16	NIST photoelectric pyrometer filter transmittances. The wavelength of the pyrometer is selected by using a 1 nm and a 10 nm spectral bandwidth (SBW) filter	22
Figure 17	The linearity correction factors for the PEP	25
Figure 18	Relative spectral response of the NIST photoelectric pyrometer	28
Figure 19	The size of source corrections for the variable temperature blackbody to the working standard lamp comparison and the working standard lamp to the gold-point blackbody.....	30
Figure 20	Calibration history of the ribbon filament lamp check standards.....	35
Figure 21	Calibration history of the radiation thermometer check standard (Minolta Cyclops 52).....	36
Figure 22	Calculation of I_{corr} from T_{mean} and $\Delta I/\Delta T$	39
Figure 23	Graph of T_{mean} and $\Delta I/\Delta T$	40
Figure 24	Modified NIST photoelectric pyrometer	53

LIST OF TABLES

Table 1	Radiance temperature calibration services	2
Table 2	Calibration uncertainties	3
Table 3	Current amplifier normalization factors	24
Table 4	Detector system correction factors and linearity coefficients	25
Table 5	Default hardware settings used by calibration program	26
Table 6	Corrections (°C) to the radiance temperature at 900 nm and at 1000 nm	29
Table 7	Typical values of WS variables and parameters.....	33
Table 8	Uncertainty budget for the NIST WS lamp spectral radiance realization	34
Table 9	Typical values of TL variables and parameters.....	42
Table 10	Uncertainty budget for the TL spectral radiance calibration	43
Table 11	Typical values of DF variables and parameters.....	47
Table 12	Uncertainty budget for the DF spectral radiance calibration.....	48
Table 13	Typical values of RT variables and parameters	51
Table 14	Uncertainty budget for the RT spectral radiance calibration.....	52

1. Introduction

The National Institute of Standards and Technology (NIST) is responsible for the establishment, maintenance, and dissemination of the U.S. National Scale of Radiance Temperature through its Radiance Temperature Calibration Laboratory (RTCL), which provides temperature measurements above 700 °C. The 1990 NIST Radiance Temperature Scale [1] is based upon the NIST realization of the International Temperature Scale 1990 [2] using principles from blackbody physics. Beginning with the determination of the freezing point of gold ($T_{90} = 1064.18$ °C), the scale is realized for the range from 700 °C to 4200 °C using radiometric methods. The radiance temperature (also referred to as the brightness or apparent temperature) of a radiator is equivalent to the temperature of a blackbody with the same radiant intensity as that of the radiator's surface at a chosen wavelength.

This document describes the determination of the ratio of the spectral radiance of a tungsten ribbon filament lamp used as a temperature standard to that of a gold fixed-point blackbody at a wavelength of 655.3 nm. A description of the NIST RTCL along with the use and calibration of radiance temperature standards therein is presented in detail. Issues relating to the wavelength calibration, size of source, and linearity of the NIST Photoelectric Pyrometer (PEP) are then discussed. Stability and uncertainties in the scale are considered in detail within the framework of the NIST policy on uncertainty statements described in reference [3].

The NIST disseminates the radiance temperature scale by two methods. The first method involves calibrating ribbon filament lamps as radiance temperature standards in the temperature range from 800 °C to 2300 °C with expanded uncertainties ($k = 2$)¹ of 0.6 °C at 800 °C, and 1.4 °C at 2300 °C. The second method involves the calibration of customer pyrometers and radiation thermometers. Disappearing filament optical pyrometers are calibrated from 800 °C to 4200 °C with expanded uncertainties ($k = 2$) of 4 °C at 800 °C, 7 °C at 2300 °C, and 25 °C at 4200 °C. Radiation thermometers are measured from 800 °C to 2700 °C with expanded uncertainties ($k = 2$) of 2 °C at 800 °C, and 3 °C at 2700 °C.

The radiance temperature calibration services shown in Table 1 are offered through the NIST Calibration Services Program and are listed in the Radiation Thermometry section of the Thermodynamics Quantities Chapter of the NIST Calibration Services Users Guide [4]. A summary of the calibration uncertainties is given in Table 2. This laboratory's quality system is based on the ANSI/NCSL Z540-1-1994 standard and the ISO/IEC Guide 25.

¹ Throughout this paper, all uncertainty values are given as an expanded uncertainty with coverage factor $k = 2$, thus a two standard deviation estimate. Uncertainties of fundamental units given as a combined standard uncertainty in other documents are restated as an expanded uncertainty ($k = 2$).

Table 1. Radiance temperature calibration services

Service ID No.	Measurement Description
Calibration reports are issued giving the radiance temperature of the blackbody at 655.3 nm versus the scale reading, output current, or output voltage.	
35010C	Radiance Temperature Standard, Disappearing Filament Optical Pyrometer (800 °C to 2400 °C, 4 to 12 points, 1 range)
35020C	Radiance Temperature Standard, Disappearing Filament Optical Pyrometer (each additional range up to 4200 °C, only available with No. 35010C)
35030C	Additional Interpolated Values
35040C	Radiance Temperature Standard, Disappearing Filament Optical Pyrometer (800 °C to 4200 °C, 3 or fewer points, 1 range)
Calibration reports are issued giving the radiance temperature of the lamp at 655.3 nm versus the lamp current.	
35050C	Radiance Temperature Standard, Ribbon Filament Lamp (800 °C to 2300 °C, 6 to 16 points)
35060C	Radiance Temperature Standard, Ribbon Filament Lamp (800 °C to 2300 °C, 5 or fewer points)
Test reports are issued giving the radiance temperature of the blackbody at 655.3 nm, 900 nm, or 1000 nm versus the indicated reading, output current, or output voltage.	
35070S	Special Tests of Radiation Thermometers

Table 2. Calibration uncertainties

Standard	Temperature Range	Expanded Uncertainty
Ribbon filament lamp	800 °C to 1600 °C	0.7 °C
	1600 °C to 1900 °C	1.0 °C
	1900 °C to 2300 °C	1.5 °C
Leeds & Northrup Model 8000 series	800 °C to 1600 °C	4 °C
	1600 °C to 2100 °C	5 °C
	1900 °C to 2300 °C	7 °C
Disappearing filament optical pyrometer	2400 °C to 2700 °C	8 °C
	2700 °C to 3200 °C	17 °C
	3200 °C to 4200 °C	25 °C
Pyrometer Instrument Model 95	800 °C to 1400 °C	3 °C
	1400 °C to 1800 °C	4 °C
	1800 °C to 2400 °C	5 °C
Disappearing filament optical pyrometer	2400 °C to 2700 °C	8 °C
	2700 °C to 3200 °C	12 °C
Infrared radiation thermometer	800 °C to 2000 °C	2 °C
	2000 °C to 2700 °C	3 °C

The expanded uncertainty (coverage factor $k = 2$) is equivalent to a two standard deviation estimate.

2. Measurement theory

Temperatures above the freezing point of silver (1234.93 K or 961.78 °C) are defined on the International Temperature Scale of 1990 (ITS-90) [2] in terms of the ratio of spectral radiances of two blackbody sources, one of which is maintained at the temperature of freezing silver, gold (1337.33 K or 1064.18 °C), or copper (1357.77 K or 1084.62 °C). The 1990 NIST Scale of Radiance Temperature (1990 NIST) is a realization of the ITS-90 using a gold fixed-point blackbody. In this section, the blackbody temperature will be defined in terms of the spectral radiance. Using the signal measurement equation, the measurement equation for the calibration of a transfer standard will be derived.

The signal measurement equation, defined by Nicodemus and Kostkowski in 1978 [5], relates the detector signal output, S [V]², to the source flux input parameters through a detector responsivity term, R_λ [A·W⁻¹], by the integral relationship,

$$S = \int_A \int_\omega \int_\lambda R_\lambda \cdot L_\lambda \cdot d\lambda \cdot \cos \theta \cdot d\omega \cdot dA, \quad (1)$$

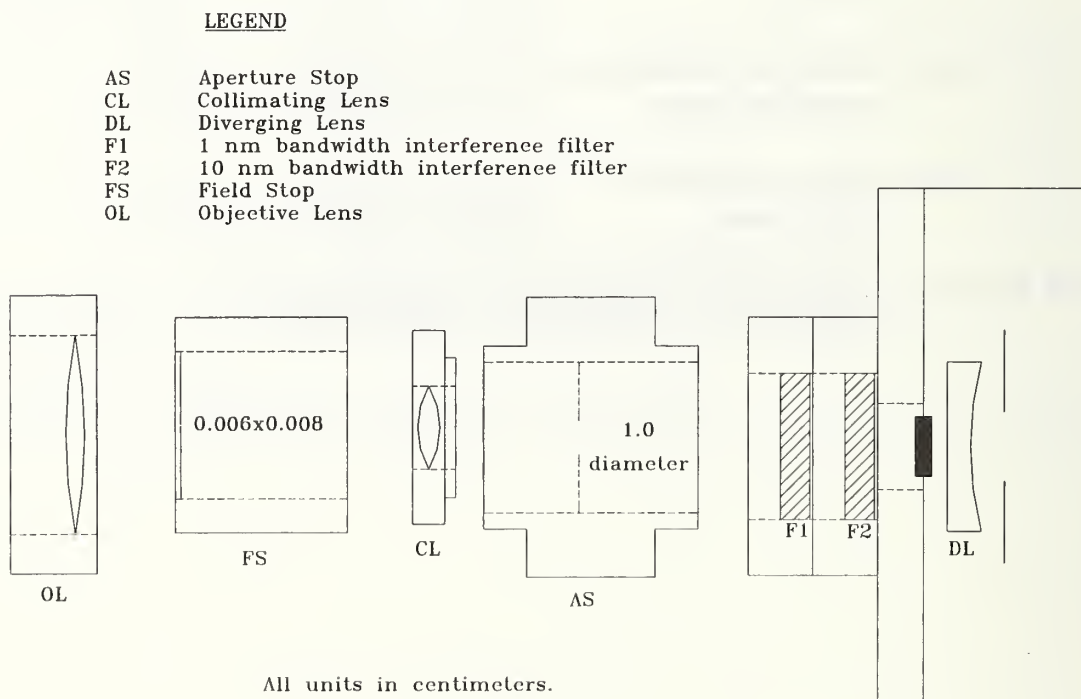


Figure 1. Schematic of optics for the NIST photoelectric pyrometer.

² As an aid to the reader, the appropriate coherent SI unit in which a quantity should be expressed is indicated in brackets when the quantity is first introduced.

where L_λ [$\text{W} \cdot \text{m}^{-3} \cdot \text{sr}^{-1}$] is the spectral radiance, $d\lambda$ [m] is the wavelength band, θ [rad] is the angle between the aperture normal and the line connecting the aperture centers, $d\omega$ [sr] is the differential solid angle originating at the source aperture as defined by the detector aperture, and dA [m^2] is the differential source aperture area. For the NIST PEP in figure 1, the spectral responsivity R_λ includes the spectral transmittance of the interference filters, the spectral transmittance of all other optical elements, and the spectral responsivity of the detector. In terms of its specific components, the spectral responsivity is

$$R_\lambda = \tau_{\lambda,\text{OL}} \cdot \tau_{\lambda,\text{CL}} \cdot \tau_{\lambda,\text{F1}} \cdot \tau_{\lambda,\text{F2}} \cdot \tau_{\lambda,\text{DL}} \cdot \tau_{\lambda,\text{EC}} \cdot R_{\lambda,\text{PEP}}, \quad (2)$$

where $\tau_{\lambda,\text{OL}}$ is the spectral transmittance of the objective lens, $\tau_{\lambda,\text{CL}}$ is the spectral transmittance of the collimating lens, $\tau_{\lambda,\text{F1}}$ is the spectral transmittance of the 1 nm bandpass interference filter, $\tau_{\lambda,\text{F2}}$ is the spectral transmittance of the 10 nm bandpass interference filter, $\tau_{\lambda,\text{DL}}$ is the spectral transmittance of the diverging lens, $\tau_{\lambda,\text{EC}}$ is the spectral transmittance of the evacuated window cell, and $R_{\lambda,\text{PEP}}$ is the detector absolute spectral responsivity.

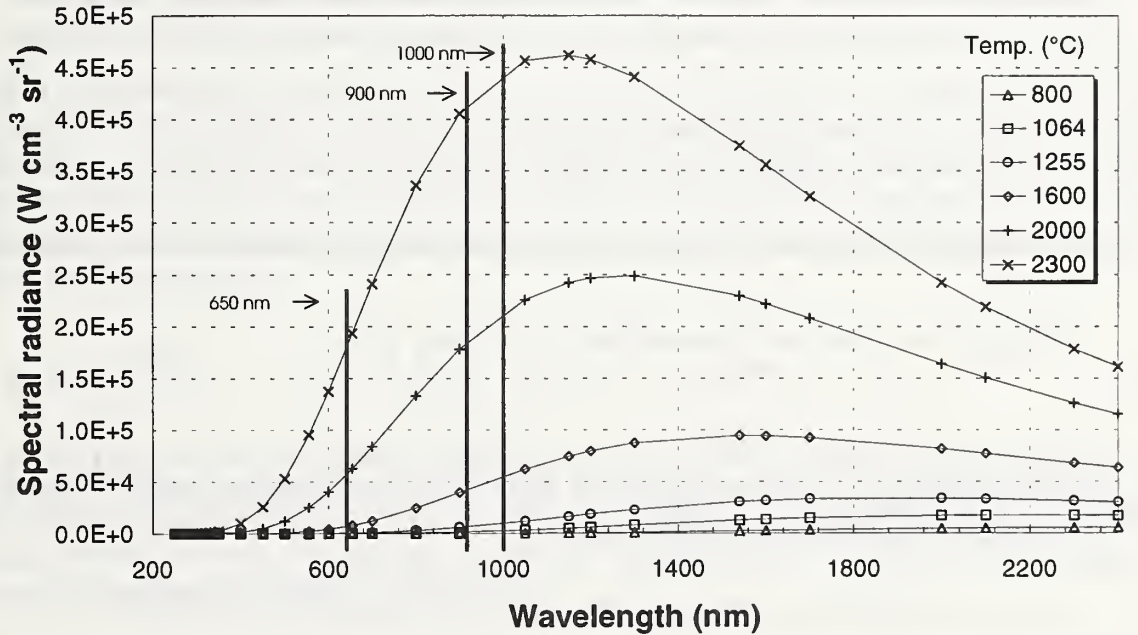


Figure 2. Blackbody spectral distribution. As the temperature increases, the peak moves towards shorter wavelengths, and the slope increases at each wavelength.

Spectral radiance L_λ is the radiant power contained in a defined area, solid angle, direction, and wavelength interval,

$$L_\lambda = d^3\Phi / dA \cdot \cos \alpha \cdot d\beta \cdot d\lambda , \quad (3)$$

where Φ is the radiant flux [W], A is the source area, α is the plane angle between the surface normal and the direction of propagation [rad], β is the solid angle about that direction [sr], and λ is the wavelength [m]. A few blackbody distributions from 250 nm to 2500 nm between 800 °C and 2300 °C are illustrated in figure 2. For an ideal blackbody in a vacuum environment, the relation between spectral radiance, wavelength, and temperature is given by the Planck equation [6],

$$L_\lambda = \frac{c_{1L}}{\lambda^5 \cdot (\exp(c_2 / (\lambda \cdot T)) - 1)} , \quad (4)$$

where c_{1L} is the first radiation constant in radiance form with a value of 119.1044 W·nm², c_2 is the second radiation constant with a value of 1.4388×10^{-7} nm·K [7], λ is the wavelength in vacuum, and T is the temperature [K]. The Planck equation in the form of eq (4) is the definition for an ideal blackbody radiator. When using a non-ideal blackbody such as a fixed-point blackbody to realize the temperature scale, the following form of the Planck equation is used,

$$L_\lambda = \frac{\epsilon_\lambda \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot (\exp(c_2 / (n_\lambda \cdot \lambda \cdot T)) - 1)} , \quad (5)$$

and temperature is defined as a function of spectral radiance using the following equation

$$T(\lambda, L_\lambda) = \frac{c_2}{n_\lambda \cdot \lambda \cdot \ln(1 + \epsilon_\lambda \cdot c_{1L} / n_\lambda^2 \cdot \lambda^5 \cdot L_\lambda)} , \quad (6)$$

where ϵ_λ is the spectral emissivity of the blackbody (which is equal to unity in the case of an ideal blackbody), λ is the wavelength in air, and n_λ is the refractive index of air at 15 °C and 101,325 Pa. From the Cauchy [8] formula,

$$n_\lambda = 1 + (2726.43 + 12.288 \text{ nm}^2 / (\lambda^2 \times 10^{-6}) + 0.3555 \text{ nm}^4 / (\lambda^4 \times 10^{-12})) \times 10^{-7} , \quad (7)$$

the refractive index of air at 655.3 nm is about 1.00028. The temperatures described in this document are radiance temperatures. The radiance temperature of a radiator is equivalent to the temperature of a blackbody with the same radiant intensity of the radiator's surface at a specified wavelength. The relationship between the radiance temperature and the true or thermodynamic temperature of a blackbody is given by:

$$\frac{1}{T} = \frac{1}{T_\lambda} + \frac{\lambda}{c_2} \cdot \ln \varepsilon_\lambda, \quad (8)$$

where T is the thermodynamic temperature, T_λ is the radiance temperature, λ is the mean effective wavelength of the NIST PEP (655.3 nm), and ε_λ is the estimated emissivity of the blackbody (0.99).

Determination of the spectral radiance temperature of a working standard (WS) lamp requires measurement of the ratio r_1 ,

$$r_1 = \frac{S_{\text{WS}}}{S_{\text{Au}}}, \quad (9)$$

of the signals from the transfer standard and the goldpoint blackbody (Au) with the NIST PEP. From eq (2), this measured ratio is actually

$$r_1 = \frac{\int_{\lambda} \int_{A_s} \int_{\omega} R_\lambda \cdot (L_\lambda \cdot \cos \theta_s \cdot d\omega \cdot dA_s)_{\text{WS}} \cdot d\lambda}{\int_{\lambda} \int_{A_s} \int_{\omega} R_\lambda \cdot (L_\lambda \cdot \cos \theta_s \cdot d\omega \cdot dA_s)_{\text{Au}} \cdot d\lambda}. \quad (10)$$

To simplify the complex expression in eq (10), it is assumed that the spectral radiances $L_{\lambda,\text{WS}}$ and $L_{\lambda,\text{Au}}$, spectral responsivities $R_{\lambda,\text{WS}}$ and $R_{\lambda,\text{Au}}$, and amplifier gains G_{WS} and G_{Au} are independent of both direction and spatial location. Furthermore, these three variables can be defined by unique for the transfer source and the gold-point blackbody at equivalent wavelengths (to be defined later) over the same small wavelength band $d\lambda$. The solid angle terms can be replaced with the definition of the solid angle,

$$\omega = \frac{A_d \cdot \cos \theta_d}{D^2}, \quad (11)$$

where A_d is the detector area [m^2], θ_d is the angle between the optical axis and the normal to the detector surface [rad], and D is the distance between the detector area and the source area [m]. Assuming that the areas are independent of direction and that the solid angles are independent of area or spatial location, this ratio then becomes

$$r_1 = \frac{(R_\lambda \cdot L_\lambda \cdot A_s \cdot \cos \theta_s \cdot A_d \cdot \cos \theta_d)_{\text{WS}}}{(R_\lambda \cdot L_\lambda \cdot A_s \cdot \cos \theta_s \cdot A_d \cdot \cos \theta_d)_{\text{Au}}} \cdot \frac{D_{\text{Au}}^2}{D_{\text{WS}}^2}, \quad (12)$$

where A_s is the source area [m^2], θ_s is the angle between the optical axis and the normal to the source surface [rad], $R_{\lambda,\text{WS}}$ is the detector responsivity when viewing the transfer source, and $R_{\lambda,\text{Au}}$ is the detector responsivity when viewing the gold-point blackbody. The spectral responsivities $R_{\lambda,\text{WS}}$ and $R_{\lambda,\text{Au}}$ are assumed to be equal. If the source aperture areas dA_{TS} and

dA_{Au} , source solid angles $d\omega_{WS}$ and $d\omega_{Au}$, and the source inclination angles θ_{WS} and θ_{Au} are the same for measuring both the transfer standard and the gold-point blackbody, the measured ratio simplifies to the expression,

$$r_1 = \frac{L_{\lambda, WS}}{L_{\lambda, Au}} = \frac{L_{\lambda}(T_{WS})}{L_{\lambda}(T_{Au})}. \quad (13)$$

This relationship is the defining equation for the ITS-90 above 1337.33 K. In terms of eq (4), it can be written at a discrete wavelength, λ , as

$$r_1 = \frac{L_{\lambda}(T_{WS})}{L_{\lambda}(T_{Au})} = \frac{\exp(c_2 / (n_{\lambda} \cdot \lambda \cdot T_{Au})) - 1}{\exp(c_2 / (n_{\lambda} \cdot \lambda \cdot T_{WS})) - 1}, \quad (14)$$

where $L_{\lambda}(T_{WS})$ and $L_{\lambda}(T_{Au})$ are the spectral radiances of the two blackbodies at temperatures T and T_{Au} , T_{Au} is the temperature of freezing gold defined as 1337.33 K, and r_1 is their ratio. In principle, a measurement of the ratio at a discrete wavelength with a linear response instrument yields the value of T .

The radiance temperature scale is typically maintained and disseminated on tungsten ribbon filament lamps, which possess a repeatable lamp current versus radiance temperature relationship. At the NIST, a pyrometer system is presently being used with a mean effective wavelength of 655.3 nm. This method requires that the pyrometer relative spectral response extends only over a small spectral range, or is known accurately enough to determine the wavelength at which the integrands of eq (10) have the same ratio as the integrals. Equation (13) above is an approximation, and, in practice, corrections, which will be presented in the next section, are used.

3. Realization of the U.S. national scale of radiance temperature

3.1 1990 NIST scale of radiance temperature

The reference temperature standard, a gold fixed-point blackbody (Au) with a temperature (T_{Au}) of 1064.18 °C (1337.33 K), and the Planck radiation law are used to realize and disseminate the 1990 NIST Radiance Temperature Scale. Equation (5) is used to calculate the spectral radiance $L_{\lambda,Au}(\lambda, T_{Au})$ of the fixed-point blackbody for $\lambda = 655.3$ nm in all the measurements of this calibration facility. Measurements are performed from 800 °C to 2300 °C for lamps, from 800 °C to 2700 °C for radiation thermometers, and extrapolated to 4200 °C for some disappearing filament optical pyrometers.

The NIST PEP is the transfer device used to compare the spectral radiances of the sources by the direct substitution method. The signals are corrected for size of source, amplifier gain, and linearity. With these corrections, the spectral radiance and the radiance temperature can be determined from eqs (12) and (13). The NIST PEP is a filtered radiometer that uses two interference filters to select the bandpass. The spectral bandwidth is 5 nm with a mean effective wavelength of 655.3 nm. A photomultiplier tube with an S-20 spectral response is used in the dc mode. The measurement spot size is a 0.6 mm by 0.8 mm rectangle.

A high stability vacuum lamp operated at a single radiance temperature of approximately 1255 °C is the working standard (WS). By equating eqs (9) and (13), the spectral radiance ratio is

$$r_1 = \frac{L_{\lambda}(T_{WSI})}{L_{\lambda}(T_{Au})} = \frac{S_{WSI}}{S_{Au}}. \quad (15)$$

After applying correction factors to the signals in eq (15) for amplifier calibration (C_A), linearity (C_L), and size of source (C_S), the spectral radiance of the WS1 lamp can be written as

$$L_{WS} = \frac{\epsilon_{Au} \cdot c_{1L}}{n_{\lambda}^2 \cdot \lambda^5 \cdot \left[\exp\left(\frac{c_2}{n_{\lambda} \cdot \lambda \cdot T_{Au}}\right) - 1 \right]} \cdot r_1 \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{Au}}. \quad (16)$$

References regarding linearity issues are presented in Section 3.3.5. The uncertainty analysis in Section 3.6 can then be derived from eq (15) for the NIST radiance temperature scale up to the calibration of the spectral radiance of the working standard lamp.

3.2 Temperature standards

3.2.1 Gold-point blackbody

In the RTCL, a gold fixed-point blackbody with a calculated emissivity of 0.9999, designed and built by the NIST Optical Technology Division, is the primary standard used to realize the 1990 NIST Radiance Temperature Scale. The blackbody in figure 3 consists of a graphite cavity, a crucible of gold, and a cylindrical heat-pipe furnace. The cavity, which is

76 mm in length and 6 mm in diameter and has a 60° conical end shown, is made from Ultra “F” grade graphite (spectrographic purities of 10 ppm or less). Surrounding this cavity is a crucible containing 0.99999% pure gold. The cavity, along with graphite rings and silica glass spacers, is placed in an alumina tube. The front rings define a solid angle with a $f/6$ field of view, while the back rings support the thermocouple. A furnace (see fig. 3), which consists of a sodium heat-pipe heated by two semi-cylindrical ceramic heater elements inside of a mullite tube, is enclosed in a water-cooled housing (631 mL/min) and is operated in an argon environment (37 mL/min with furnace door closed and 235 mL/min with furnace door open).

The gold-point blackbody liquid-to-solid phase transition is shown in figure 4. The duration of a melt or freeze plateau is approximately 40 min, and the time delay between these observation periods is about 45 min. Measurements during the freeze cycle show a negative slope of 20 mK in 30 min. The blackbody is slowly heated over about 8 h before reaching the melting point and is typically ramped up over night so that it is held just below the melting point the next morning. After the initial heating at 8 A, the melt cycle is begun by increasing the current to 8.5 A until the temperature reaches 1071 °C. The freeze cycle is begun by lowering the blackbody current to 7.95 A. Then the blackbody current is raised to 8.5 A at 1050 °C to begin the melt cycle again.

- 1 cavity
- 2 crucible
- 3 heat pipe
- 4 semicylindrical heater element
- 5 water cooling coil
- 6 graphite rings
- 7 silica glass spacers
- 8 shutter opening device

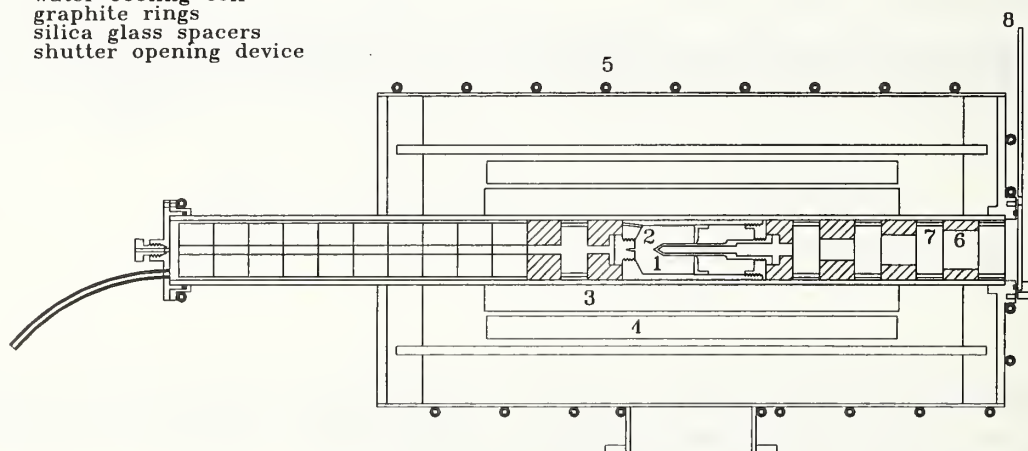


Figure 3. Schematic of gold-point blackbody.


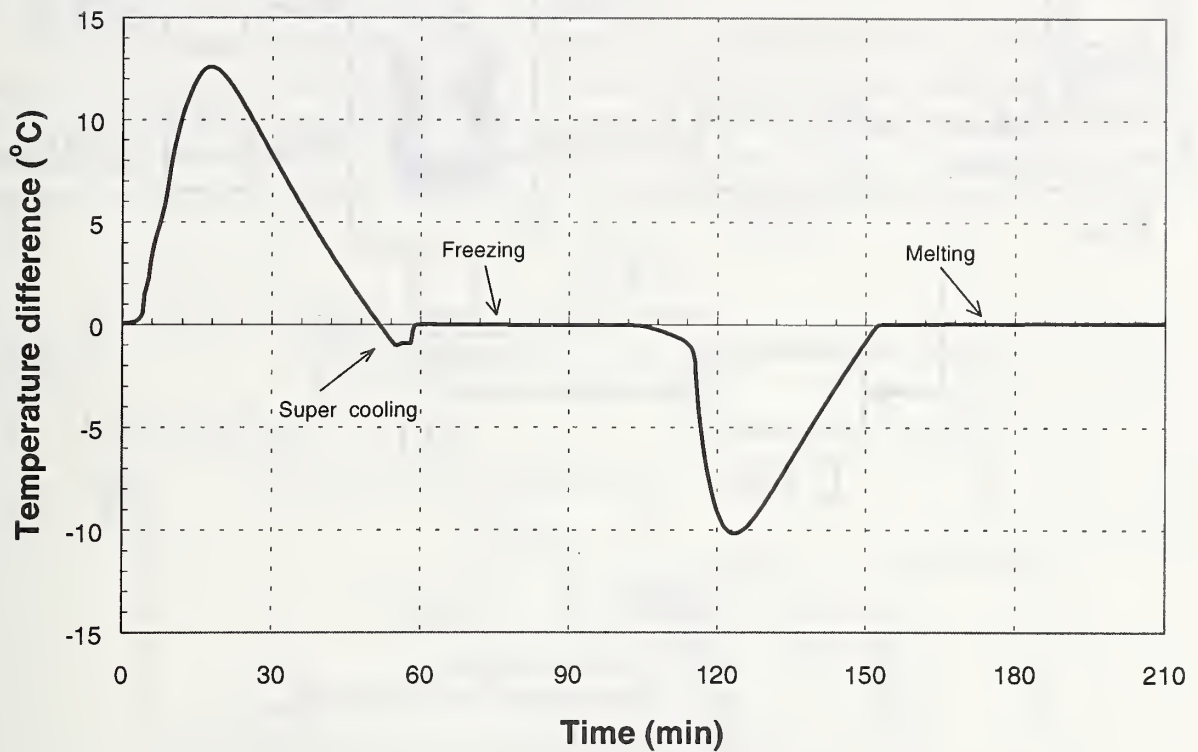


Figure 4. Gold-point blackbody liquid-to-solid phase transition.

The gold-point blackbody control system shown in figure 5 consists of a 286-based personal computer, which regulates the power supply and reads the thermocouple measurements from the digital voltmeter. To monitor the set points, platinum and platinum/10% rhodium (type S) thermocouples are inserted into round double bore alumina tubing and are placed in a silica glass sheathing which makes contact with the rear end of the crucible.

3.2.2 Tungsten ribbon filament lamps

Because of a reproducible radiance temperature versus lamp current relationship, tungsten ribbon filament lamps are used in the radiance temperature scale realization and are issued as radiance temperature standards. Vacuum lamps are typically used for temperatures from 800 °C to 1700 °C, and argon gas-filled lamps for temperatures from 1300 °C to 2300 °C. Gas is introduced into the high temperature lamps to slow down the evaporation rate of the tungsten. The evaporated tungsten atoms collide with the gas atoms, lose their energy, and re-condense on the filament surface before being swept away by convection [9].

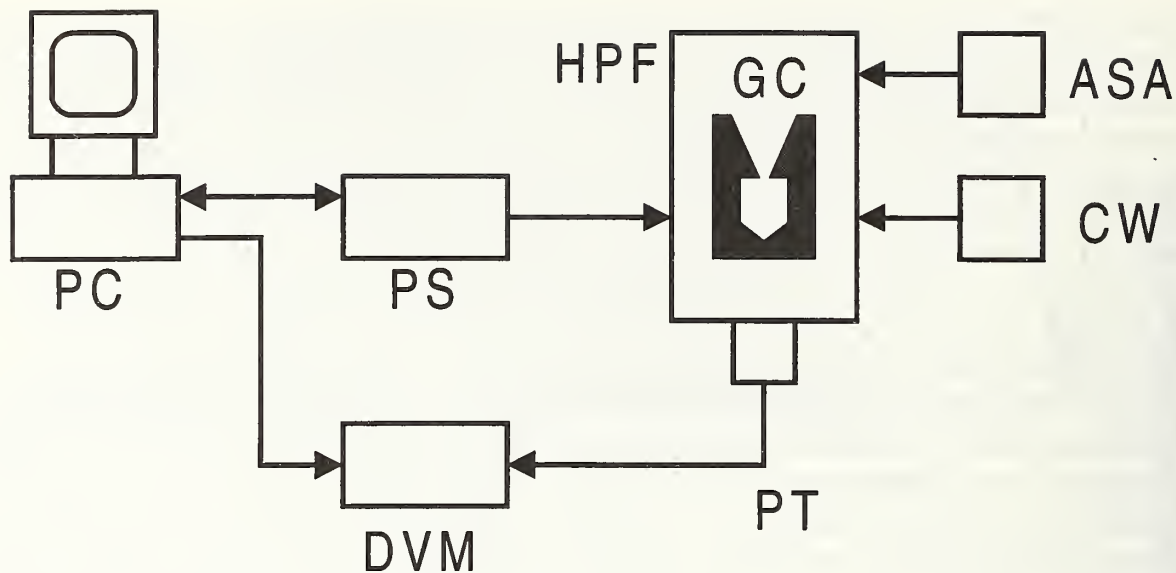

LEGEND

ASA	Argon Supply and Alarm
CW	55 °C Chilled Water/Cold Water Switch
DVM	Digital Voltmeter
GC	Gold-filled Crucible
HPF	Heat Pipe Furnace
PC	Personal Computer
PS	Power Supply
PT	Pt/Pt-10% Rh Thermocouple (Type S)


Figure 5. Gold-point blackbody control system.

The lamps have cylindrical glass envelopes and are operated on direct current with the base facing downward. A computer-stabilized power supply regulates the lamp current to within ± 0.2 mA. A schematic of the lamp current monitoring system is shown in figure 6. The filament is connected to the screw base by wires, which conduct the current and thereby minimize contact problems associated with the screw bases. To achieve the required precise alignment, the lamps are rigidly fastened in source mounts that allow translation along and rotation about three mutually perpendicular axes. Located at the intersection of the center of the lamp filament at the height of the notch opening, the rectangular calibration area on the surface of the lamp filament is 0.6 mm wide by 0.8 mm high. The notch is about midway along one edge of the filament. An alignment mark or an etched arrow is placed on the envelope opposite the side viewed by the pyrometer to permit reproducible angular positioning. Then the lamp is rotated so that the arrowhead is centered at the mouth of the notch.

Some lamps use a pointer that is connected to the filament support to indicate the calibration area. The advantage to this approach is a more uniform filament. However, one

disadvantage is the requirement of predictable filament expansion and contraction for reproducible calibration area indication by the pointer.

Until 1995, the NIST issued General Electric³ type 30A/6V/T24 ribbon filament lamps as radiance temperature standards. However, General Electric stopped manufacturing these lamps in 1990. The lamp had a glass cylindrical envelope that was 300 mm long and 75 mm in diameter and a filament that was 3 mm wide and 50 mm long. Figure 7 shows typical lamp currents as a function of temperature for the General Electric ribbon filament lamps.

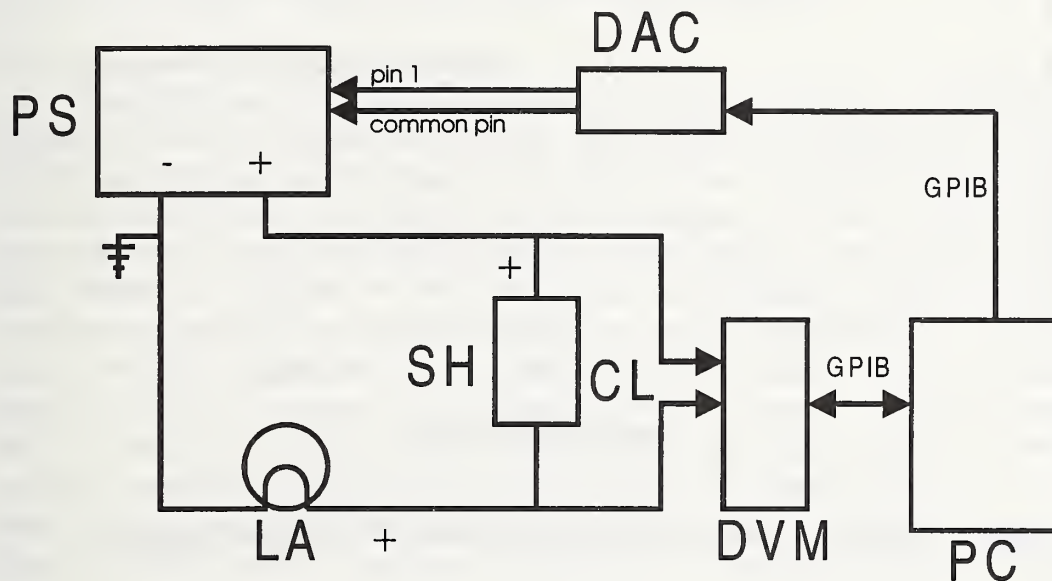

LEGEND

CL	Current Sensing Leads
DAC	Two-channel 12-bit D/A Converter
DVM	Digital Voltmeter
GPIB	GPIB IEEE-488 Interface
LA	Lamp
PC	Personal Computer
PS	Power Supply
SH	0.1 Ω Shunt

Figure 6. Schematic of lamp current monitoring system.

³ Certain commercial equipment, instruments, or materials are identified in this paper for the sole purpose of adequately describing experimental or test procedures. Such identification does not imply recommendation or endorsement by the National Institute of Standards and Technology of a particular product, nor does it imply that the materials or equipment identified are necessarily the best available for the purpose it serves.

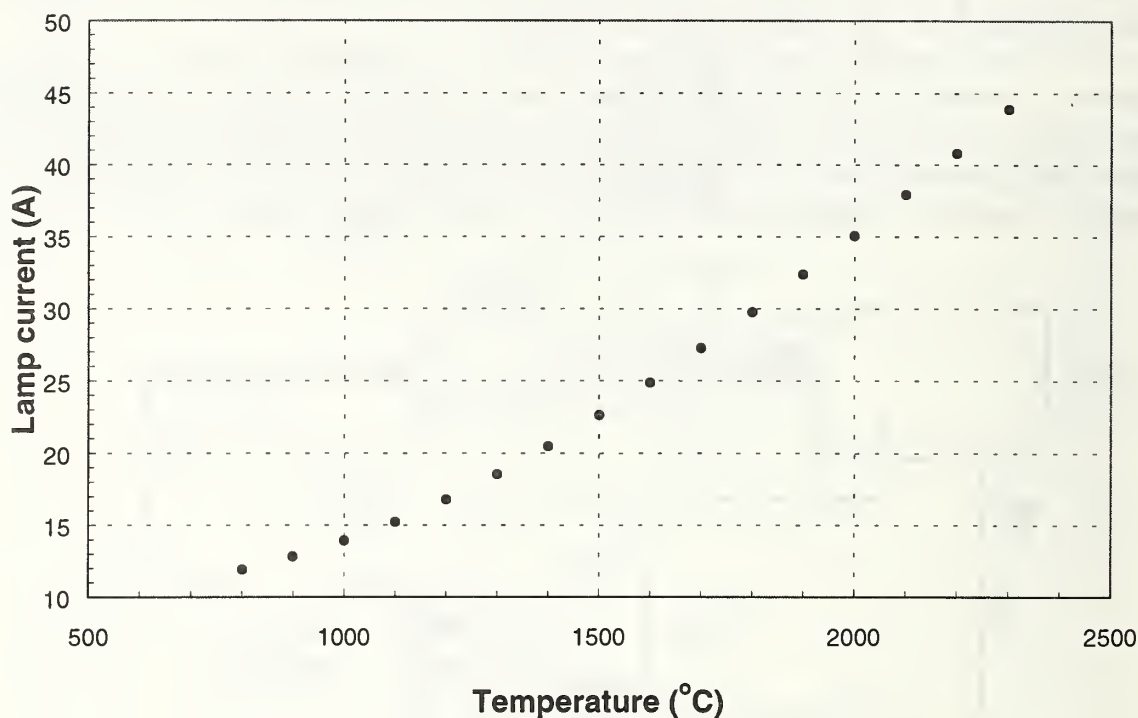


Figure 7. Typical lamp currents for the General Electric 30A/T24/6 ribbon filament lamp.

The NIST is currently issuing U.K. General Electric Company (GEC) lamps, although the lamps are no longer available from this source. The GEC lamps (type 20/G and 20/V) have glass cylindrical envelopes that are 235 mm long with a 64 mm diameter. The 20/V is a vacuum lamp with a filament 1 mm wide by 50 mm long and requires about 10 A dc at 1500 °C. The 20/G is a gas-filled lamp with a 2 mm wide by 30 mm long filament and requires about 20 A dc at 2300 °C.

In 1992, Type TRU 1100-2350 lamps were purchased from the Moscow Lamp Factory (Russia) for testing as radiance temperature standards. Due to a decline in requests for radiance temperature standard lamps, the testing was not done. At 2300 °C, about 25 A dc current is required from the gas-filled lamp, which has a cylindrical envelope, a 2.5 mm wide by 20 mm long filament, and an alignment pointer. The TRU lamps, which have a 35 mm silica glass window that projects about 5 mm from the envelope, were used in addition to the GE lamp in the NIST/VNIOFFI spectral radiance intercomparison [10].

The NIST has also bought type 24/G lamps with double silica glass windows from Polaron Engineering LTD (U. K.) for use as working and check standards. These lamps are similar in design and construction to the GEC lamps.

3.2.3 Working standard lamps

A vacuum tungsten ribbon filament lamp of the Quinn-Lee type [11,12] is used in the temperature scale realization as the secondary temperature standard which maintains the

temperature scale between scale realizations and as the transfer standard for calibration measurements. The temperature of the working standard lamp (serial number SL20) is determined by spectral comparison with the gold-point blackbody. A drawing of this lamp is shown in figure 8. This lamp is operated at a single current (7.7788 A dc) to produce a spectral radiance about eight times higher than that of the gold-point blackbody at 655.3 nm (a radiance temperature of about 1530 K). This lamp is stable to better than 0.1 °C over 100 h when operated under these single-level conditions. A graph of the calibration history of the working standard from August 1989 to February 1996 is shown in figure 9.

3.2.4 Variable-temperature blackbody

The NIST uses a commercial variable-temperature blackbody (VTBB) for its radiance temperature transfer standard (See fig. 10). This VTBB was manufactured by Thermogage Inc. in Frostburg, MD and was supplied with the Type II dual blackbody assembly (2.54 cm or 1 in ID cavity), the 48 kW power supply, a control program, the model 7000-1 (4 range) optical pyrometer, and a digital temperature computer control module. Modifications requested by the NIST include the addition of the water-cooled semi-cylindrical mirrors, and enlarging the extension tube opening to accommodate measurements with the NIST PEP.

The VTBB is operated between 700 °C and 2700 °C. The ranges, which are selected by placing different size apertures on the optical control pyrometer, are the low range (700 °C to 1300 °C), the medium range (1300 °C to 1800 °C), the high range (1800 °C to 2500 °C), and the extra high range (2500 °C to 2700 °C). The electrodes are water-cooled by using a 13 °C (55 °F) chilled water source. Before the blackbody is turned on, argon gas displaces the air in the cavity. When the VTBB cavity is operating, the argon exits from both the front and back extension tubes.

The variable-temperature blackbody control system regulates the blackbody temperature to within ± 0.1 °C (see fig. 11). The NIST measures the uniformity of the center partition of the Thermogage heater elements, and uses those with a spatial uniformity of better than ± 0.2 °C at 2000 °C over an area of 10 mm diameter (see fig. 12). The blackbody heats up to the operating point in less than 30 min and responds quickly to temperature changes (see fig. 13). The day-to-day stability of the blackbody is better than ± 0.3 °C, and the 6-month stability is within ± 2 °C.

The VTBB operates between 800 °C and 2700 °C with an estimated emissivity of 0.99 from analytical modeling. The temperature distribution within the VTBB cavity is less than 50 °C. The blackbody cavity is composed of a single piece of graphite, which is specially tapered for achieving high temperature uniformity. This graphite tube is cylindrically hollow on both ends to form two 2.54 cm diameter cavities with a common partition in the center. One cavity is used for temperature control, and the other cavity is the blackbody source.

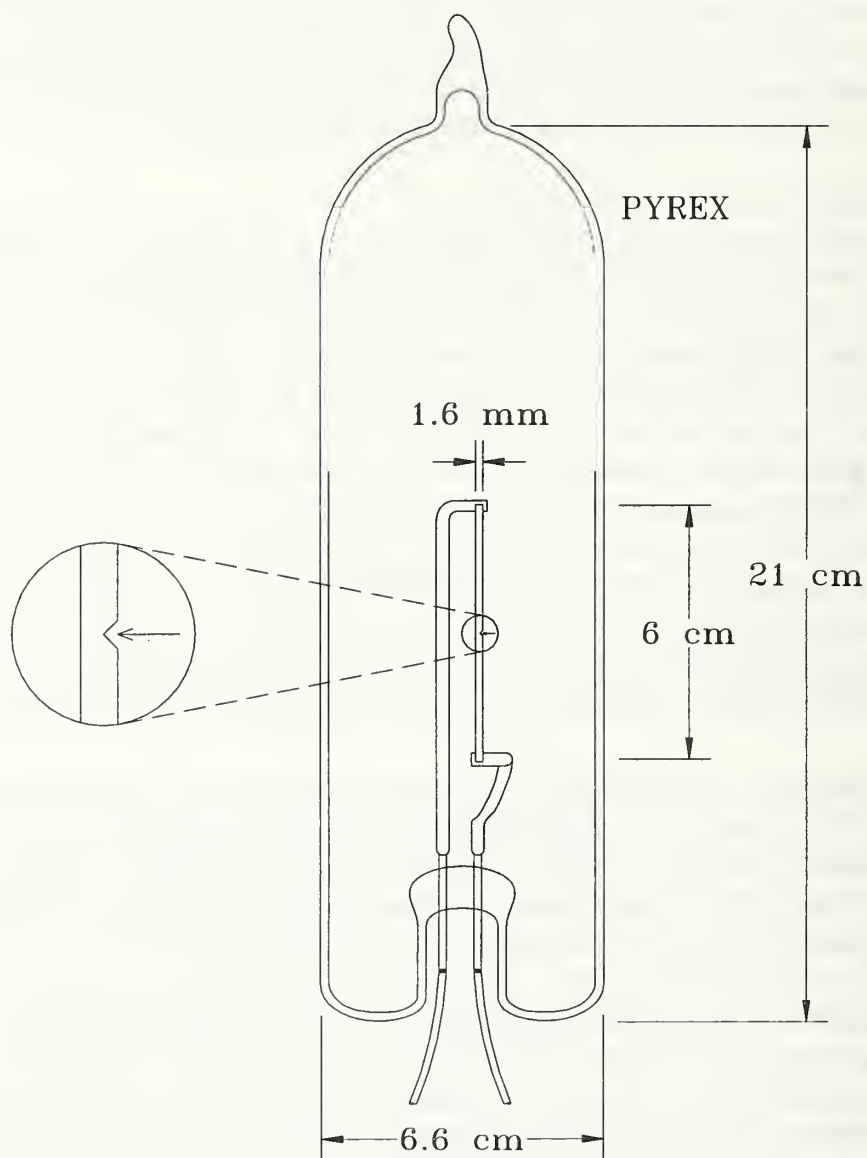


Figure 8. Ribbon filament lamp.

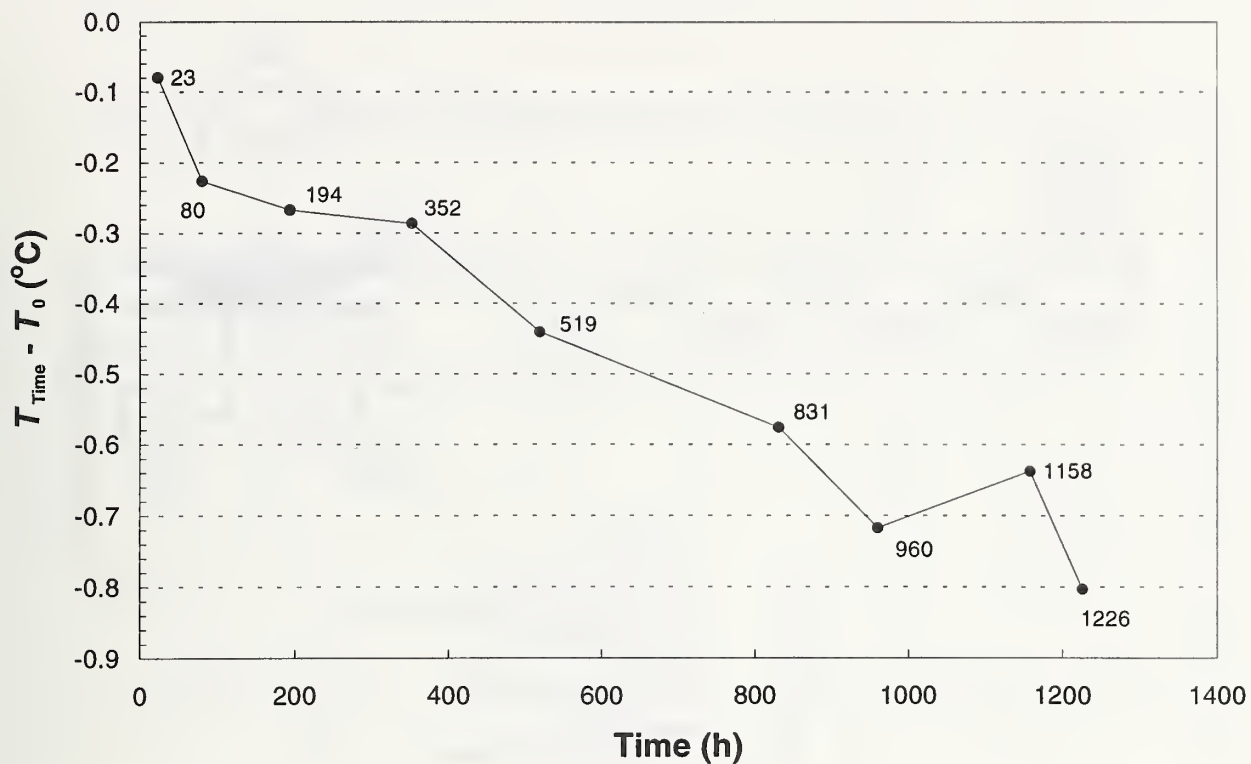


Figure 9. Calibration history of the WS (no. SL20).


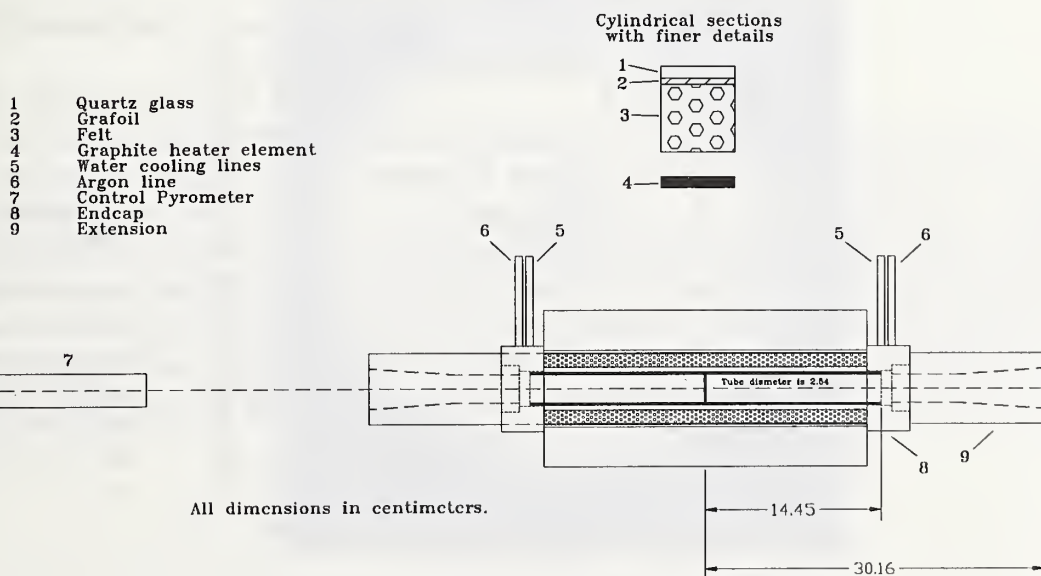


Figure 10. Schematic of variable-temperature blackbody.

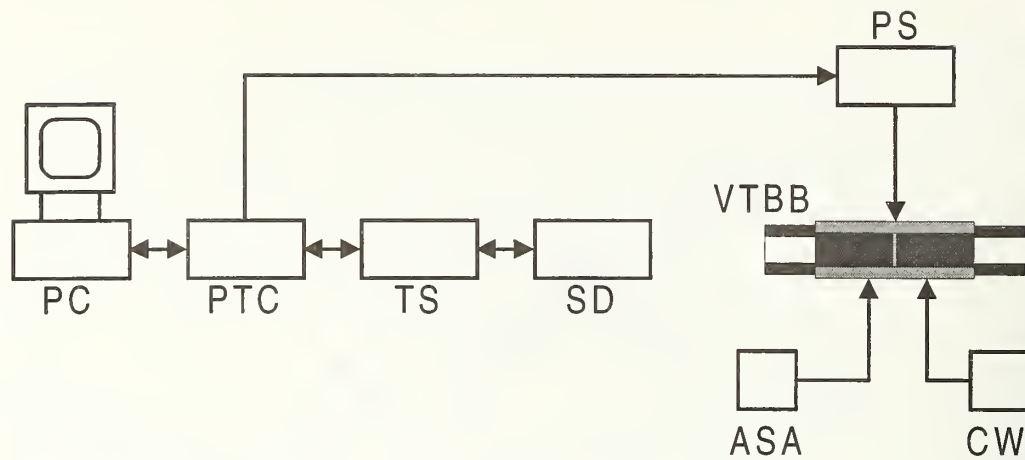

LEGEND

ASA	Argon Supply and Alarm
CW	55 °C Chilled Water/Cold Water Switch
PC	Personal Computer
PS	Power Supply
PTC	PID Temperature Controller
SD	Silicon Detector
TS	Thermal Stabilizer
VTBB	Variable Temperature Blackbody

Figure 11. Variable-temperature blackbody control system.

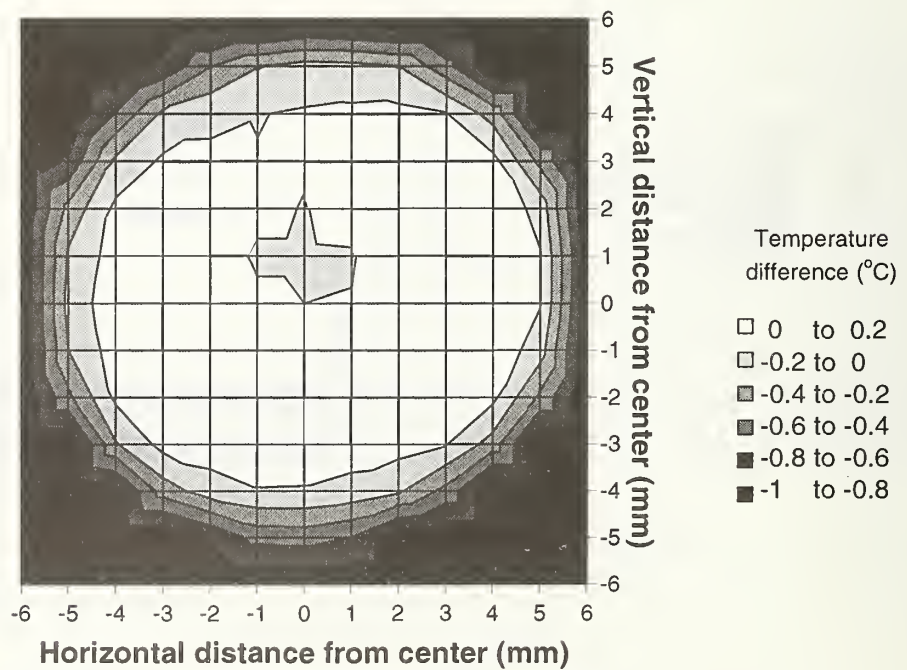


Figure 12. Spatial scan of the variable-temperature blackbody at 2000 °C. The temperature differences are from the center position (0,0).

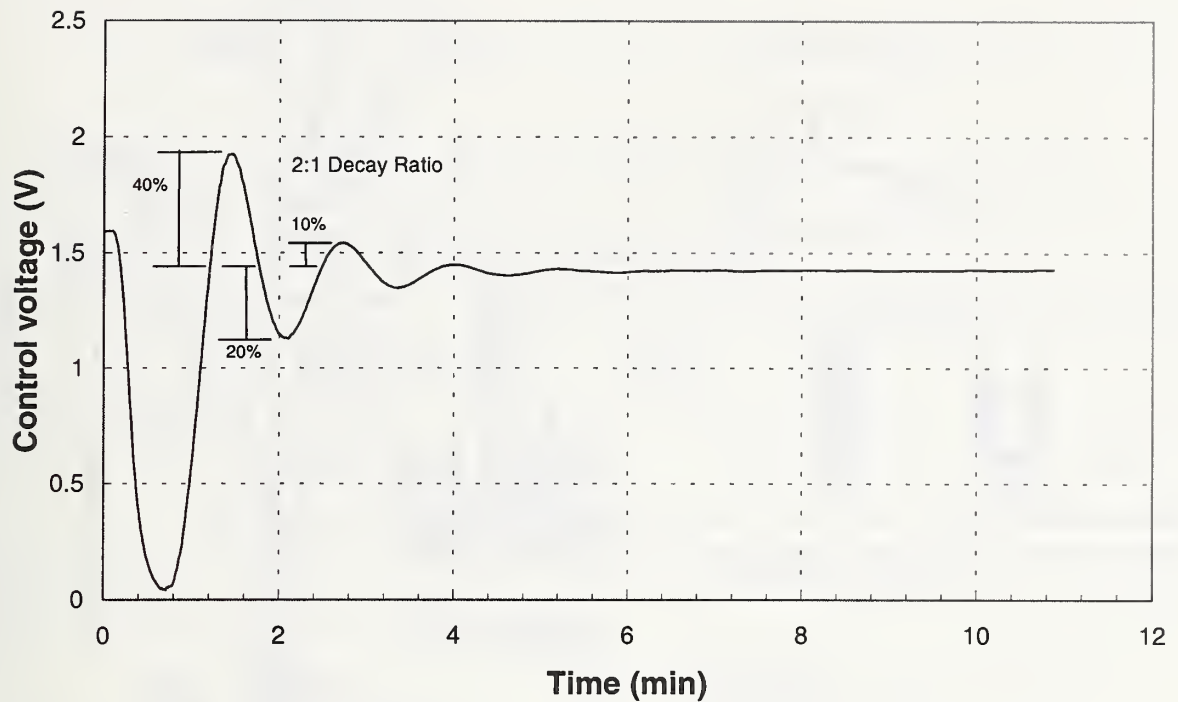


Figure 13. Variable temperature blackbody control. After decreasing the set point by 100 °C, the blackbody stabilizes in less than 10 min.

3.3 NIST photoelectric pyrometer

3.3.1 Measurement system

The PEP is a NIST designed transfer radiometer, which uses refractive optics to image the source onto the detector. The schematic of the NIST PEP is shown in figure 14. A drawing of the measurement system is shown in figure 15. The measurement system is completely automated and controlled by a personal computer, while the laboratory environment is monitored by temperature and relative humidity sensors. Lamps and blackbodies are positioned onto the optical axis of the PEP using a closed-loop motor controller system that allows positioning to within 0.01 mm. The components of the pyrometer are discussed in the next three sections followed by a discussion of the method used to characterize the system in the last three sections.

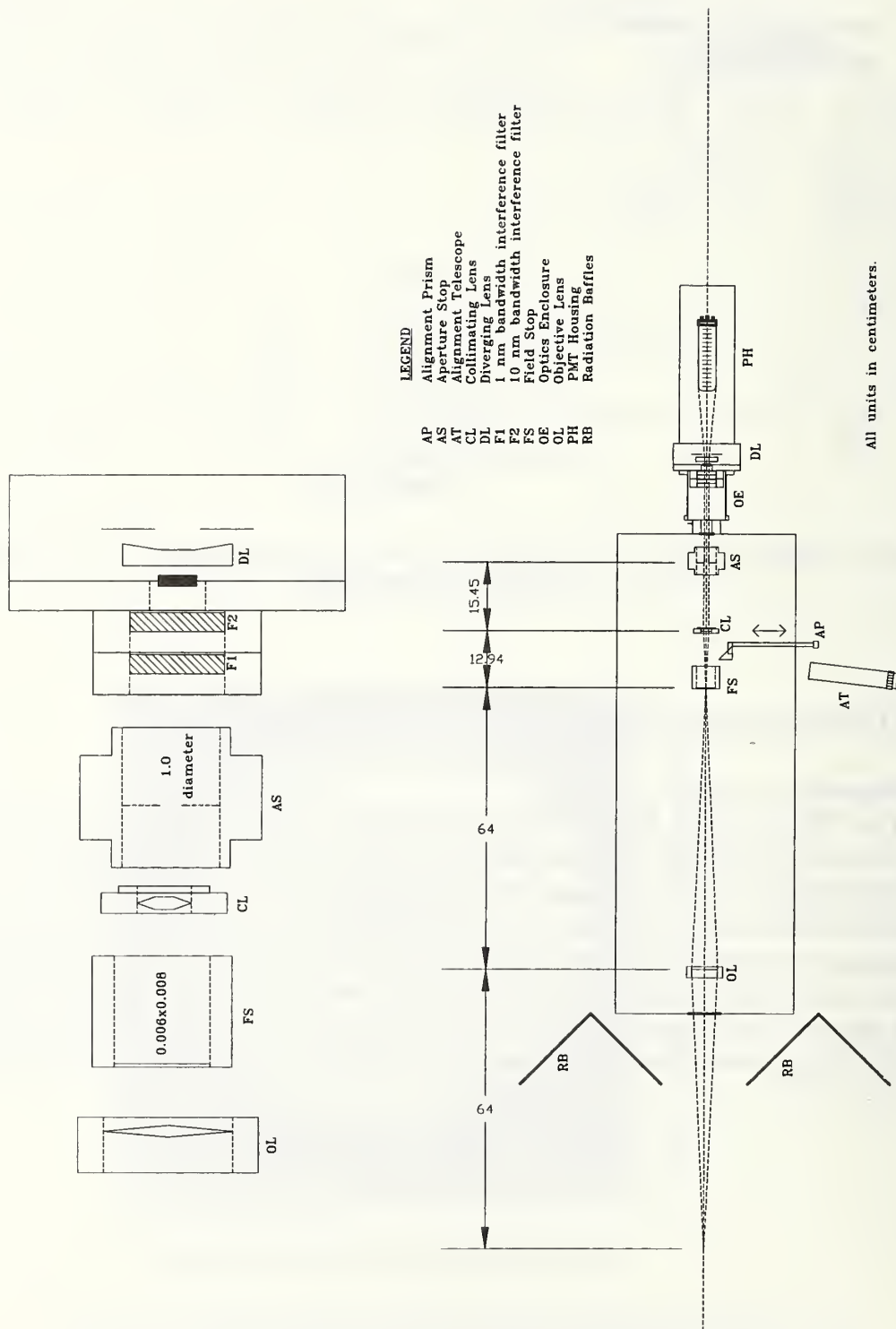


Figure 14. NIST photoelectric pyrometer.

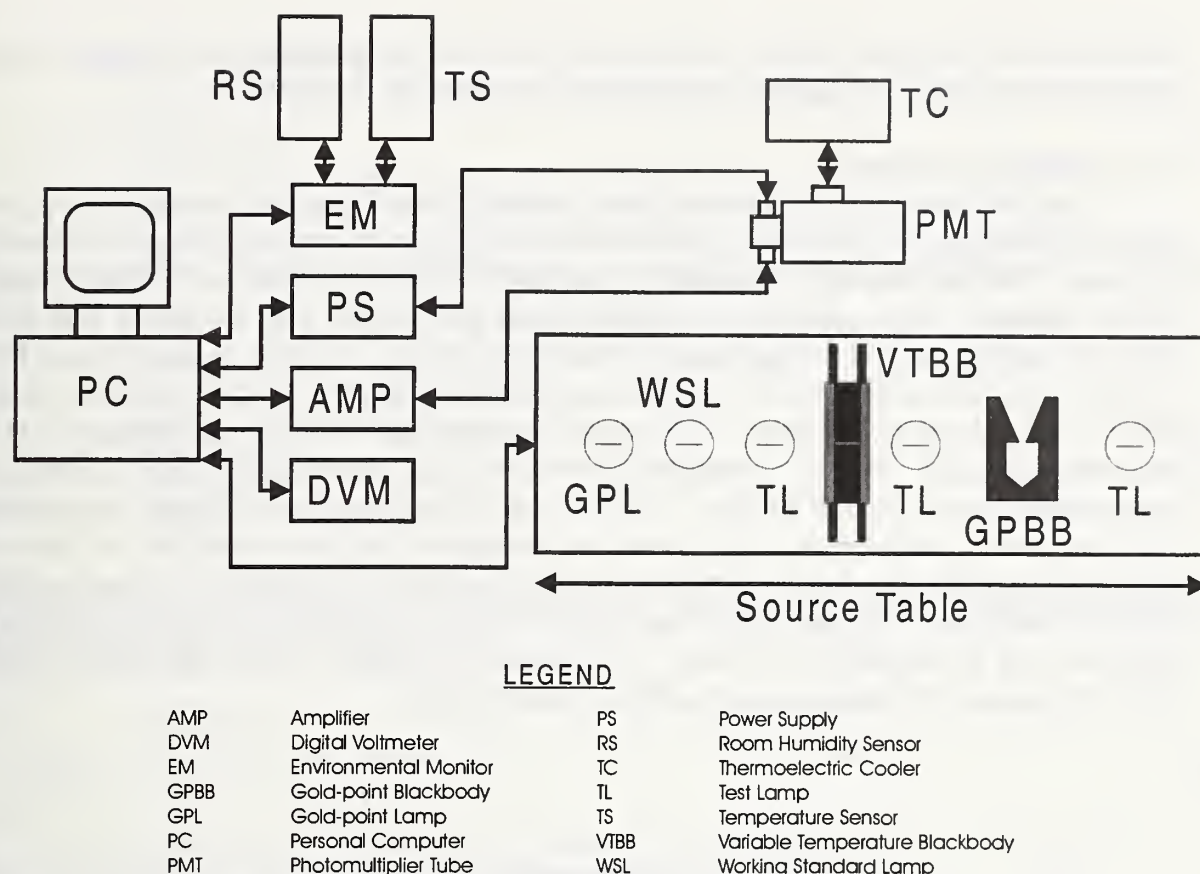


Figure 15. NIST radiance temperature laboratory measurement system.

3.3.2 Optical system

The first element is the objective lens (OL) which is a 5.8 cm diameter bi-convex lens with a focal length of 32 cm. The source is placed at twice the focal length from the lens and is imaged with unit magnification onto the field stop (FS), which is located at the back radius of curvature. The dimensions of the rectangular field stop are 0.6 mm wide by 0.8 mm high. The field stop defines the size and shape of the calibration target area. The next element is the collimating lens (CL), which is a 2.45 cm diameter bi-convex lens with a focal length of 12.94 cm. The field stop is placed at the front focal point of this lens with the result being a collimated beam with a diameter of approximately 1.2 cm. Next the aperture stop (AS) reduces the beam diameter to 1 cm. The beam continues through two interference filters (F1 & F2) through a plano-concave diverging lens that spreads the radiant flux over the photomultiplier tube cathode. The beam completely covers the surface of the cathode.

In all cases, the source overfills the objective lens and the aperture stop that corresponds to a system with an f -number of 13. The solid angle (4.67×10^{-3} sr) used is a cone having a half angle of 0.0772 rad (4.4°), the apex of which is at the center of the field stop.

To align sources, the field stop is removed, and the alignment prism (AP) is placed in the

path behind the field stop mount. The operator views the lamp filament or blackbody aperture and moves the filament or aperture to the target area viewed by the PEP.

3.3.3 Interference filters

The PEP uses two interference filters centered on 655.3 nm to produce a very narrow spectral bandwidth filter pack with transmittances in the wings six orders of magnitude lower than the peak. The first filter is a three-cavity type, and the second is a four-cavity type. Both are 5.08 cm diameter filters made from polished optical quality glass and are coated with a metal film. The regular spectral transmittances of the filters (see fig. 16) were measured from 400 nm to 1000 nm in 0.2 nm steps (with a relative expanded uncertainty of 0.3 % at 655 nm) using the NIST Transfer Spectrophotometer for Regular Transmittance described in reference [13]. The expanded uncertainty of the wavelength calibration of the spectrometer used to make these measurements was 0.5 nm at 655 nm. The first and second filters have spectral bandwidths of 1 nm and 10 nm, respectively. The central wavelength of the filters measured in series was determined to be 654.9 nm with a bandwidth of 2.4 nm when using transmittance data around the peak from 644 nm to 667 nm. For the full spectral range, the measured central wavelength was 836.6 nm with a bandwidth of 9.4 nm. The central wavelength λ_c was determined from the following equation for determining the effective wavelength,

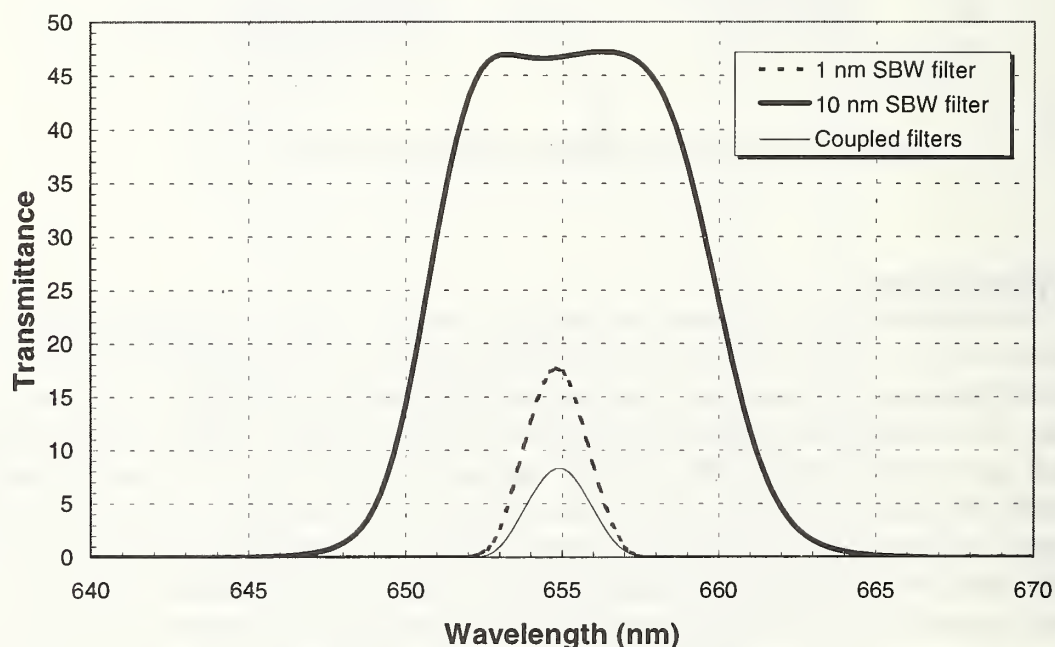


Figure 16. NIST photoelectric pyrometer filter transmittances. The wavelength of the pyrometer is selected by using a 1 nm and a 10 nm spectral bandwidth (SBW) filter.

$$\lambda_c = \frac{\int_{\lambda} \lambda \cdot \tau_{\lambda} \cdot d\lambda}{\int_{\lambda} \tau_{\lambda} \cdot d\lambda}. \quad (17)$$

By approximating the integral with finite sums over equal wavelength bandwidths $d\lambda$, λ_c is found to be

$$\lambda_c = \frac{\sum_{i=400}^{800} (\lambda_i \cdot \tau_i)}{\sum_{i=400}^{800} \tau_i}, \quad (18)$$

where the value of i indicates the wavelength in nm. The spectral bandwidth $\Delta\lambda_s$ was determined from the following integral relationship,

$$\Delta\lambda_s = \frac{\int_{\lambda} \tau_{\lambda} \cdot d\lambda}{\tau_{\text{peak}}}. \quad (19)$$

When the integral is converted into finite sums, the SBW becomes

$$\Delta\lambda_s = \frac{\sum_{i=400}^{800} (\tau_i \cdot \Delta\lambda_i)}{\tau_{\text{peak}}}. \quad (20)$$

3.3.4 Detectors

The PEP uses an 11-stage photomultiplier tube with a quartz window and an extended S-20 spectral response to measure spectral radiances. The detector is housed in a thermoelectrically cooled (-25 °C) chamber. The anode current is amplified and converted to a 0.1 V to 10 V signal by a programmable dc current amplifier and measured using a high accuracy digital voltmeter, capable of integration times ranging from 0.2 s to several minutes.

The signals for the test source and the WS lamp can vary by as much as a factor of 1100. Because several amplifier gain settings are required to cover this dynamic range, all signals are normalized to the 10^8 V/A range. The amplifier normalization factors C_A are determined by measuring the values of the four amplifier feedback resistors. Electrical currents I in the range from 1 nA to 1000 nA are supplied to the current amplifier, and the output voltages V are measured using a digital voltmeter. The resistances R , determined using Ohm's law, $V = I \cdot R$, are R_7 , R_8 , R_9 , and R_{10} , for the (10^7 , 10^8 , 10^9 , and 10^{10}) V/A ranges, respectively. The

normalization factors for each range are shown in table 3.

Table 3. Current amplifier normalization factors

Range [V/A]	Normalization factor
10^7	$C_{A,R7} = R_8/R_7 = 9.9776$
10^8	$C_{A,R8} = R_8/R_8 = 1$
10^9	$C_{A,R9} = R_8/R_9 = 9.9858 \times 10^{-2}$
10^{10}	$C_{A,R10} = R_8/R_{10} = 9.9846 \times 10^{-3}$

3.3.5 Linearity of response

The degree of linearity of the PEP response is determined with a NIST designed automated beam conjoiner [14]. A beam from a constant source is split into two branches, whose fluxes are independently attenuated or blocked before recombination and further attenuation. The flux contribution from both branches is equal to the sum of the fluxes from each branch when measured separately (additivity). The device provides 100 levels of flux ranging over a factor of about 500. The levels are presented in random order to avoid systematic errors, and are interspersed with 25 zero flux levels. A personal computer controls the attenuating filters and records the filter positions and radiometer signals. The data is least-squares fitted to a polynomial response function to determine a correction factor (C_L) by which the radiometer output signal must be multiplied to obtain a quantity proportional to radiant flux. The correction factor is given by

$$C_L = \text{flux} / \text{signal} = (A_0 + A_1 \cdot S + A_2 \cdot S^2) / S, \quad (21)$$

where A_0 , A_1 , and A_2 are the coefficients found from the least squares approximation and S is the signal as defined in eq (1). The current values for the coefficients are listed in table 4. The magnitude of linearity correction factors is graphically shown in figure 17.

The measured instrument response is linear to within $\pm 0.2\%$ for a range of photomultiplier anode currents from 0.1 nA to 500 nA. For currents much less than 0.1 nA, the signal is limited by noise. In order to ensure linearity of response, the current is restricted to a maximum of 500 nA by selection of the appropriate photomultiplier tube voltage. Typical default hardware settings used by the calibration program are listed in table 5.

Table 4. Detector system correction factors and linearity coefficients

Gain [V/A]	C_L	C_L	A_0 [V]	A_1	A_2 [V ⁻¹]
10^7	0.9987 for 500 nA	0.9973 for 1000 nA	-2.1274×10^{-8}	1	-2.6687×10^{-4}
10^8	0.9998 for 50 nA	0.9996 for 100 nA	-2.1308×10^{-6}	1	-3.8995×10^{-5}
10^9	1.0002 for 5 nA	1.0003 for 10 nA	-8.9698×10^{-9}	1	-3.1787×10^{-5}

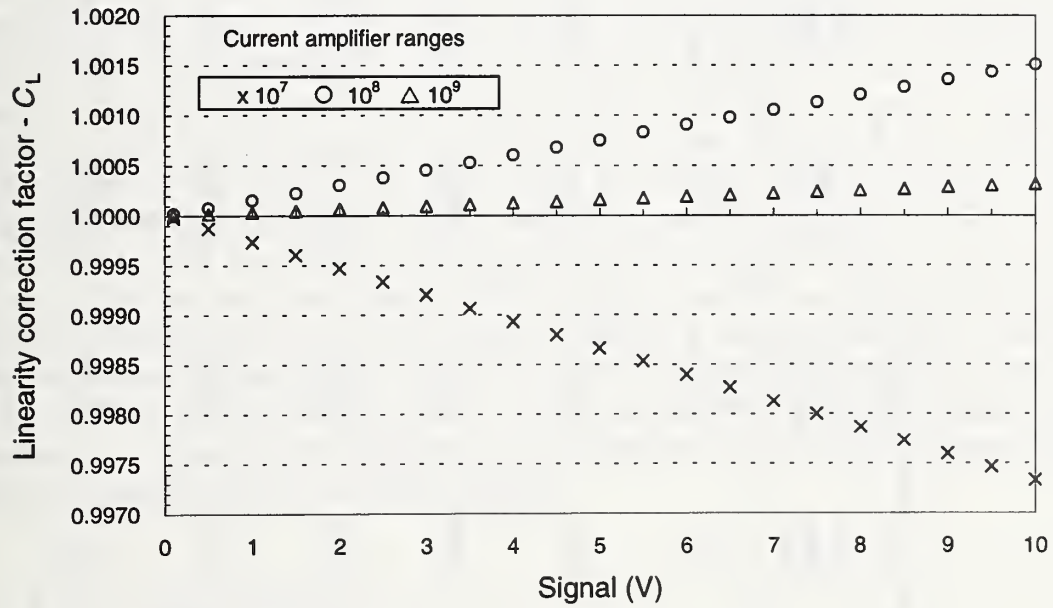


Figure 17. The linearity correction factors for the PEP.

Table 5. Default hardware settings used by calibration program

Temp [°C]	SL20 Signal [V]	Gain	Range Factor	SL20 Range Correction	L_1 (SL20) [W·cm ⁻³ ·sr ⁻¹]	L_2 (TEST) [W·cm ⁻³ ·sr ⁻¹]	Ratio Range Correction	Test Range Correction	Gain	Range Factor	Test Signal [V]	PMT HV [V]	Wait Time [min]
800	0.95	7	9.9776	9.4787	575.055	1.3015×10^0	2.2633×10^{-3}	0.0215	9	0.09986	0.2148	730	20
1100	0.95	7	9.9776	9.4787	575.055	1.1324×10^2	1.9693×10^{-1}	1.8666	8	1	1.8666	730	10
1500	1	9	0.099858	0.0999	575.055	4.1608×10^3	7.2355×10^0	0.7225	9	0.09986	7.2355	439	10
1900	1	9	0.099858	0.0999	575.055	4.0567×10^4	7.0544×10^0	7.0444	8	1	7.0444	439	5
2300	0.3	9	0.099858	0.0300	575.055	1.9487×10^5	3.3887×10^0	10.1517	7	9.9776	1.0174	385	5
3000	0.3	9	0.099858	0.0300	575.055	1.2079×10^6	2.1004×10^0	62.9232	7	9.9776	6.3064	385	5

3.3.6 Wavelength calibration

The mean effective wavelength of the PEP is calculated for temperature from 800 °C to 3000 °C by measuring its relative spectral response in the NIST Visible/Near Infrared Spectral Comparator Facility [15]. The collimating lens, aperture stop, optics enclosure, and PMT system of the PEP were set up in the Spectral Comparator Facility and the relative spectral response was measured while maintaining the same system geometry as in the RTCL. The PEP was measured from 400 nm to 900 nm in 5 nm steps and from 630 nm to 670 nm in 0.2 nm steps. The relative spectral response was calculated from the following equation,

$$R_{\text{PEP}}(\lambda) = \frac{\frac{S_{\text{PEP}}(\lambda)}{S_{\text{m1}}(\lambda)}}{\frac{S_{\text{W}}(\lambda)}{S_{\text{m2}}(\lambda)}} \cdot R_{\text{W}}(\lambda) \cdot \frac{G_{\text{W}}}{G_{\text{PEP}}}, \quad (23)$$

where $R_{\text{PEP}}(\lambda)$ is the relative spectral response of the PEP in [A/W] as defined in eq (2), $S_{\text{PEP}}(\lambda)$ is the signal from the PEP in [V], $S_{\text{m1}}(\lambda)$ is the signal from the monitor detector simultaneous with $S_{\text{PEP}}(\lambda)$, $S_{\text{W}}(\lambda)$ is the signal from the detector working standard (DWS), $S_{\text{m2}}(\lambda)$ is the signal from the monitor detector simultaneous with $S_{\text{W}}(\lambda)$, $R_{\text{W}}(\lambda)$ is the absolute spectral response of the DWS, G_{W} is the calibrated amplifier gain for the DWS in [V/A], and G_{PEP} is the calibrated amplifier gain for the PEP. The PEP results are relative spectral responses because the measurement was performed in the overfill mode that resulted in a slightly different spot size. From eq (18), the CW was determined to be 655.0 nm and the SBW was calculated to be 5.0 nm from eq (20) for data over the 400 nm to 900 nm range. The CW was 655.3 nm with a SBW of 4.2 nm for the 630 nm to 670 nm data. Relative to the peak value, the measured relative spectral responsivity values decrease to about 10^{-4} at 5 nm from the central wavelength and to about 10^{-6} at 50 nm as illustrated in figure 18. Spectral radiance as defined in eq (4) is a function of a single wavelength, λ . However, the PEP has a finite bandwidth, and the PEP measures radiance as defined by an integral form of the Wien equation,

$$L = \int_{400}^{900} \frac{c_1}{\lambda^5 \cdot \exp(c_2 / (\lambda \cdot T))} \cdot R_{\lambda} \cdot d\lambda, \quad (24)$$

which could be used in the ITS-90 definition of eq (14). In the RTCL, the single wavelength $\lambda_{T_1-T_2}$ used in eq (4) is defined as follows,

$$\lambda_{T_1-T_2} = \frac{c_2}{\ln(L_2 / L_1)} \cdot \left(\frac{1}{T_1} - \frac{1}{T_2} \right), \quad (25)$$

where $\lambda_{T_1-T_2}$ is the mean effective wavelength, L_1 is the radiance of the test lamp at temperature T_1 , and L_2 is the radiance of the WS at temperature T_2 .

All of the temperature measurements are performed for a mean effective wavelength of 655.3 nm, and these corrections are calculated for the VTBB with an emissivity of 0.99. The spectral radiance $L_\lambda(\lambda, \epsilon_\lambda, T)$ of the blackbody is defined by eq (5) to be a function of the wavelength, spectral emissivity, and the temperature, while the temperature $T(\lambda, L_\lambda, \epsilon_\lambda)$ is defined by eq (6) to be a function of the wavelength, spectral emissivity, and the spectral radiance. The following steps are performed to calculate the corrections at 900 nm for $T_0 = 1073$ K.

- ⇒ Calculate $L_0 = 1.2856 \text{ W} \cdot \text{cm}^{-3} \cdot \text{sr}^{-1}$, where $\lambda = 655.3 \text{ nm}$, $\epsilon_\lambda = 1$, and $T = T_0$
- ⇒ Calculate $T_1 = 1073.53 \text{ K}$, where $\lambda = 655.3 \text{ nm}$, $\epsilon_\lambda = 0.99$, and $L_\lambda = L_0$
- ⇒ Calculate $L_1 = 68.31 \text{ W} \cdot \text{cm}^{-3} \cdot \text{sr}^{-1}$, where $\lambda = 900 \text{ nm}$, $\epsilon_\lambda = 0.99$, and $T = T_1$
- ⇒ Calculate $T_2 = 1072.80 \text{ K}$, where $\lambda = 900 \text{ nm}$, $\epsilon_\lambda = 1$, and $L_\lambda = L_1$

The corrections $T_2 - T_0$ that are added to the radiance temperature at 655.3 nm to determine the radiance temperatures at 900 nm and 1000 nm are given in table 6. These corrections are equivalent to a spectral radiance correction of -0.27 % at 900 nm and -0.34 % at 1000 nm.

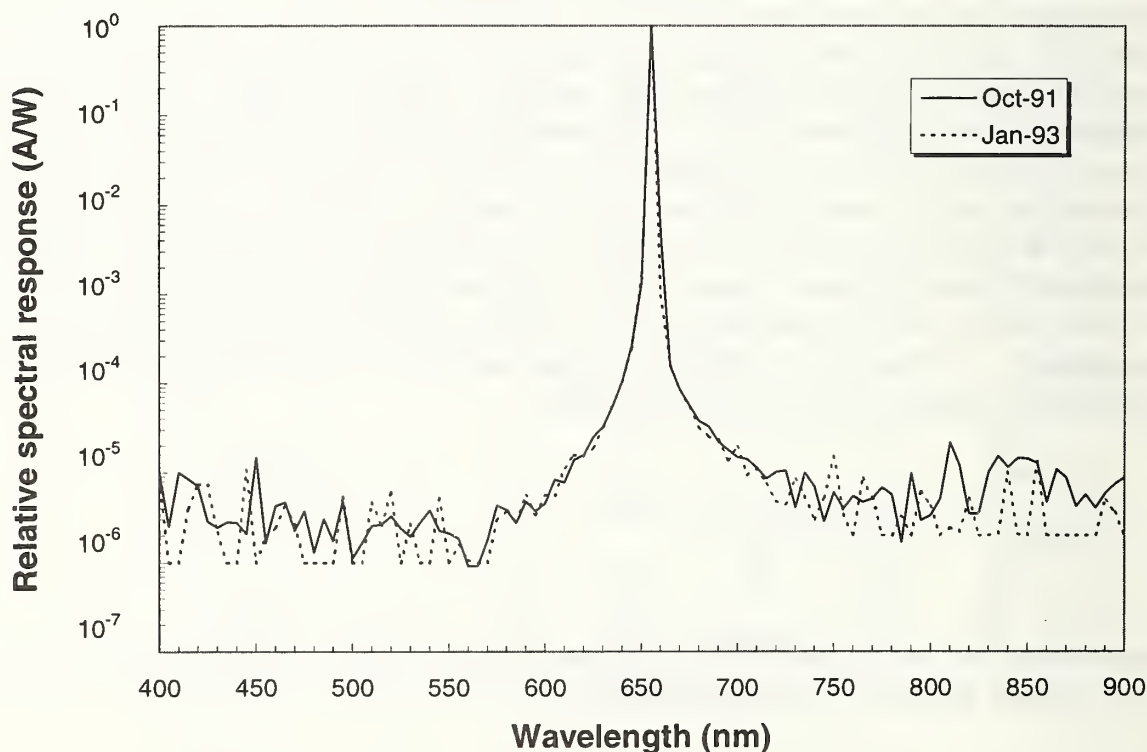


Figure 18. Relative spectral response of the NIST photoelectric pyrometer.

Table 6. Corrections (°C) to the radiance temperature at 900 nm and at 1000 nm

Blackbody temperature (°C)	900 nm	1000 nm
800	-0.20	-0.27
1100	-0.32	-0.45
1500	-0.53	-0.75
1900	-0.80	-1.13
2300	-1.12	-1.58
3000	-1.82	-2.55

3.3.7 Size of source effect

Measurements are performed to determine how much of the measured flux is from sources other than the target. This additional flux from outside the target area is scattered into the measured beam by the optical system. The measurement is performed to determine the size of source effect or the sensitivity of the PEP when viewing a 0.6 mm wide by 0.8 mm high target area on sources of various sizes. The magnitude of the size of source correction is determined by measuring the spectral radiance of a large uniform diffuse source with various apertures in front of the source. The large uniform diffuse source designed at NIST is used to measure the source of source effect. A 1 kW frosted quartz-halogen lamp is placed in a 20 cm by 23 cm by 20 cm vented housing. A lens focuses the lamp onto an opening in the housing that is covered by a diffuser. The apertures are on an aperture slide for quick positioning and reproducibility. The aperture sizes measured are a 1.4 mm by 25.4 mm slit, a 3 mm by 25.4 mm slit, a 15 mm diameter hole, and 25.4 mm diameter hole which approximate the sizes of the WS filament, the TL filament, the Au blackbody opening, and the BB opening. The ratio r_1 is the measurement of the ratio of the WS to the AuBB in eq (9). The size of source correction is calculated from the ratio of the measurements of the 15 mm diameter hole to the 1.4 mm by 25.4 mm slit. The spectral radiance correction was 0.09 %, which corresponds to a correction factor of 1.0009 ± 0.0006 and a temperature correction of 0.10 °C. A negligible difference was observed for the comparison of the measurements of the 1.4 mm by 25.4 mm slit and the 3 mm by 25.4 mm slit; therefore, the correction factor for ratio r_2 is 1. The size of source correction for ratios r_3 and r_4 is -0.13 % which corresponds to a correction factor of 0.9987 ± 0.0002 . The correction was calculated from the ratio of the measurements of the 1.4 mm by 25.4 mm slit to the 25.4 mm diameter hole. The magnitude of the temperature correction varies from -0.07 °C at 800 °C to -0.52 °C at 2700 °C. The size of source corrections are shown in figure 19. Because of flux contributions from outside the target area, the larger area BB results in a higher temperature. A negative correction is added to subtract the additional flux.

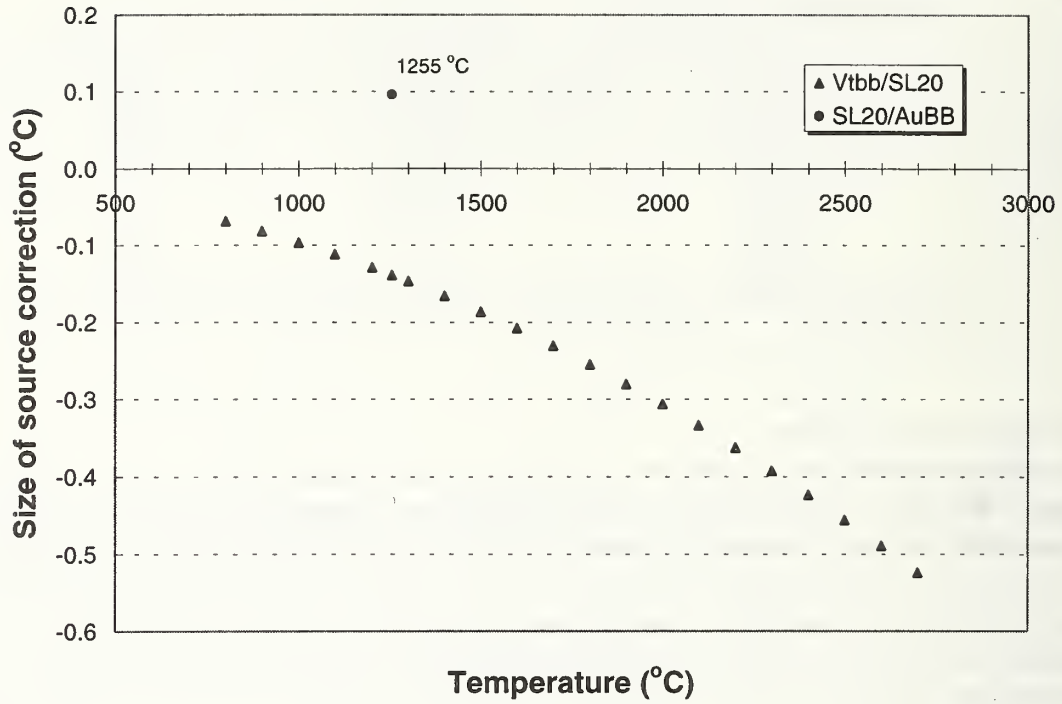


Figure 19. The size of source corrections for the variable temperature blackbody to the working standard lamp comparison and the working standard lamp to the gold-point blackbody.

3.4 Radiance temperature scale uncertainties

From the laws of statistical theory, when the covariances of the independent variables are assumed to be negligible, the propagation of standard uncertainty $u(y)$ for a variable y ,

$$y = f(x_1, x_2, x_3, \dots, x_N), \quad (26)$$

can be defined in terms of the following sum,

$$u(y) = \left[\sum_{i=1}^N \left(\frac{\partial f}{\partial x_i} \cdot u(x_i) \right)^2 \right]^{1/2}, \quad (27)$$

where $u(x_i)$ are the standard uncertainties. The uncertainty in the WS lamp spectral radiance,

$$\begin{aligned}
u_0(L_{\text{WS}}) = & \left[\left(\frac{\partial L_{\text{WS}}}{\partial n_\lambda} \cdot u(n_\lambda) \right)^2 + \left(\frac{\partial L_{\text{WS}}}{\partial \lambda} \cdot u(\lambda) \right)^2 + \left(\frac{\partial L_{\text{WS}}}{\partial T_{\text{Au}}} \cdot u(T_{\text{Au}}) \right)^2 \right. \\
& \left. + \left(\frac{\partial L_{\text{WS}}}{\partial c_2} \cdot u(c_2) \right)^2 + \sum_{i=1}^{11} \left(\frac{\partial L_{\text{WS}}}{\partial x_i} \cdot u(x_i) \right)^2 \right]^{1/2},
\end{aligned} \tag{28}$$

can be calculated by the root sum of squares (RSS) of the products of partial derivatives with their respective uncertainties. This uncertainty $u_0(L_{\text{WS}})$ can be derived directly from the measurement equation in eq (16) by using the propagation of standard uncertainty relationship in eq (27). The partial derivatives are those of the spectral radiance with respect to the independent variables in eq (16) and are listed in eqs (29) through (33). The eleven variables in the summation of eq (28) are the following variables: ϵ_{Au} , $c_{1\text{L}}$, r_1 , $C_{\text{A,WS}}$, $C_{\text{L,WS}}$, $C_{\text{S,WS}}$, G_{WS} , $C_{\text{A,Au}}$, $C_{\text{L,Au}}$, $C_{\text{S,Au}}$, or G_{Au} . Since all of the measurements in the RTCL are based on the WS lamp, the uncertainty analysis of each calibration service (tungsten ribbon filament lamps, disappearing filament optical pyrometers, infrared thermometers) will be completed in Sections 4, 5, and 6 using the uncertainty in $u_0(L_{\text{WS}})$.

The measurement equation for the radiance temperature scale from the gold-point blackbody to the WS lamp can be represented by eq (16). The partial derivatives of the WS spectral radiance in eq (16) with respect to the variables, n_λ , λ , T_{Au} , and c_2 , are

$$\frac{\partial L_{\text{WS}}}{\partial n_\lambda} = \frac{L_{\text{WS}}}{n_\lambda} \cdot \left[\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}} \cdot \frac{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}}\right)}{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}}\right) - 1} - 2 \right], \tag{29}$$

$$\frac{\partial L_{\text{WS}}}{\partial \lambda} = \frac{L_{\text{WS}}}{\lambda} \cdot \left[\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}} \cdot \frac{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}}\right)}{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}}\right) - 1} - 5 \right], \tag{30}$$

$$\frac{\partial L_{\text{WS}}}{\partial T_{\text{Au}}} = \frac{L_{\text{WS}}}{T_{\text{Au}}} \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}} \cdot \frac{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}}\right)}{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{Au}}}\right) - 1}, \text{ and} \tag{31}$$

$$\frac{\partial L_{WS}}{\partial c_2} = \frac{-L_{WS}}{c_2} \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{Au}} \cdot \frac{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{Au}}\right)}{\exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{Au}}\right) - 1}. \quad (32)$$

The expression in eq (28),

$$\left| \frac{\partial L_{WS}}{\partial x} \right| = \frac{L_{WS}}{x}, \quad (33)$$

represents the partial derivatives of L_{WS} with respect to x , where x is one of the eleven variables in eq (28).

Table 7 details typical values of the variables in eq (16), while table 8 summarizes the uncertainties in eq (28). The uncertainty $u(T_{Au})$ in table 8 is derived from the NIST radiometric determination of the freezing point of gold [16]. The combined uncertainty $u(L_{WS})$ in the spectral radiance of the WS lamp can then be calculated by calculating the RSS from the following expression,

$$\frac{u(L_{WS})}{L_{WS}} = \left[\left(\frac{u_0(L_{WS})}{L_{WS}} \right)^2 + \left(\frac{u(\text{DMM})}{L_{WS}} \right)^2 + \left(\frac{u(\text{LTC})}{L_{WS}} \right)^2 + \left(\frac{u(\text{LD})}{L_{WS}} \right)^2 \right]^{1/2}, \quad (34)$$

where $u(\text{DMM})$ is the digital multimeter uncertainty, $u(\text{LTC})$ is the lamp current uncertainty, and $u(\text{LD})$ is the lamp drift uncertainty. From tables 7 and 8, a typical relative expanded uncertainty in the WS spectral radiance is equal to $(2.870 \text{ W} \cdot \text{m}^{-2} \cdot \mu\text{m}^{-1} \cdot \text{sr}^{-1}) / (569.9 \text{ W} \cdot \text{m}^{-2} \cdot \mu\text{m}^{-1} \cdot \text{sr}^{-1})$, or 0.50 %. The final value of the WS spectral radiance will be utilized in Sections 4 through 6 to continue the uncertainty analyses for each type of calibration service.

3.5 Quality control

The Radiance Temperature Measurements Calibration Service ensures quality measurements through quality system implementation, internal quality control procedures, periodic audits, and interlaboratory comparisons. The Quality System Manual for the Radiance Temperature Measurements Calibration Service is based on the ANSI/NCSL Z540-1-1994 standard [17] and the ISO/IEC Guide 25 [18]. The Quality System Manual is organized into three sections that form a pyramidal shape. The top section is the Optical Technology Division Quality System Policies, the middle section contains the Division Procedures that uniformly apply to all the calibration services, and the bottom section is the Calibration Service Procedures for each calibration service. Additional information on the Divisions Quality System can be found in reference [19].

Table 7. Typical values of WS lamp variables and parameters

Variable	Symbol	Value
Refractive index	n_λ	1.00028
Wavelength	λ	655.3 nm
Freezing temperature of gold	T_{Au}	1337.33 K
Second radiation constant	c_2	14387.69 $\mu\text{m}\cdot\text{K}$
Emissivity of Au	ϵ_{Au}	0.9999
First radiation constant	c_{1L}	$1.191 \times 10^8 \text{ W}\cdot\mu\text{m}^4\cdot\text{m}^{-2}$
Ratio of WS signal to Au signal	r_1	7.764
WS amplifier calibration correction	$C_{\text{A,WS}}$	1
WS linearity correction	$C_{\text{L,WS}}$	1
WS size of source correction	$C_{\text{S,WS}}$	1.0009
WS amplifier gain	G_{WS}	$1 \times 10^8 \text{ V}\cdot\text{A}^{-1}$
Au amplifier calibration correction	$C_{\text{A,Au}}$	1
Au linearity correction	$C_{\text{L,Au}}$	1
Au size of source correction	$C_{\text{S,Au}}$	1
Au amplifier gain	G_{Au}	$1 \times 10^8 \text{ V}\cdot\text{A}^{-1}$
WS spectral radiance	L_{WS}	$569.9 \text{ W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$

Table 8. Uncertainty budget for the NIST WS lamp spectral radiance realization

Uncertainty component	Symbol	Expanded Uncertainty ($k = 2$)	
		Type A	Type B
Refractive index	$u(n_\lambda)$		0.00002
Wavelength	$u(\lambda)$		0.2 nm
Freezing temperature of gold	$u(T_{\text{Au}})$		0.23 K
Second radiation constant	$u(c_2)$		0.24 $\mu\text{m}\cdot\text{K}$
Emissivity of Au	$u(\varepsilon_{\text{Au}})$		0.0002
First radiation constant	$u(c_{1\text{L}})$		440 $\text{W}\cdot\mu\text{m}^4\cdot\text{m}^{-2}$
Ratio of WS signal to Au signal	$u(r_1)$	0.006	
WS amplifier calibration correction	$u(C_{\text{A,WS}})$	0.0001	
WS linearity correction	$u(C_{\text{L,WS}})$	0.001	
WS size of source correction	$u(C_{\text{S,WS}})$	0.0006	
WS amplifier gain	$u(G_{\text{WS}})$		0 $\text{V}\cdot\text{A}^{-1}$
Au amplifier calibration correction	$u(C_{\text{A,Au}})$	0.0001	
Au linearity correction	$u(C_{\text{L,Au}})$	0.001	
Au size of source correction	$u(C_{\text{S,Au}})$	0.0002	
Au amplifier gain	$u(G_{\text{Au}})$		0 $\text{V}\cdot\text{A}^{-1}$
WS lamp spectral radiance	$u_0(L_{\text{WS}})$		2.754 $\text{W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
Digital voltmeter	$u(\text{DMM})$		0.000 $\text{W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
WS current	$u(\text{LTC})$		0.570 $\text{W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
WS drift	$u(\text{LD})$	0.570 $\text{W}/\text{m}^2/\mu\text{m}/\text{sr}$	
NIST WS lamp spectral radiance	$u(L_{\text{WS}})$	2.870 $\text{W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$	
NIST WS lamp temperature	$u(T_{\text{WS}})$	0.536 K	

The calibrations resulting from the processes described in this document are derived in principle from a complete realization of the radiance temperature scale, and therefore are not dependent upon the maintenance of in-house standards. In practice, the realization often omits the GPBB comparison and relies upon the spectral radiance that is maintained on the stable working standard lamp. The calibration history of the lamp compared to the GPBB is monitored to confirm the stability of the lamp (see fig. 9). The working standard lamp is compared to the gold point blackbody at least twice a year. Lamp and radiation thermometer check standards are periodically measured to ensure the accuracy of the measurements. The calibration history of the ribbon filament lamp check standards (Serial numbers P51 and P74) are shown in figure 20, while the calibration history of the radiation thermometer check standard (Minolta Cyclops 52) is shown in figure 21. The calibration process is further monitored by using control limits in the data acquisition and data reduction programs. A check-list is used to track the status of the test unit throughout the calibration process. A calibration summary sheet is created for each calibrated item to record customer information, measurements performed, test unit information, environmental conditions, filenames, and computer programs used for data acquisition and reduction. The calibration report is carefully reviewed by a party who was not involved in producing the values.

Audits are performed annually by the Division's quality control team. This calibration service participates in bilateral and multilateral laboratory comparisons of radiance temperature with various national laboratories usually involving guest scientists [20,21].

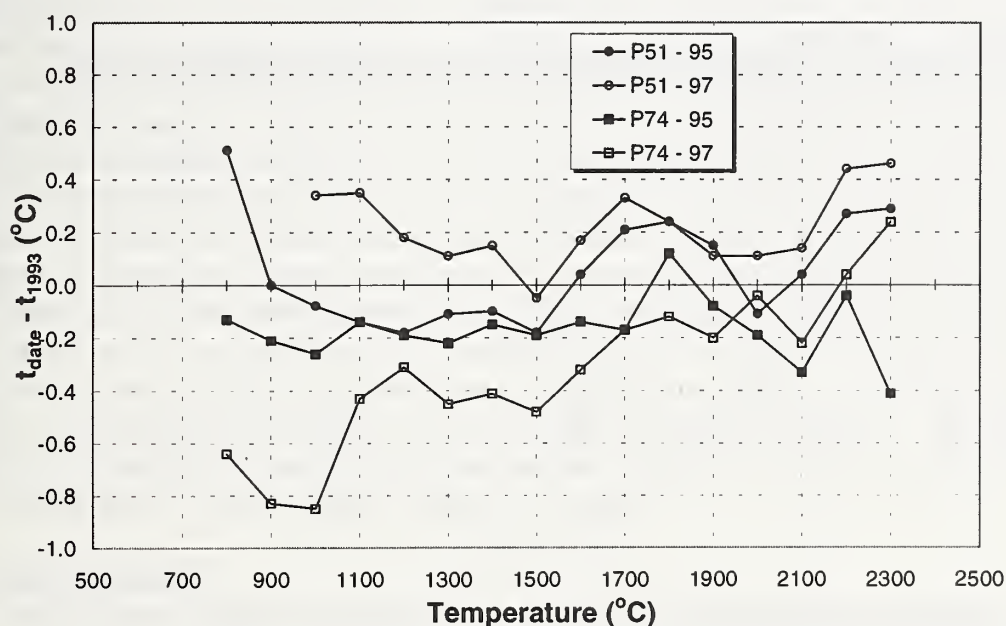


Figure 20. Calibration history of the ribbon filament lamp check standards.

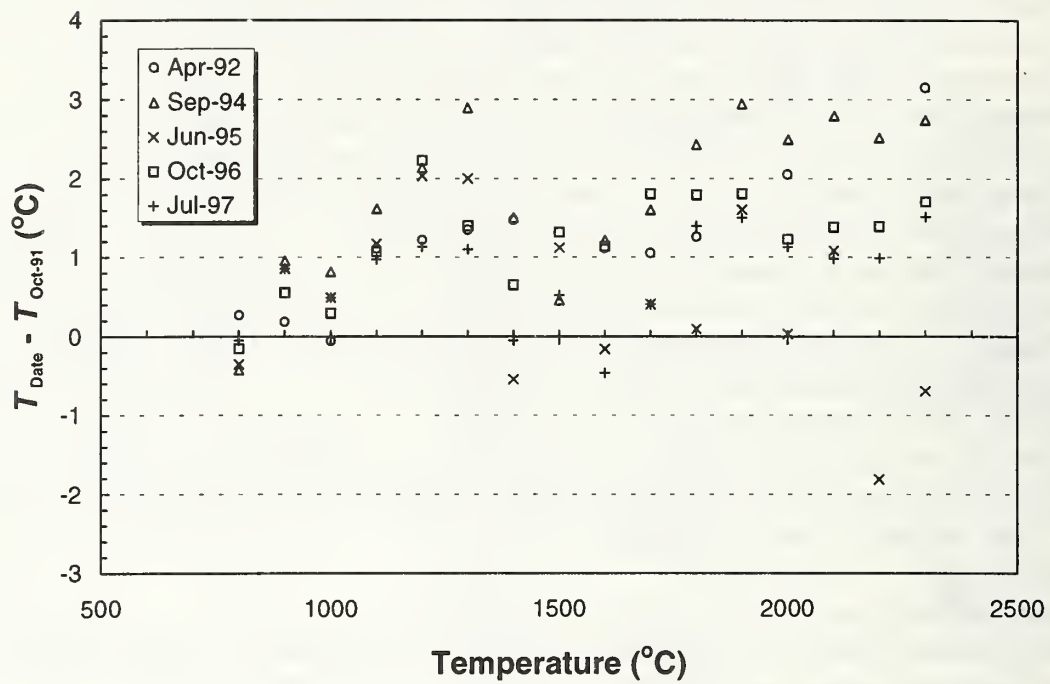


Figure 21. Calibration history of the radiation thermometer check standard (Minolta Cyclops 52).

4. Tungsten Ribbon Filament Lamp Calibrations

4.1 Lamp preparation

Tungsten ribbon filament lamps are calibrated by NIST for the temperature range 800 °C to 2300 °C in terms of the 1990 NIST Radiance Temperature Scale by spectral comparison to a standard temperature lamp calibrated at 655.3 nm. New lamps are inspected for filament continuity, envelope and window clarity, and filament flatness and centering. A digital multimeter is used to check the filament continuity. The envelope or window is checked for defects and striations in the optical path. The filament is checked to determine that it is flat, parallel to the envelope walls, and within 2 mm of being centered in the envelope or window. The lamps are cleaned with water and laboratory glassware detergent.

Test lamps (TL) are operated base down on direct current with the center contact (longer filament support for a bipost base lamp) at the positive potential. The test lamps are tested for air leaks by turning the lamp on, then off, and then on again after 15 min. If air has leaked into the lamp envelope because of bad or broken seals, then the lamp will fail when it is turned on. Lamps that pass the leak test are annealed to improve lamp stability. New gas lamps are annealed at a radiance temperature of 2350 °C at 655.3 nm for 2 h on direct current and vacuum lamps are annealed at 1450 °C for about 14 h. The gas lamps are aged for an additional 50 h at 2300 °C.

The lamp orientation chosen for the calibration of a new lamp minimizes the variation in lamp output while maintaining the optical axis of the measuring instrument approximately normal to the lamp filament. This orientation is determined with the lamp operating at approximately 1700 °C for a gas lamp, or 1200 °C for a vacuum lamp. An arrow is etched onto the rear surface of the lamp envelope to allow reproducible alignment of this orientation. The alignment is performed with the lamp operating base down, the filament vertical, and the optical axis of the PEP passing through the lamp envelope and intersecting the center of the filament at the height of the notch. The sides of the target area, area 0.6 mm wide by 0.8 mm high, are approximately parallel to the sides of the lamp filament. The center of the target area is located at the intersection of two orthogonal lines on the filament surface. One line bisects the filament lengthwise, and the other passes through the point centered at the mouth of the notch. The lamp is positioned so that the etched arrow on the lamp envelope is to the rear, as viewed from the PEP. The center of the target area is viewed along the horizontal optical axis of the PEP. A plumb line is used to make the notch side of the filament vertical. The image of the lamp filament is focused onto the field stop of the PEP to within ± 2.5 mm. The lamp is then successively rotated about the horizontal and vertical centerlines through the target area until the tip of the arrowhead is centered at the mouth of the notch.

4.2 Lamp calibration

The values of radiance temperature apply when the lamp has been aligned to a specified orientation while operating at a designated radiance temperature and after the lamp has reached thermal stability at each specified operating current. The WS and TL are aligned as described in Section 4.1 at one temperature. At all other temperatures, the TL is translated vertically and/or horizontally so that the target area viewed by the PEP is always centered on the lamp filament at the height of the notch. No additional rotational alignments are performed

The initial lamp current that corresponds to 1700 °C is selected from previous data for the lamp type. The TL is turned on and set to approximately 1700 °C and aligned after 30 min. The spectral radiance ratio of the TL to the WS, r_2 , is calculated, the WS signal, S_{WS2} , is measured, and the TL signal, S_{TL} , that corresponds to 1700 °C is calculated. The lamp current is adjusted until the measured signal is equal to S_{TL} . After waiting 10 min, the final alignment of the TL and the WS is performed and the TL is spectrally compared to the WS to determine its radiance temperature. The spectral radiance ratio of the TL to the WS is measured three times by alternately translating the WS and the TL to the PEP. If the percent standard deviation of the ratio is less than the control limit, then the next calibration point is measured; otherwise, the point is repeated. The control limits are equal to the relative standard uncertainty ($k = 1$) of the ratio and were determined from previous calibration data. The control limits are the expected uncertainties for the spectral radiance ratio measurement. Next, the TL is aligned and measured at 2300 °C after waiting 20 min. The additional calibration temperatures are aligned and measured in decreasing order after waiting 10 min. A typical calibration interval is 100 °C. The following temperatures are repeated at the previously determined lamp currents to estimate the lamps stability: 1700 °C, 2300 °C, 1900 °C, 1500 °C, 1100 °C, and 800 °C. The following data is stored in a computer file: nominal temperature, measured temperature, lamp current, and spectral radiance ratio. See the “Ribbon Filament Lamp Calibration Procedure” in Appendix D for a more details on lamp calibrations in the RTLC.

4.3 Lamp data analysis

The repeated calibration temperatures are averaged and the standard error of the mean is calculated and compared to the absolute control limits. These results are used to determine the expanded uncertainties listed in the calibration report.

Calibration data is reported at the nominal temperatures 800 °C, 900 °C, 1000 °C, ... 2300 °C. One method to determine the nominal temperatures is to measure it directly. Through an iterative process, the TL is repeatedly measured and the lamp current is adjusted until the nominal temperature is actually measured. However, the RTCL calculates the lamp current that corresponds to the nominal temperature from measurements of the measured temperature (T_m) near the nominal temperature (T_n).

For a small ΔT , where $\Delta T = T_n - T_m$, the small change in lamp current is approximated by the slope of the tangent at T and given by

$$\frac{dI(T)}{dT} = \lim_{\Delta T \rightarrow 0} \frac{\Delta I}{\Delta T} = \lim_{\Delta T \rightarrow 0} \frac{I(T + \Delta T) - I(T)}{\Delta T} . \quad (35)$$

This limit simplifies to a working equation of the form

$$I(T + \Delta T) = I(T) + \frac{dI}{dT} \cdot \Delta T , \quad (36)$$

where $I(T)$ is the measured lamp current I_m at T_m and $I(T + \Delta T)$ is the measured lamp current I_{corr} . The derivative dI/dT can be approximated by $\Delta I/\Delta T$, which can be modeled as

$$\frac{\Delta I}{\Delta T} = \sum_{j=0}^5 a_j T_{\text{mean}}^j . \quad (37)$$

The coefficients a_j are calculated from a least squares fit of two other parameters, T_{mean} and $\Delta I/\Delta T$. The first parameter is defined by

$$T_{\text{mean}} = \frac{T_{m,k} + T_{m,k+1}}{2} , \quad (38)$$

and the slope of the line between the calibration points P_k and P_{k+1} in figure 22 is given by

$$\frac{\Delta I}{\Delta T} = \frac{I_{m,k} - I_{m,k+1}}{T_{m,k} - T_{m,k+1}} . \quad (39)$$

Figure 22 shows the shape of the curve for T_m versus I_m . $\Delta I/\Delta T(T_{\text{mean}})$ is fit to a fifth degree polynomial giving the coefficients a_j . Figure 23 shows the shape of the curve with of T_{mean} and $\Delta I/\Delta T$ used to calculate the coefficients a_j . From eq (36), the corrected lamp current is

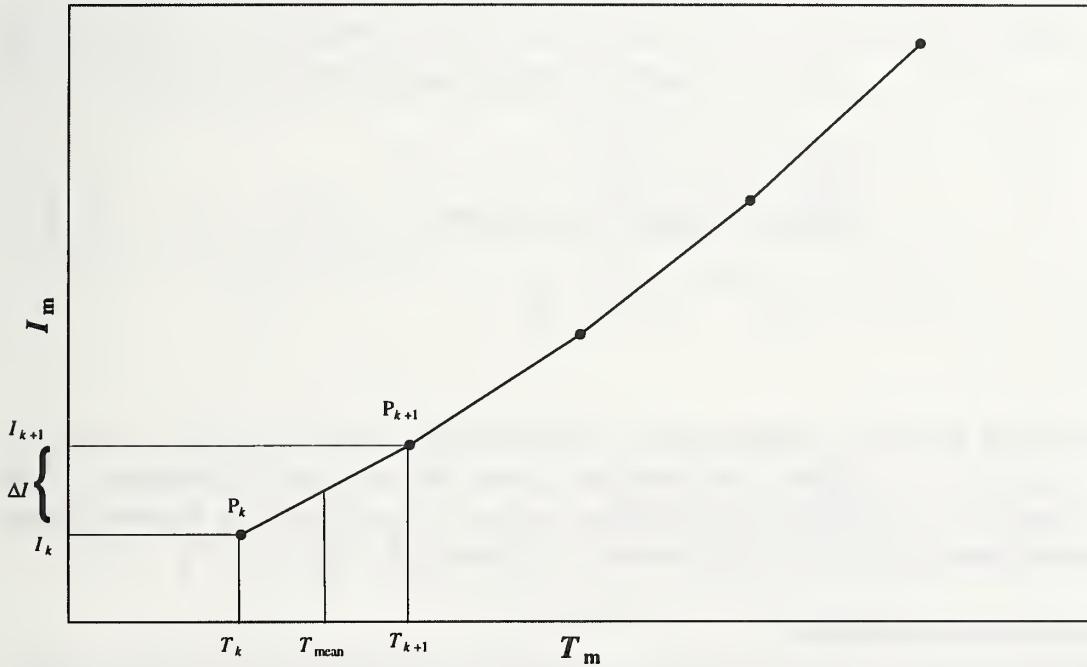


Figure 22. Calculation of I_{corr} from T_{mean} and $\Delta I/\Delta T$.

determined to be

$$I_{\text{corr}} = I_m + \frac{dI}{dT} (T_n - T_m) . \quad (40)$$

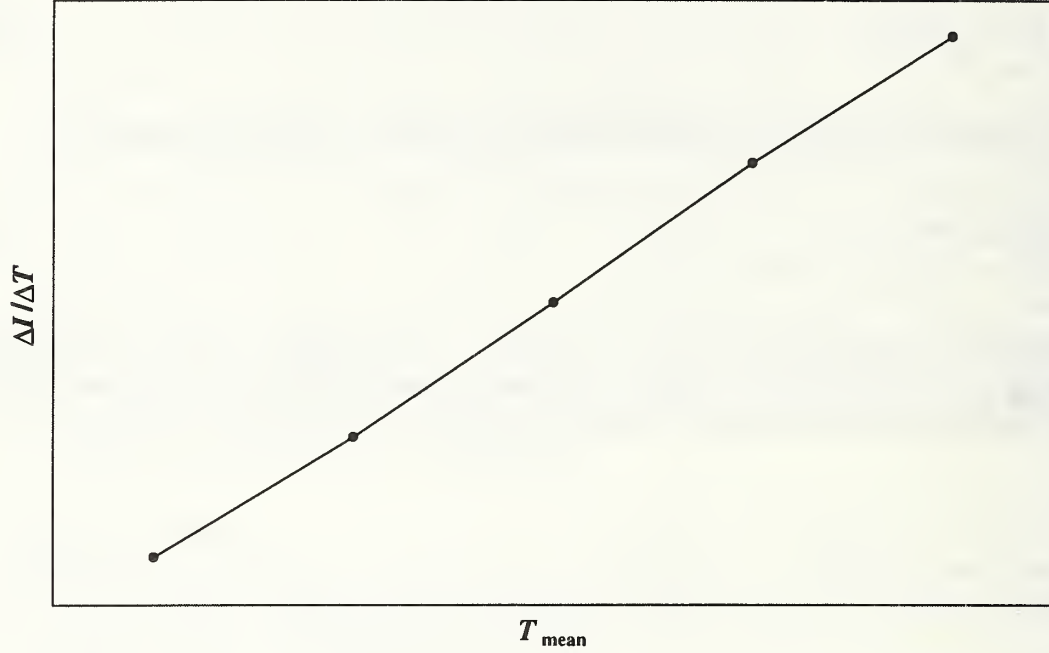


Figure 23. Graph of T_{mean} and $\Delta I/\Delta T$.

A least squares fit of the entire calibration data sets finds the best curve for the data set but not necessarily for each data point. This method ensures that the direction of the correction is right and has been verified by comparison to the iterative method. Sample data can be found in the “Ribbon Filament Lamp Data Reduction Procedure” in Appendix D.

4.4 Lamp calibration uncertainty

To calibrate test lamps (TL), the ratio (r_2) of the spectral radiance of the WS lamp to that of the TL,

$$r_2 = \frac{L_\lambda(T_{\text{TL}})}{L_\lambda(T_{\text{WS2}})} = \frac{S_{\text{TL}}}{S_{\text{WS2}}} , \quad (41)$$

is first measured. Taking into account the correction factors for the signals, the overall measurement equation for the 1990 NIST Scale of Radiance Temperature for calibration of the TL temperature is

$$T_{\text{TL}} = \frac{c_2}{n_\lambda \cdot \lambda \cdot \ln \left[1 + \frac{\varepsilon_{\lambda, \text{TL}} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{\text{WS2}} \cdot r_2} \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{WS}}}{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{TL}}} \right]}. \quad (42)$$

In a manner similar to eq (28), the uncertainty in the spectral radiance of the WS can be represented by the RSS of products of partial derivatives,

$$u_0(T_{\text{TL}}) = \left[\left(\frac{\partial T_{\text{TL}}}{\partial n_\lambda} \cdot u(n_\lambda) \right)^2 + \left(\frac{\partial T_{\text{TL}}}{\partial \lambda} \cdot u(\lambda) \right)^2 + \left(\frac{\partial T_{\text{TL}}}{\partial c_2} \cdot u(c_2) \right)^2 + \sum_{i=1}^{12} \left(\frac{\partial T_{\text{TL}}}{\partial x_i} \cdot u(x_i) \right)^2 \right]^{1/2}, \quad (43)$$

where x_i is one of the variables: ε_{TL} , c_{1L} , L_{WS} , r_2 , $C_{A, \text{WS}}$, $C_{L, \text{WS}}$, $C_{S, \text{WS}}$, G_{WS} , $C_{A, \text{TL}}$, $C_{L, \text{TL}}$, $C_{S, \text{TL}}$, or G_{TL} . The partial derivatives of T_{TL} with respect to n_λ , λ , and c_2 are

$$\frac{\partial T_{\text{TL}}}{\partial n_\lambda} = \frac{T_{\text{TL}}}{n_\lambda} \cdot \left[\frac{\varepsilon_{\lambda, \text{TL}} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{\text{WS2}}} \cdot \frac{2 \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{WS}}}{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{TL}}}}{r_2 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{TL}}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{TL}}}\right)} - 1 \right], \quad (44)$$

$$\frac{\partial T_{\text{TL}}}{\partial \lambda} = \frac{T_{\text{TL}}}{\lambda} \cdot \left[\frac{\varepsilon_{\lambda, \text{TL}} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{\text{WS2}}} \cdot \frac{5 \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{WS}}}{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{TL}}}}{r_2 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{TL}}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{\text{TL}}}\right)} - 1 \right], \text{ and} \quad (45)$$

$$\frac{\partial T_{\text{TL}}}{\partial c_2} = \frac{T_{\text{TL}}}{c_2}. \quad (46)$$

The partial derivatives of T_{TL} with respect to x are

$$\left| \frac{\partial T_{\text{TL}}}{\partial x} \right| = \frac{T_{\text{TL}}}{x} \cdot \frac{\varepsilon_{\lambda, \text{TL}} \cdot c_{1L}}{n_{\lambda}^2 \cdot \lambda^5 \cdot L_{\text{WS2}}} \cdot \frac{\frac{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{WS}}}{(C_A \cdot C_L \cdot C_S \cdot G)_{\text{TL}}}}{r_2 \cdot \frac{c_2}{n_{\lambda} \cdot \lambda \cdot T_{\text{TL}}} \cdot \exp\left(\frac{c_2}{n_{\lambda} \cdot \lambda \cdot T_{\text{TL}}}\right)}, \quad (47)$$

where x is one of the variables labeled x_i in eq (43).

Typical values of the variables in eq (42) are shown in table 9, and the uncertainties of eq (43) are displayed in table 10. The total expanded relative uncertainty $u(T_{\text{TL}})/T_{\text{TL}}$ in the TL temperature has contributions from $u_0(T_{\text{TL}})$ in eq (43), from the uncertainty $u(\text{TLC})$ in the TL current, and from the uncertainty $u(\text{AL})$ in alignment. These contributions are shown in table 10 and are summed by the RSS technique as follows,

$$\frac{u(T_{\text{TL}})}{T_{\text{TL}}} = \left[\left(\frac{u_0(T_{\text{TL}})}{T_{\text{TL}}} \right)^2 + \left(\frac{u(\text{TLC})}{T_{\text{TL}}} \right)^2 + \left(\frac{u(\text{AL})}{T_{\text{TL}}} \right)^2 \right]^{1/2}. \quad (48)$$

From tables 9 and 10, a typical relative expanded uncertainty in the TL temperature realization is 1.10 K/1973.52 K, or about 0.059 %. For other nominal TL temperatures besides 1700 °C, the user is referred to table 2.

Table 9. Typical values of TL variables and parameters

Variable	Symbol	Value
Refractive index	n_{λ}	1.00028
Wavelength	λ	655.3 nm
Second radiation constant	c_2	14387.69 $\mu\text{m} \cdot \text{K}$
Emissivity of TL	ε_{TL}	1
First radiation constant	c_{1L}	$1.191 \times 10^8 \text{ W} \cdot \mu\text{m}^4 \cdot \text{m}^{-2}$
WS spectral radiance	L_{WS}	$569.9 \text{ W} \cdot \text{m}^{-2} \cdot \mu\text{m}^{-1} \cdot \text{sr}^{-1}$
Ratio of TL signal to WS signal	r_2	2.551
WS amplifier calibration correction	$C_{A, \text{WS}}$	0.09986
WS linearity correction	$C_{L, \text{WS}}$	1
WS size of source correction	$C_{S, \text{WS}}$	1
WS amplifier gain	G_{WS}	$1 \times 10^9 \text{ V} \cdot \text{A}^{-1}$
TL amplifier calibration correction	$C_{A, \text{TL}}$	1
TL linearity correction	$C_{L, \text{TL}}$	1
TL size of source correction	$C_{S, \text{TL}}$	1
TL amplifier gain	G_{TL}	$1 \times 10^8 \text{ V} \cdot \text{A}^{-1}$
TL temperature	T_{TL}	1973.52 K

Table 10. Uncertainty budget for the TL temperature calibration

Uncertainty component	Symbol	Expanded Uncertainty ($k = 2$)	
		Type A	Type B
Refractive index	$u(n_\lambda)$		0.00002
Wavelength	$u(\lambda)$		0.2 nm
Second radiation constant	$u(c_2)$		0.24 $\mu\text{m}\cdot\text{K}$
Emissivity of TL	$u(\varepsilon_{\text{TL}})$		0.0002
First radiation constant	$u(c_{1\text{L}})$		440 $\text{W}\cdot\mu\text{m}^4\cdot\text{m}^{-2}$
WS spectral radiance	$u_0(L_{\text{WS}})$		2.998 $\text{W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
Ratio of TL signal to WS signal	$u(r_2)$	0.0051	
WS amplifier calibration correction	$u(C_{\text{A,WS}})$	0.00001	
WS linearity correction	$u(C_{\text{L,WS}})$	0.001	
WS size of source correction	$u(C_{\text{S,WS}})$	0.0002	
WS amplifier gain	$u(G_{\text{WS}})$		0 $\text{V}\cdot\text{A}^{-1}$
TL amplifier calibration correction	$u(C_{\text{A,TL}})$	0.0001	
TL linearity correction	$u(C_{\text{L,TL}})$	0.001	
TL size of source correction	$u(C_{\text{S,TL}})$	0.0002	
TL amplifier gain	$u(G_{\text{TL}})$		0 $\text{V}\cdot\text{A}^{-1}$
TL temperature	$u_0(T_{\text{TL}})$		1.10 K
TL current	$u(\text{TLC})$	0.09 K	
Alignment	$u(\text{AL})$	0.30 K	
TL temperature calibration	$u(T_{\text{TL}})$		1.10 K

5. Disappearing filament optical pyrometer calibrations

5.1 Optical pyrometer preparation

Disappearing filament optical pyrometers (DFP) are calibrated from 800 °C to 4200 °C by comparison to the VTBB. When the test unit arrives, the batteries are checked, the objective lens and eyepiece are cleaned, and the eyepiece is focused onto the DFP lamp filament. The lamp is checked to ensure that it is operating and that matching is possible. If an operating manual is not on file, then one is requested from the customer because test pyrometers are operated according to the manufacturer's instructions unless other instructions are given. The pyrometers and associated electronics are initialized according to manufacturer's specifications.

5.2 Optical pyrometer calibration

The DFP is aligned onto the geometric center of the VTBB opening. The DFP is measured at an appropriate distance from the front surface of the objective lens holder to the blackbody cavity partition and is focussed on the center of the partition. After the VTBB is turned on, the pyrometer centering is checked and adjusted as necessary.

For the 800 °C to 2700 °C temperature range, the PEP is used to spectrally compare the VTBB to the WS to determine the radiance temperature of the blackbody at each calibration temperature. The VTBB and the WS are compared three times and the mean is used. Once the temperature of the blackbody has been determined, the test DFP lamp filament is matched in brightness to that of the blackbody, and the DFP lamp current is measured using a standard resistor and a digital voltmeter. The matching is performed from bright to disappearance and from dark to disappearance twice and the mean of the four readings is used. The indicated temperature is recorded for pyrometers with either a temperature scale or a digital panel indicator. The output voltage is measured for pyrometers with internal current measuring resistors. The DFP measurements are not bracketed by VTBB measurements because the VTBB has been shown to be stable to within ± 0.1 °C during the time required to complete the measurements.

For temperatures above 2700 °C, the method used to extrapolate the calibration to 4200 °C involves measuring temperatures below 2700 °C on two or more ranges. For a three-range DFP with a low range (L) from 1100 °C to 1700 °C, a medium range (M) from 1700 °C to 2800 °C, and a high range (H) from 2900 °C to 4200 °C, the following measurements are performed. The DFP lamp currents are measured for VTBB temperatures from 1100 °C to 1700 °C on the L range, from 1600 °C to 2700 °C on the M range, and from 1700 °C to 2700 °C on the H range.

5.3 Optical pyrometer data analysis

The various data reduction methods used are determined by how the DFP indicates temperature. In all cases, the data is sorted by range and averaged and the standard deviation is calculated. For DFPs with ranges below 2700 °C, the data is analyzed and the blackbody temperature versus the scale or panel reading, output current, or output voltage is reported. The pyrometer correction, blackbody temperature minus DFP reading, that is added to the pyrometer reading is reported for scale and panel indicator DFPs.

For the three-range DFP case when temperatures are determined above 2700 °C, more calculations are required before the data is reported. The current measurement averages are fit to

the curve $I(T) = a_0 + a_1 \cdot T + a_2 \cdot T^2 + a_3 \cdot T^3$ calculating a least squares approximation for each DFP range. For the L range the blackbody temperatures are calculated at the nominal temperatures 1100 °C, 1200 °C, ..., 1700 °C and reported versus the pyrometer lamp currents. For the M range, the blackbody temperatures are calculated at the nominal temperatures 1700 °C, 1800 °C, ..., 2700 °C and reported versus the pyrometer lamp currents. The 2800 °C calibration point is determined in two steps. The first step is to calculate A_M , the “A” value [22] for the M range, which is defined by the relationship,

$$A = \frac{1}{T_{APP}} - \frac{1}{T_{BB}} = \frac{\lambda}{c_2} \cdot \ln \tau, \quad (44)$$

where T_{APP} is the apparent temperature [K], T_{BB} is the blackbody temperature [K], τ is the filter transmittance, and A is the “A” value [K⁻¹]. From eq (44), it becomes apparent that $\ln \tau$ must be approximately equal to $1/\lambda$ in order for A to be constant. Typically, the value of A is on the order of $173 \times 10^{-6} \text{ K}^{-1}$. T_{BB} is determined during the M range measurements. T_{APP} is the temperature on the L range that corresponds to the current measured on the M range which is determined through iteration of the L range function $I(t)$. The second step is to use the “A” value to calculate the pyrometer lamp current that corresponds to the extrapolated blackbody temperature at 2800 °C. Use eq (44) to calculate T_{APP2} when $T_{BB2} = 3073.15 \text{ K}$ and $A = A_M$. Next calculate the pyrometer lamp current using the L range function $I(t)$.

This procedure is repeated to determine the extrapolated blackbody temperatures for the H range. The H range “A-value”, A_H , is determined where T_{APP} is the temperature on the L range that corresponds to the current measured on the H range and T_{BB} is the blackbody temperature during the H range measurements. A_H is then used to calculate the apparent temperature T_{APP2} for each extrapolated blackbody temperature from 2900 °C to 4200 °C. The pyrometer lamp currents are calculated from the L range function $I(t)$. The extrapolated blackbody temperatures for the M and H ranges are reported versus the pyrometer lamp currents.

5.4 Optical pyrometer calibration uncertainty

To calibrate disappearing filament optical pyrometers (DFP), the ratio (r_3) of the spectral radiance of the WS lamp to that of the variable-temperature blackbody (BB),

$$r_3 = \frac{L_\lambda(T_{BB})}{L_\lambda(T_{WS3})} = \frac{S_{BB}}{S_{WS3}}, \quad (45)$$

is first measured to determine the radiance temperature of the BB, $T(\lambda, L_{\lambda, BB})$, at each disappearing filament optical pyrometer temperature by the equation,

$$T_{BB} = \frac{c_2}{n_\lambda \cdot \lambda \cdot \ln \left[1 + \frac{\epsilon_{\lambda, BB} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS3} \cdot r_3} \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}} \right]}. \quad (46)$$

The uncertainty in the spectral radiance of the BB can be calculated by the RSS of products of partial derivatives of eq (46) with their respective uncertainties as follows,

$$u_0(T_{BB}) = \left[\left(\frac{\partial T_{BB}}{\partial n_\lambda} \cdot u(n_\lambda) \right)^2 + \left(\frac{\partial T_{BB}}{\partial \lambda} \cdot u(\lambda) \right)^2 + \left(\frac{\partial T_{BB}}{\partial c_2} \cdot u(c_2) \right)^2 + \sum_{i=1}^{12} \left(\frac{\partial T_{BB}}{\partial x_i} \cdot u(x_i) \right)^2 \right]^{1/2}, \quad (47)$$

where x_i is one of the following variables: ε_{BB} , c_{1L} , r_3 , $C_{A,WS}$, $C_{L,WS}$, $C_{S,WS}$, G_{WS} , $C_{A,BB}$, $C_{L,BB}$, $C_{S,BB}$, or G_{BB} . The partial derivatives of n_λ , λ , and c_2 can be derived from eq (46) as

$$\frac{\partial T_{BB}}{\partial n_\lambda} = \frac{T_{BB}}{n_\lambda} \cdot \left[\frac{\varepsilon_{\lambda,BB} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS3}} \cdot \frac{2 \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}}}{r_3 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}}\right)} - 1 \right], \quad (48)$$

$$\frac{\partial T_{BB}}{\partial \lambda} = \frac{T_{BB}}{\lambda} \cdot \left[\frac{\varepsilon_{\lambda,BB} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS3}} \cdot \frac{5 \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}}}{r_3 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}}\right)} - 1 \right], \text{ and} \quad (49)$$

$$\frac{\partial T_{BB}}{\partial c_2} = \frac{T_{BB}}{c_2}. \quad (50)$$

The expression,

$$\left| \frac{\partial T_{BB}}{\partial x} \right| = \frac{T_{BB}}{x} \cdot \frac{\varepsilon_{\lambda,BB} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS3}} \cdot \frac{\frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}}}{r_3 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}}\right)}, \quad (51)$$

represents the partial derivatives of T_{BB} with respect to x , where x is one of the variables, ϵ_{BB} , c_{1L} , r_3 , $C_{A,WS}$, $C_{L,WS}$, $C_{S,WS}$, G_{WS} , $C_{A,BB}$, $C_{L,BB}$, $C_{S,BB}$, or G_{BB} .

A summary of the DFP variable values in eq (46) in table 11. Uncertainty values from eq (47) are given in table 12. Table 12 also shows the standard uncertainties leading up to the total standard uncertainty in the DFP temperature,

$$\frac{u(T_{BB})}{T_{BB}} = \left[\left(\frac{u_0(T_{BB})}{T_{BB}} \right)^2 + \left(\frac{u(DMM)}{T_{BB}} \right)^2 \right]^{1/2}. \quad (52)$$

The uncertainty in the pyrometer temperature is related to the uncertainty in the variable-temperature blackbody temperature by the following relationship,

$$\frac{u(T_{DFP})}{T_{DFP}} = \left[\frac{u^2(T_{BB})}{T_{BB}^2} + \left(\frac{u(BU)}{T_{DFP}} \right)^2 + \left(\frac{u(PR)}{T_{DFP}} \right)^2 \right]^{1/2}, \quad (53)$$

where $u(BU)$ is the blackbody uniformity uncertainty and $u(PR)$ is the pyrometer reading uncertainty. The relative expanded uncertainty in the DFP temperature from tables 11 and 12 is 7.57 K/2579.07 K, or 0.29 %. The uncertainty values for other DFP temperatures are listed in table 2.

Table 11. Typical values of DFP variables and parameters

Variable	Symbol	Value
Refractive index	n_λ	1.00028
Wavelength	λ	655.3 nm
Second radiation constant	c_2	14387.69 $\mu\text{m}\cdot\text{K}$
Emissivity of BB	ϵ_{BB}	0.99
First radiation constant	c_{1L}	$1.191 \times 10^8 \text{ W}\cdot\mu\text{m}^4\cdot\text{m}^{-2}$
WS spectral radiance	L_{WS}	$569.9 \text{ W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
Ratio of BB signal to WS signal	r_3	3.445
WS amplifier calibration correction	$C_{A,WS}$	0.09986
WS linearity correction	$C_{L,WS}$	1
WS size of source correction	$C_{S,WS}$	1
WS amplifier gain	G_{WS}	$1 \times 10^9 \text{ V}\cdot\text{A}^{-1}$
BB amplifier calibration correction	$C_{A,BB}$	9.978
BB linearity correction	$C_{L,BB}$	0.9997
BB size of source correction	$C_{S,BB}$	0.9987
BB amplifier gain	G_{BB}	$1 \times 10^7 \text{ V}\cdot\text{A}^{-1}$
BB temperature	T_{BB}	2579.07 K

Table 12. Uncertainty budget for the DFP pyrometer temperature calibration

Uncertainty factor	Symbol	Expanded Uncertainty ($k = 2$)	
		Type A	Type B
Refractive index	$u(n_\lambda)$		0.00002
Wavelength	$u(\lambda)$		0.2 nm
Second radiation constant	$u(c_2)$		0.24 $\mu\text{m}\cdot\text{K}$
Emissivity of BB	$u(\epsilon_{\text{BB}})$		0.0002
First radiation constant	$u(c_{1\text{L}})$		440 $\text{W}\cdot\mu\text{m}^4\cdot\text{m}^{-2}$
WS spectral radiance	$u(L_{\text{WS}})$		2.998 $\text{W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
Ratio of BB signal to WS signal	$u(r_3)$	0.00682	
WS amplifier calibration correction	$u(C_{\text{A,WS}})$	0.00001	
WS linearity correction	$u(C_{\text{L,WS}})$	0.001	
WS size of source correction	$u(C_{\text{S,WS}})$	0.0002	
WS amplifier gain	$u(G_{\text{WS}})$		0 $\text{V}\cdot\text{A}^{-1}$
BB amplifier calibration correction	$u(C_{\text{A,BB}})$	0.001	
BB linearity correction	$u(C_{\text{L,BB}})$	0.001	
BB size of source correction	$u(C_{\text{S,BB}})$	0.0002	
BB amplifier gain	$u(G_{\text{BB}})$		0 $\text{V}\cdot\text{A}^{-1}$
Digital voltmeter	$u(\text{DMM})$	0 K	
BB temperature	$u(T_{\text{BB}})$		2.34 K
Blackbody uniformity	$u(\text{BU})$	0.20 K	
DFP reading	$u(\text{PR})$	7.20 K	
DFP temperature calibration	$u(T_{\text{DFP}})$		7.57 K

6. Infrared thermometer calibrations

6.1 Radiation pyrometer preparation

Radiation thermometers (RT) are measured for the temperature range from 800 °C to 2700 °C. The RTs are operated according to the manufacturer's operating instructions unless special instructions are given. First, the lens, eyepiece, and filters are cleaned; then the emissivity control is set to unity. The RT is aligned so that its optical axis, when viewing the VTBB, coincides with the geometrical center of the blackbody.

6.2 Radiation pyrometer calibration

The VTBB is turned on and set to the highest temperature measured. After waiting 30 min, the PEP is used to determine the radiance temperature of the VTBB by spectrally comparing it to the WS lamp at each calibration temperature from 2700 °C to 800 °C in 100 °C steps. The measurements are repeated at 2300 °C, 1900 °C, 1500 °C, 1100 °C, and 800 °C to determine the stability of the RT. The VTBB is compared to the WS three times, and the average of the measurements is used. After the temperature of the VTBB is determined, three RT measurements of the VTBB are made, and the output current, output voltage, or indicated temperature is recorded. The spectral radiance temperatures of the VTBB used in the calibration are determined at wavelengths of 655.3 nm, 900 nm, and 1000 nm. For 100 °C temperature changes, the VTBB is measured after waiting 5 min for the VTBB to stabilize.

6.3 Radiation pyrometer data analysis

For RTs with output currents or voltages, the nominal temperature, the VTBB temperature, and the output current or voltage are reported in the test report. For RTs with temperature indicators, the nominal temperature, the VTBB temperature, the temperature indicated by the RT, and the RT correction are reported in the test report. The RT correction is equal to the VTBB temperature minus the RT indicated temperature.

6.4 Radiation pyrometer calibration uncertainty

To calibrate radiation thermometers (RT), the ratio (r_4) of the spectral radiance of the WS lamp to that of the variable-temperature blackbody (BB),

$$r_4 = \frac{L_\lambda(T_{RT})}{L_\lambda(T_{WS4})} = \frac{S_{RT}}{S_{WS4}}, \quad (54)$$

is first measured to determine the radiance temperature of the BB, $T(\lambda, L_{\lambda, BB})$, at each test temperature by the equation,

$$T_{BB} = \frac{c_2}{n_\lambda \cdot \lambda \cdot \ln \left[1 + \frac{\epsilon_{\lambda, RT} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS4} \cdot r_4} \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}} \right]}. \quad (55)$$

The uncertainty in the spectral radiance of the BB can be calculated by the RSS of the products of the partial derivatives in eqs (57) through (60) with their respective uncertainties,

$$u_0(T_{BB}) = \left[\left(\frac{\partial T_{BB}}{\partial n_\lambda} \cdot u(n_\lambda) \right)^2 + \left(\frac{\partial T_{BB}}{\partial \lambda} \cdot u(\lambda) \right)^2 + \left(\frac{\partial T_{BB}}{\partial c_2} \cdot u(c_2) \right)^2 + \sum_{i=1}^{12} \left(\frac{\partial T_{BB}}{\partial x_i} \cdot u(x_i) \right)^2 \right]^{1/2}, \quad (56)$$

where x_i is one of the following variables: ϵ_{RT} , c_{1L} , r_4 , $C_{A,WS}$, $C_{L,WS}$, $C_{S,WS}$, G_{WS} , $C_{A,BB}$, $C_{L,BB}$, $C_{S,BB}$, or G_{BB} . From eq (55), the partial derivatives of n_λ , λ , and c_2 are

$$\frac{\partial T_{BB}}{\partial n_\lambda} = \frac{T_{BB}}{n_\lambda} \cdot \left[\frac{\epsilon_{\lambda,RT} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS4}} \cdot \frac{2 \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}}}{r_4 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}}\right)} - 1 \right], \quad (57)$$

$$\frac{\partial T_{BB}}{\partial \lambda} = \frac{T_{BB}}{\lambda} \cdot \left[\frac{\epsilon_{\lambda,RT} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS4}} \cdot \frac{5 \cdot \frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}}}{r_4 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}}\right)} - 1 \right], \text{ and} \quad (58)$$

$$\frac{\partial T_{BB}}{\partial c_2} = \frac{T_{BB}}{c_2}. \quad (59)$$

The expression,

$$\left| \frac{\partial T_{BB}}{\partial x} \right| = \frac{T_{BB}}{x} \cdot \frac{\epsilon_{\lambda,RT} \cdot c_{1L}}{n_\lambda^2 \cdot \lambda^5 \cdot L_{WS4}} \cdot \frac{\frac{(C_A \cdot C_L \cdot C_S \cdot G)_{WS}}{(C_A \cdot C_L \cdot C_S \cdot G)_{BB}}}{r_4 \cdot \frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}} \cdot \exp\left(\frac{c_2}{n_\lambda \cdot \lambda \cdot T_{BB}}\right)}, \quad (60)$$

represents the partial derivatives of T_{BB} with respect to x , where x is one of the variables, ϵ_{RT} , c_{1L} ,

r_4 , $C_{A,WS}$, $C_{L,WS}$, $C_{S,WS}$, G_{WS} , $C_{A,BB}$, $C_{L,BB}$, $C_{S,BB}$, or G_{BB} .

Table 13 shows the typical values of the variables in eq (55). The uncertainties of eq (56) are displayed in table 14. Finally, the total expanded uncertainty shown in table 14 is

$$\frac{u(T_{BB})}{T_{BB}} = \left[\left(\frac{u_0(T_{BB})}{T_{BB}} \right)^2 + \left(\frac{u(\text{DMM})}{T_{BB}} \right)^2 \right]^{1/2} \quad (61)$$

The uncertainty in the RT temperature is related to the uncertainty in the variable-temperature blackbody temperature by the following relationship,

$$\frac{u(T_{RT})}{T_{RT}} = \left[\frac{u^2(T_{BB})}{T_{BB}^2} + \left(\frac{u(\text{BU})}{T_{RT}} \right)^2 + \left(\frac{u(\text{RTR})}{T_{RT}} \right)^2 \right]^{1/2} \quad (62)$$

From tables 13 and 14, the relative expanded uncertainty in the T_{RT} is 3.05 K/2579.07 K, or 0.12 %. For other RT temperatures, table 2 can be used to approximate the uncertainties.

Table 13. Typical values of RT variables and parameters

Variable	Symbol	Value
Refractive index	n_λ	1.00028
Wavelength	λ	655.3 nm
Second radiation constant	c_2	14387.69 $\mu\text{m}\cdot\text{K}$
Emissivity of BB	ϵ_{BB}	0.99
First radiation constant	c_{1L}	$1.191 \times 10^8 \text{ W}\cdot\mu\text{m}^4\cdot\text{m}^{-2}$
WS spectral radiance	L_{WS}	$569.9 \text{ W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
Ratio of BB signal to WS signal	r_4	3.445
WS amplifier calibration correction	$C_{A,WS}$	0.09986
WS linearity correction	$C_{L,WS}$	1
WS size of source correction	$C_{S,WS}$	1
WS amplifier gain	G_{WS}	$1 \times 10^9 \text{ V}\cdot\text{A}^{-1}$
BB amplifier calibration correction	$C_{A,BB}$	9.978
BB linearity correction	$C_{L,BB}$	0.9997
BB size of source correction	$C_{S,BB}$	0.9987
BB amplifier gain	G_{BB}	$1 \times 10^7 \text{ V}\cdot\text{A}^{-1}$
BB temperature	T_{BB}	2579.07 K

Table 14. Uncertainty budget for the RT temperature calibration

Uncertainty component	Symbol	Expanded Uncertainty ($k = 2$)	
		Type A	Type B
Refractive index	$u(n_\lambda)$		0.00002
Wavelength	$u(\lambda)$		0.2 nm
Second radiation constant	$u(c_2)$		0.24 $\mu\text{m}\cdot\text{K}$
Emissivity of BB	$u(\varepsilon_{\text{BB}})$		0.0002
First radiation constant	$u(c_{1\text{L}})$		440 $\text{W}\cdot\mu\text{m}^4\cdot\text{m}^{-2}$
WS spectral radiance	$u(L_{\text{WS}})$		2.998 $\text{W}\cdot\text{m}^{-2}\cdot\mu\text{m}^{-1}\cdot\text{sr}^{-1}$
Ratio of BB signal to WS signal	$u(r_4)$	0.00579	
WS amplifier calibration correction	$u(C_{\text{A,WS}})$	0.00001	
WS linearity correction	$u(C_{\text{L,WS}})$	0.001	
WS size of source correction	$u(C_{\text{S,WS}})$	0.0002	
WS amplifier gain	$u(G_{\text{WS}})$		0 $\text{V}\cdot\text{A}^{-1}$
BB amplifier calibration correction	$u(C_{\text{A,BB}})$	0.001	
BB linearity correction	$u(C_{\text{L,BB}})$	0.001	
BB size of source correction	$u(C_{\text{S,BB}})$	0.0002	
BB amplifier gain	$u(G_{\text{BB}})$		0 $\text{V}\cdot\text{A}^{-1}$
Digital voltmeter	$u(\text{DMM})$		0 K
BB temperature	$u(T_{\text{BB}})$		2.29 K
Blackbody uniformity	$u(\text{BU})$	0.20 K	
RT reading	$u(\text{RTR})$	2.00 K	
RT temperature calibration	$u(T_{\text{RT}})$		3.05 K

7. Modifications to the NIST photoelectric pyrometer

The radiance temperature calibrations described in this document are performed using the NIST PEP, which is a single wavelength instrument. The spectral radiance temperatures at 900 nm and 1000 nm are calculated from measurements made at 655.3 nm. Modifications to the PEP are underway to allow direct measurements at 900 nm and 1000 nm. The modified pyrometer also uses refractive optics, interference filters, and a photomultiplier tube. The modifications will include adding a filtered silicon detector and amplifier system, and a laser to simplify source alignments. A schematic of the modified PEP is shown in figure 24. The methods described in the report to realize the radiance temperature scale at 655.3 nm will be employed to realize the scale in the infrared.

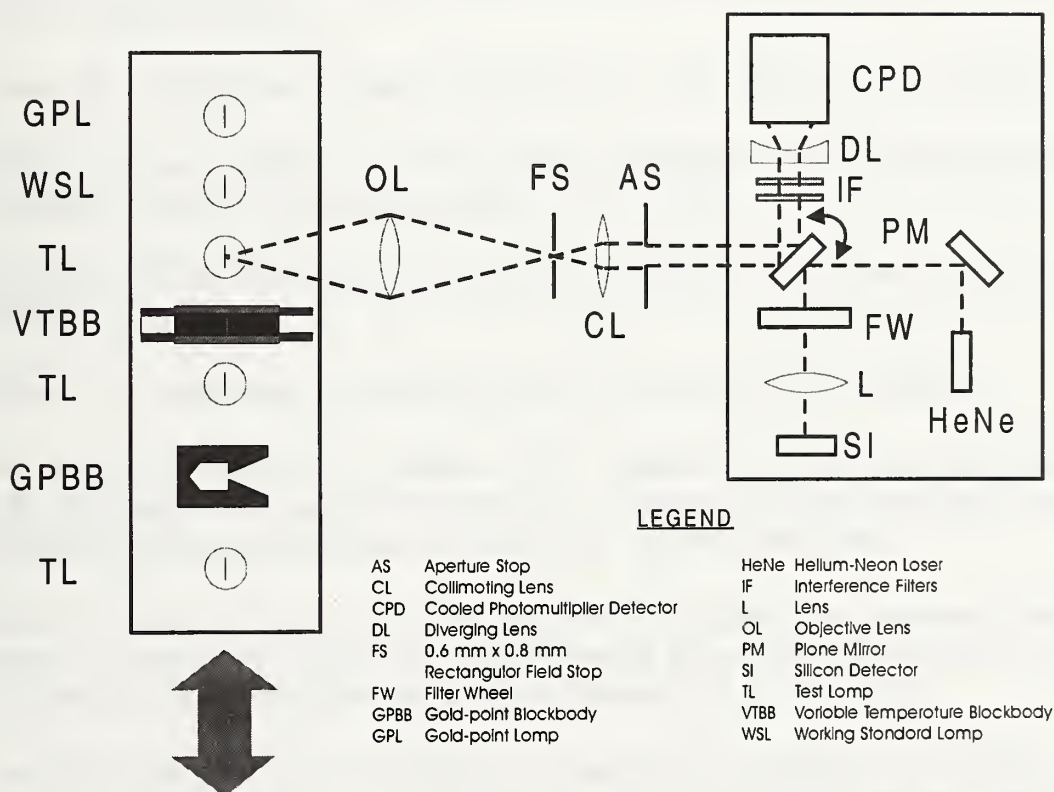


Figure 24. Modified NIST photoelectric pyrometer.

Acknowledgments

The authors would like to thank Robert D. Saunders for his technical assistance in reading the manuscript and his constant encouragement. They also like to thank Frank Lovas, Dean Ripple, and Chris Kuyatt for useful suggestions and comments.

References

- [1] Mielenz, K.D., Saunders, R.D., Parr, A.C., and Hsia, J.J., "The 1990 NIST Scales of Thermal Radiometry," *J. Res. NIST*, **95**, 621 (1990).
- [2] Preston-Thomas, H., "The International Temperature Scale of 1990 (ITS-90)," *Metrologia*, **27**, 3-10, 107 (1990).
- [3] Taylor, B.N., and Kuyatt, C.E., "Guidelines for Evaluating and Expressing the Uncertainty of NIST Measurement Results," *NIST TN 1297* (2nd ed., 1994).
- [4] Simmons, J.D. Ed., "NIST Calibration Services Users Guide," *NIST SP 250*, 69 (1991).
- [5] Kostkowski, H.J., and Nicodemus, F.E., "An Introduction to the Measurement Equation, Chapter 5 of Self-Study Manual on Optical Radiation Measurements: Part I – Concepts," *NBS TN 910-2*, 58 (1978).
- [6] Planck, M., "Distribution of Energy in the Spectrum," *Ann. Phys.*, **4**, 3, 553 (1901).
- [7] Cohen, E.R. and Taylor, B.N., "The 1986 Adjustment of the Fundamental Physical Constants," *Review of Modern Physics*, **59**, 4, 1121 (1987).
- [8] Lide, D.R., "CRC Handbook of Chemistry and Physics," CRC Press, 67th ed., E-373 (1986-1987).
- [9] Quinn, T.J., "Temperature," Academic Press Inc., San Diego, CA, 380 (1990).
- [10] Saunders, R.D., Gibson, C.E., Mielenz, K.D., Sapritsky, V.I., Sudarev, K.A., Khlevnoy, B.B., Mekhontsev, S.N, and Harchenko, G.D. , "Results of a NIST/VNIIOFI Comparison of Spectral-Radiance Measurements," *Metrologia*, **32**, 449 (1995/1996).
- [11] Quinn, T.J. and Lee, R.D., "Vacuum Tungsten Strip Lamps with Improved Stability as Radiance Temperature Standards," *Temperature, Its Measurement and Control in Science and Industry*, Instrument Society of America, Pittsburgh, **4**, 1, 395 (1972).
- [12] Lee, R.D., Kostkowski, H.J., Quinn, T.J., Chandler, P.R., Jones, T.N., Tapping, J., and Kunz, H., "Intercomparison of the IPTS68 Above 1064°C by Four National Laboratories," *Temperature, Its Measurement and Control in Science and Industry*, Instrument Society of America, Pittsburgh, **4**, 1, 377 (1972).
- [13] Eckerle, K.L., Hsia, J.J., Mielenz, K.D., and Weidner, V.R., "Regular Spectral Transmittance," *NBS SP 250-6* (1987).

- [14] Saunders, R.D. and Shumaker, J.B., "Automated Radiometric Linearity Tester," *Appl. Opt.*, **23**, 20, 3504, (1984).
- [15] Larason, T.C., Bruce, S.S., and Cromer, C.L., "The NIST High Accuracy Scale for Absolute Spectral Response from 406 nm to 920 nm," *J. Res. NIST*, **101**, 2, 133 (1996).
- [16] Mielenz, K.D., Saunders, R.D., and Shumaker, J. B., "Spectroradiometric Determination of the Freezing Temperature of Gold," *J. Res. NIST*, **95**, 49 (1990).
- [17] "American National Standard for Calibration – Calibration Laboratories and Measuring and Test Equipment – General Requirements," ANSI/NCSL Z540-1-1994.
- [18] "ISO Guide 25 General requirements for the competence of calibration and testing laboratories," 3rd edition (1990).
- [19] Bruce, S.S. and Larason, T.C., "Developing Quality System Documentation Based on ANSI/NCSL Z540-1-1994 : The Optical Technology Division's Effort," *NISTIR 5866* (1996).
- [20] Johnson, B.C., Gibson, C., Machin, G., and Rusby, R.L., "Intercomparison of the ITS-90 radiance temperature scales of the national physical laboratory and NIST," *J. Res. NIST*, **99**, 731 (1994).
- [21] Johnson, B.C., Sakuma, F., Sakate, H., Gibson, C., Machin, G., Ricolfi, T., Battuello, M., Fischer, J., and Jung, H.J., "International comparison of radiation temperature scales among five national metrological laboratories using a transfer standard radiation thermometer," *Metrolgia*, **33**, 241 (1996).
- [22] Kostkowski, H.J., and Lee, R.D. "Theory and Methods of Optical Pyrometry," *NBS Monograph 41*, (1962).

Appendix A: Ribbon filament lamp sample calibration report

In reply refer to: 844/123456-95

Radiometric Systems Inc.
Attn: John Doe
123 Calibration Court
Measurement City, MD 00000-0000

Subject: Ribbon Filament Lamp Report of Calibration

Order No.: AB1234 dated January 1, 1995

Dear Mr. Doe:

Enclosed are results of the calibration you requested for Test Number 35050C Radiance Temperature Standard, Ribbon Filament Lamp. Please refer to the above file number in any further communication concerning this calibration.

Sincerely,

Albert C. Parr, Chief
Optical Technology Division
Physics Laboratory

Enclosures:

One Report of Calibration

ACP/jah

NIST Cost Center 8443600

REPORT OF CALIBRATION

35050C Radiance Temperature Standard,
Ribbon Filament Lamp

for

Model # General Electric 30A/T24/6, Serial # P000

Submitted by:

John Doe
Radiometric Systems Inc.
123 Calibration Court
Measurement City, MD 00000-0000

(See your Purchase Order No. AB1234 dated January 1, 1995)

1. Description of Calibration Item

A gas-filled tungsten ribbon filament lamp with a mogul screw base was calibrated by the National Institute of Standards and Technology (NIST) as a standard of radiance temperature from 1000 °C to 2200 °C. The serial number is located on the rear of the lamp envelope opposite the side viewed by the pyrometer.

2. Description of Calibration

The lamp was calibrated in the NIST Radiance Temperature Calibration Facility using the equipment and procedures described in Ref. [1]. The NIST Photoelectric Pyrometer, which has a mean effective wavelength of 655.3 nm, was used to spectrally compare the test lamp to working standard SL20 to determine its radiance temperature at each test lamp current. The NIST Photoelectric Pyrometer is described in the Ref. [1].

The orientation and operation of the test lamp is described in the enclosed NIST SP250-43 entitled, "Radiance Temperature Calibrations." The calibration was performed with the lamp operating on direct current and with the center contact of the lamp base at positive potential. The lamp was aligned base down while operating at a radiance temperature of 1700 °C.

Laboratory Environment:

Temperature: 24 °C \pm 1 °C

Relative Humidity: 25% \pm 5%

Calibration Date: January 2, 1995

NIST Test No.: 844/123456-95

Page 1 of 5

REPORT OF CALIBRATION

35050C Radiance Temperature Standard, Ribbon Filament Lamp
Radiometric Systems Inc.

Model #: GE 30A/T24/6
Serial #: P000

3. Results of Calibration

Table 1 gives the radiance temperature of the test lamp at 655.3 nm versus the test lamp current.

Table 2 gives the calibration uncertainties in relative to the International System of Units (SI Units). The expanded uncertainties $U = ku_c(T)$ (coverage factor $k = 2$) are two standard deviation estimates. Details on the estimation of these uncertainties are given in Ref. [1]. The NIST policy on uncertainty statements is described in Ref. [3].

4. General Information

The temperatures described in this report are radiance temperatures. The temperatures are stated with respect to the 1990 NIST Radiance Temperature Scale and correspond to temperatures when sighting on a Planckian radiator. See Ref. [2] for further details on the 1990 NIST Radiance Temperature Scale.

To maintain the highest accuracy, keep the lamp envelope clean and have the lamp recalibrated periodically. Appropriate calibration schedules vary with lamp and application and are best determined by the user.

The following was taken from Barber [4]. The radiance temperature below 1500 °C for gas-filled lamps and 900 °C for vacuum lamps will change as the room or ambient temperature changes. Barber observed a change in radiance temperature for a gas-filled lamp of between 0.3 °C and 0.4 °C per 1 °C change in room temperature at 1000 °C and 0.05 °C per 1 °C change in room temperature at 1400 °C. Changes of the same magnitude were observed at 700 °C and 900 °C for a vacuum lamp.

The results of this calibration apply only to the lamp referenced in this report. This report shall not be reproduced, except in full, without the written approval of the Radiance Temperature Measurements Calibration Service.

Prepared by:

Approved by:

Charles E. Gibson
Optical Technology Division
Physics Laboratory
(301) 975-2329

Robert D. Saunders
For the Director,
National Institute of
Standards and Technology
(301) 975-2355

Calibration Date: January 2, 1995
NIST Test No.: 844/123456-95

Page 2 of 5

REPORT OF CALIBRATION

35050C Radiance Temperature Standard, Ribbon Filament Lamp
Radiometric Systems Inc.

Model #: GE 30A/T24/6
Serial #: P000

References

- [1] C. E. Gibson, B. K. Tsai, and A. C. Parr, "Radiance Temperature Calibrations," *NIST Special Publication* **250-43** (1997).
- [2] K. D. Mielenz, R. D. Saunders, A. C. Parr, and J. J. Hsia, "The 1990 NIST Scales of Thermal Radiometry," *J. Res. Natl. Inst. Stand. Technol.*, **95**, 621-629, (1990).
- [3] B. N. Taylor and C. E. Kuyatt, "Guidelines for Evaluating and Expressing the Uncertainty of the NIST Measurement Results," *NIST Technical Note* **1297** (2nd ed., 1994).
- [4] C. R. Barber, "Factors affecting the reproducibility of brightness of tungsten strip lamps for pyrometry standardization," *J. Sci. Instr.* **23**, 238 (1946).

REPORT OF CALIBRATION

35050C Radiance Temperature Standard, Ribbon Filament Lamp
Radiometric Systems Inc.

Model #: GE 30A/T24/6
Serial #: P000

TABLE 1
Radiance Temperature of Lamp P000 at 655.3 nm

Radiance Temperature (1990 NIST) [°C]	Lamp Current [A]
1000	12.945
1100	14.220
1200	15.698
1300	17.358
1400	19.191
1500	21.155
1600	23.246
1700	25.449
1800	27.757
1900	30.175
2000	32.669
2100	35.254
2200	37.922

REPORT OF CALIBRATION35050C Radiance Temperature Standard, Ribbon Filament Lamp
Radiometric Systems Inc.Model #: GE 30A/T24/6
Serial #: P000

TABLE 2
Radiance Temperature Calibration Uncertainties

Source of Uncertainty	Type	Expanded Uncertainties ($k = 2$) [°C]				
		800 °C	1100 °C	1500 °C	1900 °C	2300 °C
1. Calibration of the reference radiance temperature lamp relative to the 1990 NIST Radiance Temperature Scale	B	0.12	0.19	0.32	0.48	0.67
2. Test lamp temperature determination	A	0.42	0.17	0.29	0.43	0.60
3. Current measurement	B	0.29	0.19	0.15	0.14	0.14
4. Mean effective wavelength measurement for the NIST Photoelectric Pyrometer	B	0.10	0.04	0.09	0.28	0.54
5. Test lamp alignment	B	0.09	0.14	0.24	0.36	0.51
6. 1990 NIST Radiance Temperature Scale relative to the Thermodynamic Temperature Scale	B	0.15	0.24	0.40	0.61	0.85
Overall uncertainty of test lamp calibration with respect to SI units		0.6	0.4	0.7	1.0	1.5

Appendix B: Optical pyrometer sample calibration report

In reply refer to: 844/123456-95

Radiometric Systems Inc.
Attn: John Doe
123 Calibration Court
Measurement City, MD 00000-0000

Subject: Optical Pyrometer Report of Calibration

Order No.: AB1234 dated January 1, 1995

Dear Mr. Doe:

Enclosed are results of the calibration you requested for Test Number 35040C Radiance Temperature Standard, Disappearing Filament Optical Pyrometer. Please refer to the above file number in any further communication concerning this calibration.

Sincerely,

Albert C. Parr, Chief
Optical Technology Division
Physics Laboratory

Enclosure:

One Report of Calibration

ACP/jah

NIST Cost Center 8443600

REPORT OF CALIBRATION

35040C Radiance Temperature Standard,
Disappearing Filament Optical Pyrometer

for

Model # Leeds and Northrup 8631-F, Serial # 0000000

Submitted by:

John Doe
Radiometric Systems Inc.
123 Calibration Court
Measurement City, MD 00000-0000

(See your Purchase Order No. AB1234 dated January 1, 1995)

1. Description of Calibration Item

An optical pyrometer was tested by the National Institute of Standards and Technology (NIST) as a standard of radiance temperature from 1149 °C (2100 °F) to 1482 °C (2700 °F).

2. Description of Calibration

The test pyrometer was measured in the NIST Radiance Temperature Calibration Facility using the equipment and procedures described in Ref. [1]. The NIST Photoelectric Pyrometer, which has a mean effective wavelength of 655.3 nm, was used to spectrally compare the variable-temperature blackbody to working standard SL20 to determine the radiance temperature of the blackbody at each calibration temperature. The portion of the test pyrometer lamp filament below the wire marker was matched in brightness to that of the blackbody and the test pyrometer scale reading was recorded at each calibration temperature. The NIST Photoelectric Pyrometer is described in the Ref. [1]

The variable-temperature blackbody has an estimated emissivity of 0.99. The blackbody cavity is a single piece of graphite, specially tapered for temperature uniformity. This graphite tube is cylindrically hollow on both ends to form two 2.54 cm diameter cavities with a common partition in the center. One cavity is used for temperature control and the other as a blackbody source.

Laboratory Environment:

Temperature: 23 °C ± 1 °C

Relative Humidity: 46% ± 5%

Calibration Date: January 2, 1995

NIST Test No.: 844/123456-95

Page 1 of 5

REPORT OF CALIBRATION

35040C Radiance Temperature Standard, Optical Pyrometer
Radiometric Systems Inc.

Model #: L&N 8631-F
Serial #: 0000000

The test pyrometer was operated according to the manufacturer's operating instructions. The pyrometer was calibrated from 1149 °C (2100 °F) to 1482 °C (2700 °F) on the high range (H). The pyrometer was aligned so that its optical axis, when viewing the variable-temperature blackbody, coincided with the geometrical center of the blackbody. The pyrometer was measured at a distance of 81.7 cm from the front surface of the objective lens holder to the blackbody cavity partition and was focussed on the center of the partition.

3. Results of Calibration

Table 1 gives the radiance temperature at 655.3 nm versus the test pyrometer scale reading.

Table 2 gives the calibration uncertainties relative to the International System of Units (SI Units).

The expanded uncertainties $U = ku_c(T)$ (coverage factor $k = 2$) are two standard deviation estimates. Details on the estimation of these uncertainties are given in Ref. [1]. The NIST policy on uncertainty statements is described in Ref. [3].

4. General Information

The temperatures described in this report are radiance temperatures. The temperatures are stated with respect to the 1990 NIST Radiance Temperature Scale and correspond to temperatures when sighting on a Planckian radiator. See Ref. [2] for further details on the 1990 NIST Radiance Temperature Scale.

To maintain the highest accuracy, keep the optical components clean and have the pyrometer recalibrated periodically. Appropriate calibration schedules vary with instrument and application and are best determined by the user.

The results of this calibration apply only to the pyrometer referenced in this report. This report shall not be reproduced, except in full, without the written approval of the Radiance Temperature Measurements Calibration Service.

Prepared by:

Approved by:

Charles E. Gibson
Optical Technology Division
Physics Laboratory
(301) 975-2329

Robert D. Saunders
For the Director,
National Institute of
Standards and Technology
(301) 975-2355

Calibration Date: January 2, 1995
NIST Test No.: 844/123456-95

Page 2 of 5

REPORT OF CALIBRATION

35040C Radiance Temperature Standard, Optical Pyrometer
Radiometric Systems Inc.

Model #: L&N 8631-F
Serial #: 0000000

References

- [1] C. E. Gibson, B. K. Tsai, and A. C. Parr, "Radiance Temperature Calibrations," *NIST Special Publication* **250-43** (1997).
- [2] K. D. Mielenz, R. D. Saunders, A. C. Parr, and J. J. Hsia, "The 1990 NIST Scales of Thermal Radiometry," *J. Res. Natl. Inst. Stand. Technol.*, **95**, 621-629, (1990).
- [3] B. N. Taylor and C. E. Kuyatt, "Guidelines for Evaluating and Expressing the Uncertainty of the NIST Measurement Results," *NIST Technical Note* **1297** (2nd ed., 1994).

REPORT OF CALIBRATION

35040C Radiance Temperature Standard, Optical Pyrometer
Radiometric Systems Inc.

Model #: L&N 8631-F
Serial #: 0000000

TABLE 1
Radiance Temperature at 655.3 nm

Nominal Temperature [°C]	Blackbody Temperature (1990 NIST) [°C]	Blackbody Temperature [°F]	Pyrometer Scale Reading Value [°F]	Pyrometer Correction [°F]
High range (H)				
1148.9	1149.2	2100.6	2098	3
1315.6	1316.5	2401.7	2401	1
1482.2	1482.9	2701.2	2700	1

This equation was used to convert the blackbody temperature T from °C to °F:

$$T/^{\circ}\text{F} = (1.8 \times T/^{\circ}\text{C}) + 32$$

REPORT OF CALIBRATION

35040C Radiance Temperature Standard, Optical Pyrometer
Radiometric Systems Inc.

Model #: L&N 8631-F

Serial #: 0000000

TABLE 2
Radiance Temperature Calibration Uncertainties

Source of Uncertainty	Type	Expanded Uncertainties ($k = 2$) [°C (°F)]		
		1148.9 °C (2100 °F)	1315.6 °C (2400 °F)	1482.2 °C (2700 °F)
1. Calibration of the variable temperature blackbody relative to the 1990 NIST Radiance Temperature Scale	B	0.3 (0.6)	0.4 (0.7)	0.4 (0.7)
2. Mean effective wavelength measurement for the NIST Photoelectric Pyrometer	B	0.0 (0.0)	0.0 (0.0)	0.1 (0.2)
3. Blackbody uniformity	B	0.2 (0.4)	0.2 (0.4)	0.2 (0.4)
4. Test pyrometer temperature determination	A	1.1 (2.0)	1.1 (2.0)	2.2 (4.0)
5. 1990 NIST Radiance Temperature Scale relative to the Thermodynamic Temperature Scale	B	0.3 (0.5)	0.3 (0.6)	0.4 (0.7)
Overall uncertainty of test lamp calibration with respect to SI units		1.2 (2.2)	1.2 (2.2)	2.3 (4.1)

Note: The expanded uncertainties in parentheses are in °F. The following equation was used to convert the blackbody temperature T from °C to °F: $T/^{\circ}\text{F} = (1.8 \times T/^{\circ}\text{C}) + 32$

Appendix C: Radiation thermometer sample calibration report

In reply refer to: 844/123456-95

Radiometric Systems Inc.
Attn: John Doe
123 Calibration Court
Measurement City, MD 00000-0000

Subject: Radiation Thermometer Report of Test

Order No.: AB1234 dated January 1, 1995

Dear Mr. Doe:

Enclosed are results of the test you requested for Test Number 35070S Special Test of Radiation Thermometer. Please refer to the above file number in any further communication concerning this test.

Sincerely,

Albert C. Parr, Chief
Optical Technology Division
Physics Laboratory

Enclosure:

One Report of Test

ACP/jah

NIST Cost Center 8443600

REPORT OF TEST

35070S Special Test of Radiation Thermometer

for

Model # Minolta-Land Cyclops 52, Serial # 00000000

Submitted by:

John Doe
Radiometric Systems Inc.
123 Calibration Court
Measurement City, MD 00000-0000

(See your Purchase Order No. AB1234 dated January 1, 1995)

1. Description of Test Item

A radiation thermometer was tested by the National Institute of Standards and Technology (NIST) as a standard of radiance temperature from 800 °C to 2300 °C. The test unit was measured with clear filter (Serial Number 0) and dark filter (Serial Number 0).

2. Description of Test

The test thermometer was measured in the NIST Radiance Temperature Calibration Facility using the equipment and procedures described in Ref. [1]. The NIST Photoelectric Pyrometer, which has a mean effective wavelength of 655.3 nm, was used to spectrally compare the variable-temperature blackbody to working standard SL20 to determine the radiance temperature of the blackbody at each calibration temperature. The test thermometer display panel value was recorded at each calibration temperature. The NIST Photoelectric Pyrometer is described in the Ref. [1]

The variable-temperature blackbody has an estimated emissivity of 0.99. The blackbody cavity is a single piece of graphite, specially tapered for temperature uniformity. This graphite tube is cylindrically hollow on both ends to form two 2.54 cm diameter cavities with a common partition in the center. One cavity is used for temperature control and the other as a blackbody source.

Laboratory Environment:

Temperature: 23 °C ± 1 °C

Relative Humidity: 46% ± 5%

Test Date: January 2, 1995

NIST Test No.: 844/123456-95

Page 1 of 5

REPORT OF TEST

35070S Special Test of Radiation Thermometer
Radiometric Systems Inc.

Model #: Minolta-Land Cyclops 52
Serial #: 00000000

The thermometer was aligned so that its optical axis, when viewing the variable-temperature blackbody, coincided with the geometrical center of the blackbody. The test thermometer was operated according to the manufacturer's operating instructions and measured with the following conditions.

Temperatures measured:	800 °C to 1200 °C with the clear filter 1200 °C to 2300 °C with the dark filter
Distance:	100 cm from the thermometer alignment mark to the blackbody cavity partition
Operating mode:	continuous mode
Calculated target size:	5.6 mm diameter focused on the center of the blackbody cavity partition
Emissivity setting:	1

3. Results of Test

Table 1 gives the radiance temperature at 1000 nm versus the test thermometer display panel value.

Table 2 gives the test uncertainties relative to the International System of Units (SI Units). The expanded uncertainties $U = ku_c(T)$ (coverage factor $k = 2$) are two standard deviation estimates. Details on the estimation of these uncertainties are given in Ref. [1]. The NIST policy on uncertainty statements is described in Ref. [3].

4. General Information

The temperatures described in this report are radiance temperatures. The temperatures are stated with respect to the 1990 NIST Radiance Temperature Scale and correspond to temperatures when sighting on a Planckian radiator. See Ref. [2] for further details on the 1990 NIST Radiance Temperature Scale.

To maintain the highest accuracy, keep the optical components clean and have the thermometer recalibrated periodically. Appropriate calibration schedules vary with instrument and application and are best determined by the user.

The results of this test apply only to the thermometer referenced in this report. This report shall not be reproduced, except in full, without the written approval of the Radiance Temperature Measurements Calibration Service.

Test Date: January 2, 1995
NIST Test No.: 844/123456-95

Page 2 of 5

REPORT OF TEST

35070S Special Test of Radiation Thermometer
Radiometric Systems Inc.

Model #: Minolta-Land Cyclops 52
Serial #: 00000000

Prepared by:

Approved by:

Charles E. Gibson
Optical Technology Division
Physics Laboratory
(301) 975-2329

Robert D. Saunders
For the Director,
National Institute of
Standards and Technology
(301) 975-2355

References

- [1] C. E. Gibson, B. K. Tsai, and A. C. Parr, "Radiance Temperature Calibrations," *NIST Special Publication* **250-43** (1997).
- [2] K. D. Mielenz, R. D. Saunders, A. C. Parr, and J. J. Hsia, "The 1990 NIST Scales of Thermal Radiometry," *J. Res. Natl. Inst. Stand. Technol.*, **95**, 621-629, (1990).
- [3] B. N. Taylor and C. E. Kuyatt, "Guidelines for Evaluating and Expressing the Uncertainty of the NIST Measurement Results," *NIST Technical Note* **1297** (2nd ed., 1994).

REPORT OF TEST

35070S Special Test of Radiation Thermometer
Radiometric Systems Inc.

Model #: Minolta-Land Cyclops 52
Serial #: 00000000

TABLE 1
Radiance Temperature at 1000 nm

Nominal Temperature [°C]	Blackbody Temperature (1990 NIST) [°C]	Thermometer Display Panel Value [°C]	Thermometer Correction [°C]
Clear Filter			
800	798.7	801	-2.3
900	899.1	900	-0.9
1000	1000.2	1000.5	-0.3
1100	1100.6	1100.5	0.1
1200	1200.5	1199.5	1.0
Dark filter			
1200	1200.4	1200	0.4
1300	1299.4	1298	1.4
1400	1399.6	1396.5	3.1
1500	1500.3	1497	3.3
1600	1600.6	1597	3.6
1700	1700.5	1696.5	4.0
1800	1801.7	1797	4.7
1900	1899.5	1894	5.5
2000	1999.9	1992	7.9
2100	2099.9	2091.5	8.4
2200	2200.7	2191	9.7
2300	2301.2	2290	11.2

REPORT OF TEST

35070S Special Test of Radiation Thermometer
Radiometric Systems Inc.

Model #: Minolta-Land Cyclops 52
Serial #: 00000000

TABLE 2
Radiance Temperature Test Uncertainties

Source of Uncertainty	Type	Expanded Uncertainties ($k = 2$) [°C]				
		800 °C	1100 °C	1500 °C	1900 °C	2300 °C
1. Calibration of the variable temperature blackbody relative to the 1990 NIST Radiance Temperature Scale	B	0.3	0.3	0.4	0.6	0.9
2. Mean effective wavelength measurement for the NIST Photoelectric Pyrometer	B	0.1	0.0	0.1	0.3	0.5
3. Blackbody temperature correction	B	0.1	0.2	0.4	0.6	0.8
4. Blackbody uniformity	B	0.2	0.2	0.2	0.2	0.2
5. Test thermometer temperature determination	A	2.0	2.0	2.0	2.0	2.0
6. 1990 NIST Radiance Temperature Scale relative to the Thermodynamic Temperature Scale	B	0.1	0.2	0.4	0.6	0.9
Overall uncertainty of test lamp calibration with respect to SI units		2.0	2.1	2.1	2.3	2.5

Radiance Temperature Measurements Laboratory Procedures Manual



Safety Issues

Safety glasses are required in lab at all times.

The following hazards have been identified in this lab. Safety measures are in place to reduce risk of injury from hazards below.

BRIGHT SOURCE

Staring at source is to be avoided because staring can cause severe eye damage. Using correct neutral density filter when viewing is required

BURN

To avoid burns, it is wise to keep body and clothes away from lamps and blackbodies.

COLLISION

The boxcar path should be kept clear of collision with any objects, and cables should not be hung from boxcar at any time.

ELECTROCUTION

Since high voltage connections are present, electrical parts should not be touched at all.

LASER

As with other bright sources, staring into laser is to be avoided because staring can cause severe eye damage. When lasers are not in use, block laser with opaque object and turn off laser when it is not in use.

Symbols, similar to one below, are used throughout this document to alert operator to potential hazards.



Table of Contents

- i. Safety Issues
- 1. Running General Data Taking Program (RTMLP1)
- 2. NIST Photoelectric Pyrometer (PEP) Operation (RTMLP2)
- 3. Gold-point Blackbody (GPBB) Alignment (RTMLP3)
- 4. Gold-point Blackbody Operation (RTMLP4)
- 5. Ribbon Filament Lamp Preparation (RTMLP5)
- 6. Ribbon Filament Lamp Alignment (RTMLP6)
- 7. Ribbon Filament Lamp Operation (RTMLP7)
- 8. Ribbon Filament Lamp Calibration Procedure (RTMLP8)
- 9. Ribbon Filament Lamp Data Reduction (RTMLP9)
- 10. Variable Temperature Blackbody (BB) Alignment (RTMLP10)
- 11. Variable Temperature Blackbody Operation (RTMLP11)
- 12. Disappearing Filament Optical Pyrometer Alignment (RTMLP12)
- 13. Radiation Thermometer (RT) Alignment (RTMLP13)
- 14. Radiation Thermometer (Optical Pyrometer) Calibration Procedure (RTMLP14)
- 15. Radiation Thermometer Data Reduction (RTMLP15)

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Running General Data Taking Program

Written by: Charles Gibson

Procedure: RTMLP1

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Running the General Data Taking Program

1. Turn on computer
2. Enter password to work online or click on Cancel button to work offline
3. Load Visual Basic program: *c:\programs\datatak\pyro.mak* (do one of following)
 - a. Click on Pyro.mak shortcut icon
 - b. Double click on folder *c:\programs\datatak\pyro.mak* using Windows Explorer
4. Run program (do one of following)
 - a. Press F5
 - b. Click on ► tool bar icon
 - c. Click on menu item "Run" then "Start"
5. Run *Hardware.frm* to set up equipment with settings used by other programs
6. Run *Software.frm* to input data used by other programs
7. Run all other programs from these forms

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
NIST Photoelectric Pyrometer (PEP) Operation

Written by: Charles Gibson

Procedure: RTMLP2

Approved by: _____

Date: _____

Revised by: _____

Date: _____

NIST Photoelectric Pyrometer (PEP) Operation

1. Turn on equipment
 - a. PMT thermoelectric cooler
 - b. Digital multimeters
 - c. Motor controller
 - d. High voltage power supply
 - e. Current amplifier
 - f. Environment monitor

CAUTION

LASER

2. Check alignment of optics
 - a. Align laser on source table so that beam goes through center of pin hole aperture and field stop
 - b. Check that collimating lens is centered on and normal to optical axis (retroreflection checked at field stop)
 - c. Check that aperture stop is centered on and normal to optical axis (place microscope slide against aperture stop to check retroreflection at field stop)
 - d. Check scope alignment (see step 3)
 - e. Check that lamp beam is centered on aperture stop
 - f. Record aperture stop micrometer position in lab notebook

CAUTION

BRIGHT SOURCE

BURN

COLLISION

ELECTROCUTION

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP2: Photoelectric Pyrometer (PEP) Operation

3. Check scope alignment (scope is used to define optical axis of PEP)
 - a. Turn on SL20 or other lamp in focal plane of PEP (see RTMLP7)
 - b. Move lamp to PEP
 - c. Close PMT shutter
 - d. Open enclosure doors
 - e. Open scope door
 - f. Place 3/4" aperture on objective lens
 - g. Put in field stop
 - h. Put in prism
 - i. Focus reticle to operator's eye
 - j. Focus scope on field stop
 - k. Rotate reticle so that #5 vertical scale line is parallel to side of field stop opening
 - l. Position reticle center (intersection of #5 vertical scale line and horizontal scale line) on field stop opening
 - m. Remove prism
 - n. Remove field stop
 - o. Put in prism
4. Set up PEP for measurements
 - a. Remove aperture from objective lens
 - b. Remove prism
 - c. Close scope door
 - d. Put in field stop
 - e. Close all enclosure doors
 - f. Position baffles at front of PEP
 - g. Open PMT shutter
 - h. Begin calibration (see (RTMLP8) for lamps or (RTMLP14) for radiation thermometers)

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Gold-point Blackbody (GPBB) Alignment

Written by: Charles Gibson

Procedure: RTMLP3

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Gold-point Blackbody (GPBB) Alignment



1. Align laser on source table so that beam goes through center of pin hole aperture and field stop
2. Move first laser to laser on PEP table
3. Align second laser so that it coincides with first laser
4. Remove outer most graphite ring
5. Put in alignment jig
6. Translate GPBB to center jig on laser
7. Rotate GPBB to make jig normal to laser (retroreflect laser beam)
8. Turn on GPBB (see RTMLP4)
9. Radiometrically center GPBB
 - a. Center X along horizontal axis
 - b. Center Y along vertical axis
 - c. Center $\angle X$ about horizontal axis
 - d. Center $\angle Y$ about vertical axis
 - e. Make ± 8 turns along optical axis
 - f. Repeat steps a and b until centered
 - g. Check centering
 - i. Monitor BB signal on graph logger
 - ii. Hold card in front of BB opening
 - iii. Check for equal signal drop-offs when blocking top, bottom, left, and right portions of BB opening

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Gold-point Blackbody Operation

Written by: Charles Gibson

Procedure: RTMLP4

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Gold-point Blackbody (GPBB) Operation

CAUTION
BRIGHT SOURCE
BURN
ELECTROCUTION

1. Close GPBB door
2. Turn on cold water at flow rate of 10 GPH
3. Set argon at flow rate of 54 mL/min (set black ball to 50)
4. Go to rack of power supplies
5. Turn on digital multimeter
6. Check that correct cable is connected to GPBB power supply connector
7. Switch on power supply (labeled HPPS)
8. Run lamp current monitoring program
 - a. Turn on computer
 - b. Type 2 and press [Enter] at DOS prompt to load and run lamp monitoring program
 - c. Prepare GPBB
 - i. Ramp up GPBB
 - (1) Ramp up GPBB to about 1040 °C
 - (2) Set current at 2 A
 - (3) Wait 3 h
 - (4) Ramp up to 8 A in 4 h (ramping rate 1.5 A/h)
 - (5) Set GPBB temperature to about 1040 °C in 8 h to 10 h
 - ii. Press 8 to ramp up GPBB
 - iii. Type (start time) and press [Enter]
 - iv. program will start ramping up GPBB at start time

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP4: Gold-point Blackbody Operation

- d. Cycle GPBB through melt and freeze cycles
 - i. Press 3 to monitor lamps and read GPBB temperature
 - ii. Press [Esc] to return to main menu
 - iii. Set current to 8.5 A and begin GPBB melt cycle when set point falls below 1050 °C
 - iv. Set current to 7.95 A and begin GPBB freeze cycle when set point rises above 1071 °C
- 9. Set argon at flow rate of 235 mL/min (set black ball to 120)
- 10. Open GPBB door
- 11. Turn off GPBB
 - a. Close GPBB door
 - b. Set argon at flow rate of 54 mL/min (set black ball to 50)
 - c. Turn of power to GPBB
 - d. Leave water and argon on overnight until GPBB is cool

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Ribbon Filament Lamp Preparation

Written by: Charles Gibson
Approved by: _____
Revised by: _____

Procedure: RTMLP5
Date: _____
Date: _____

Ribbon Filament Lamp Preparation

CAUTION
BRIGHT SOURCE
BURN
ELECTROCUTION

1. Lamp inspection
 - a. Check filament continuity
 - b. Check clarity of envelope (and window)
 - c. Check flatness of filament
 - d. Check that filament is within 2 mm of being centered in envelope (or window)
2. Clean lamp envelope
3. Anneal lamp
 - a. Anneal gas lamp
 - i. Leak test lamp
 - (1) Set lamp to 1700 °C (PMT signal = 2.5 V on 8 range, high voltage = 449 V) (About 25 A dc)
 - (2) Turn off lamp
 - (3) Wait 15 min to allow air to enter any leaks
 - (4) Set lamp to 2350 °C (PMT signal = 1.2 V on 7 range, high voltage = 394 V)(About 45 A dc)
 - (5) Wait 15 min
 - ii. Anneal lamp for 2 h
 - iii. Age lamp for 50 h at 2300 °C (About 42 A dc)

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP5: Ribbon Filament Lamp Preparation

- b. Anneal vacuum lamp
 - i. Leak test lamp
 - (1) Set lamp to 1000 °C (PMT signal = 5.3 V on 9 range, high voltage = 730 V)
 - (2) Turn off lamp
 - (3) Wait 15 min to allow air to enter any leaks
 - (4) Set lamp to 1450 °C (PMT signal = 5.1 V on 9 range, high voltage = 449 V)
 - (5) Wait 15 min
 - ii. Anneal lamp for 14 h
- 4. Locate arrow on lamp
 - a. Set lamp to 1700 °C for gas lamps or 1200 °C for vacuum lamps (for 1200 °C, PMT signal = 5.5 V on 8 range, high voltage = 730 V)
 - b. Use plumb line to make notch side of filament vertical
 - c. Translate lamp horizontally so that reticle center bisects filament lengthwise
 - d. Translate lamp vertically so that reticle center is positioned at center of notch opening
 - e. Translate lamp along optical axis until filament is in focus
 - f. Rotate lamp about vertical axis so that both edges of filament are in focus
 - g. Use graph logger to locate angular settings with smallest change in radiance with respect to rotation about horizontal and vertical axes
 - h. First person looks through alignment scope and guides second person to notch
 - i. Second person makes dot on lamp envelope at notch
 - j. Draw picture of notch in lab notebook showing its relationship to dot
 - k. Draw arrow showing relationship of alignment arrow to dot
 - l. Repeat steps for each lamp

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Ribbon Filament Lamp Alignment

Written by: Charles Gibson
Approved by: _____
Revised by: _____

Procedure: RTMLP6
Date: _____
Date: _____

Ribbon Filament Lamp Alignment

CAUTION
BRIGHT SOURCE
BURN
COLLISION
ELECTROCUTION

1. Check PEP scope alignment (see RTMLP2)
2. Place 3/4" aperture on objective lens
3. Turn on lamp (see RTMLP7)
4. Move lamp to PEP optical axis
5. Use plumb line to make notch side of filament vertical
6. Translate lamp horizontally so that reticle center bisects filament lengthwise
7. Translate lamp vertically so that reticle center is positioned at height of notch opening
8. Translate lamp along optical axis until filament is in focus
9. Rotate lamp about vertical axis so that both edges of filament are in focus
10. Translate lamp horizontally and vertically to reticle center
11. Place 1/4" aperture on objective lens
12. Place mirror behind lamp alignment arrow
13. Rotate lamp about horizontal axis so that tip of arrowhead is positioned at height of notch opening
14. Rotate lamp about vertical axis so that tip of arrowhead is positioned at notch opening
15. Repeat steps 13 and 14 until tip of arrowhead is centered at notch opening
16. Translate lamp horizontally and vertically to reticle center at height of notch
17. Repeat steps 3 to 16 for standard and test lamps
18. Set up PEP for measurements (see RTMLP2)

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Ribbon Filament Lamp Operation

Written by: Charles Gibson

Procedure: RTMLP7

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Ribbon Filament Lamp Operation

CAUTION
BRIGHT SOURCE
BURN
ELECTROCUTION

1. Use Kimwipe to avoid touching lamp envelope
2. Place lamp in socket
 - a. Check lamp polarity (center contact for screw base, longer filament support for mogul base at positive potential)
 - b. Set alignment arrow or mark away from PEP
3. Go to rack of power supplies
4. Turn on digital multimeter
5. Turn on GP-IB Programmer(s) for supply(s) to be used
 - a. cables and programmers are numbered with corresponding power supply number
 - b. programmers are two channel 12 bit D/A converters
6. Check that correct cable is connected to power supply
7. Switch on power supply
8. Run lamp current monitoring program
 - a. Turn on computer
 - b. Type 2 and press [Enter] at DOS prompt to load and run lamp monitoring program
 - c. Set lamp current
 - i. Press 1 to set power supply current
 - ii. Type (supply #) and press [Enter]
 - iii. Type (lamp current) and press [Enter]
 - d. Monitor lamp currents
 - i. Press 3
 - ii. Press [Esc] to return to main menu

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP7: Ribbon Filament Lamp Operation

- e. Turn off supply
 - i. Press 2
 - ii. Type (supply #) and press [Enter]
- f. Change lamp current
 - i. same as (1)
- g. Set multiple lamp currents
 - i. Press A to select power supplies to turn on
 - ii. Type (lamp current) and press [Enter] to select supplies to turn on
 - iii. Type 0 and press [Enter] to skip supply
- h. Turn off all supplies
 - i. Press B to turn off all lamps
 - ii. Type OFF and press [Enter] to turn off all lamps
- i. Trouble-shoot
 - i. Check program locking
 - (1) Type Ctrl [Break]
 - (2) Type GOTO 20, press [Enter] to return to main menu
- 9. Align lamps (see RTMLP6)
- 10. Use **EMERGENCY TURN OFF** procedure when necessary
 - a. Go to supply
 - b. Turn CURRENT knob counter clockwise to turn down lamp current
 - c. Turn off supply
 - d. Turn off GP-IB Programmer for supply

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Ribbon Filament Lamp Procedure

Written by: Charles Gibson

Procedure: RTMLP8

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Ribbon Filament Lamp Calibration Procedure

1. Review calibration summary
 - a. Align and measure test lamp at 1700 °C
 - b. Measure 2300 °C to 800 °C in 100 °C steps
 - c. Repeat 1700 °C, 2300 °C, 1900 °C, 1500 °C, 1100 °C, and 800 °C
2. Run general data taking program (see RTMLP1)
3. Align standard and test lamps (see RTMLP6))
4. Run calibration program
 - a. Click on "Calibrations" from *Hardware.frm* menu to load *Rtcal.frm*
 - b. Select "Lamps and blackbodies" button
 - c. Enter data filename
 - d. Enter test filename suffix\$
 - i. Create test filenames from names on *Software.frm*
 - e. Select Yes and use suggested ranges (PMT high voltage and amplifier ranges are part of *Rtcal.frm* program)
 - f. Click on "Start" button
 - g. Enter data on *Software.frm*
 - i. Enter lamp and blackbody data
 - (1) Enter data for reference and test lamps
 - (2) Select "Use" button
 - (3) Enter lamp position
 - (4) Enter lamp name
 - (5) Select "Reference" button
 - (6) Enter lamp current (no necessary, lamp current automatically read from lamp monitoring computer)
 - (7) Enter power supply number
 - ii. Select temperature scale
 - iii. Select measurement wavelength [nm]
 - iv. Enter number of repeats (number of test/standard measurements)

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP8: Ribbon Filament Lamp Procedure

- v. Enter number of readings (number of readings per DMM trigger)
- vi. Select lamp type
- vii. Select detector channel
- h. Click on "Save Configuration" button
- i. Select "Continue" button to return to *Rtcal.frm*
- j. Enter Standard temperature in °C in input box
- k. Enter Nominal radiance temperature in input box
- l. Read lamp currents automatically
- m. Adjust test lamp current so that test lamp signal matches that indicated on *Suggest.frm*
- n. Wait for lamp to stabilize (5 min, 10 min, 15 min, or 20 min depending on measurement temperature)
- o. Check lamp alignment
- p. Move source table to reference source and measure full scale signal
- q. Move source table to dark signal position and measure dark signal
- r. Move source table to test source and measure full scale signal
- s. Move source table to dark signal position and measure dark signal on test source amplifier range, if reference and test source amplifier ranges are different
- t. Repeat steps r to s for each test source
- u. Repeat steps p to s "number of repeats" times
- v. Select next action and then select "Continue" button
 - i. Repeat measurement with no changes
 - ii. Measure with new temperature
 - iii. Measure with new software changes
 - iv. Stop measurements
- w. Repeat steps k to v for each measurement point

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP8: Ribbon Filament Lamp Procedure

5. Store calibration data in log file "CALMMDD.YY"
6. Use sample log file below:

```

Date 08-15-1995   Time 12:49:14 PM
File Name = CAL0815.95
PMT Detector   PMT HV = 449   NRDS = 5   Integration Time = 8.35 s
Wavelength = 655.3 nm   Nominal Rad Temp = 1700 °C

Name           SL20           P51

Rad. Temp      1255.07 °C      1700.37 °C
Current        7.7788 A        26.960 A
Sp. Rad.       5.696E+02      1.456E+04
Ratio          1.000E+00      2.556E+01
Std. Ratio     0.01%          0.04%
Amp Range      9             8
Full Scale     1.052E+00 ±0.00%    2.685E+00 ±0.00%
               1.052E+00 ±0.00%    2.683E+00 ±0.00%
               1.052E+00 ±0.00%    2.684E+00 ±0.00%
Dark Sig.      1.501E-04          -8.035E-05
Mean           1.052E+00      2.684E+00
Std(Mean)      0.01%          0.04%
Range Cf.      9.9858E-02      1.0000E+00
SS Cf.         1.000E+00      1.000E+00
Corr. Mean     1.050E-01      2.684E+00
Position #     1             6
Table Pos.     37500          246000
  
```

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP8: Ribbon Filament Lamp Procedure

7. Store lamp data in test file "XXXXMMDD.YY"
8. Use sample test file below:

P510815.95			
08-15-1995 12:49:14 PM			
Tnom °C	Tmeas °C	Imeas °C	Test/Standard
1700	1700.37	26.960	2.5558E+01
2300	2300.22	43.340	3.4164E+02
2200	2200.29	40.330	2.4203E+02
2100	2100.04	37.440	1.6636E+02
2000	2000.04	34.670	1.1075E+02
1900	1900.16	31.990	7.1053E+01
1800	1800.13	29.420	4.3646E+01

9. Reduce data (see RTMLP9)

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Ribbon Filament Lamp Data Reduction

Written by: Charles Gibson

Procedure: RTMLP9

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Ribbon Filament Lamp Data Reduction

1. Print out test file
2. Look in lab notebook and make changes on printout and then edit file
3. Load Visual Basic program: *c:\programs\lmpredn\Rtred.mak* (do one of following)
 - a. Click on Rtred.mak shortcut icon
 - b. Double click on folder *c:\programs\lmpredn\Rtred.mak* using Windows Explorer
4. Run program (do one of following)
 - a. Press F5
 - b. Click on ► tool bar icon
 - c. Click on menu item "Run" then on "Start"
5. Select number of lamps to reduce
6. Enter 4 character serial number for each lamp
7. Enter data file suffix
8. Select directory that contains data file
9. Select degrees of freedom for curve fit
10. Click on "Run" button
11. Print data in list box
12. Calculate mean (T_{meas}) and standard deviation are calculated
13. Compare standard deviation to control limit
14. Correct lamp current to nominal temperature using following method, because measured temperatures are different from nominal temperatures (T_{nom})
 - a. Fit data to fifth degree polynomial using least squares approximation
 - b. Define independent variable: $T_{\text{mean}} = (T_{\text{meas}}(n-1) + T_{\text{meas}}(n))/2$
 - c. Define dependent variable: $\Delta I/\Delta T = (I_{\text{meas}}(n-1) - I_{\text{meas}}(n))/(T_{\text{meas}}(n-1) - T_{\text{meas}}(n))$
 - d. Calculate dI/dT from fit
 - e. Define corrected lamp current: $I_{\text{corr}} = I_{\text{meas}} + dI/dT * (T_{\text{nom}} - T_{\text{meas}})$
15. Store reduced data in log file "REDMMDD.YY"
16. Use sample log file on next page:

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP9: Ribbon Filament Lamp Data Reduction

Reducing 'C:\PYRO\P510815.95' 08-17-1995 11:26:56 AM

P510815.95

08-15-1995 12:49:14 PM

Tnom °C	Tmeas °C	Imeas A	Test/Standard
1700	1700.37	26.96	2.5558E+01
2300	2300.22	43.34	3.4164E+02
2200	2200.29	40.33	2.4203E+02
2100	2100.04	37.44	1.6636E+02
2000	2000.04	34.67	1.1075E+02
1900	1900.16	31.99	7.1053E+01
1800	1800.13	29.42	4.3646E+01

Tnom °C	Tmean °C	Imeas A	Std Error	Limit	Flag	Nrpt
800	800.82	11.850	0.09	0.21		3
900	899.58	12.720	0.00	0.19		1
1000	1000.04	13.800	0.00	0.15		1
1100	1100.02	15.100	0.00	0.09		1
1200	1199.63	16.620	0.00	0.10		1

Tnom °C	Imeas A	Tmeas °C	Tmean °C	DelI/DelT	dI/dT	DelT	DelI	Icorr A
0800	11.850	0800.82	0850.20	0.0088	0.0079	-0.82	-0.007	11.843
0900	12.720	0899.58	0949.81	0.0108	0.0097	0.42	0.004	12.724
1000	13.800	1000.04	1050.03	0.0130	0.0119	-0.04	0.000	13.800
1100	15.100	1100.02	1149.83	0.0153	0.0141	-0.02	0.000	15.100
1200	16.620	1199.63	1249.74	0.0174	0.0163	0.37	0.006	16.626

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP9: Ribbon Filament Lamp Data Reduction

Coefficients used to calculate dI/dT from a 5th degree polynomial

$C(0) = 0.063517488043157$
 $C(1) = -2.4035523727757E-04$
 $C(2) = 3.56720919376003E-07$
 $C(3) = -2.29083987235765E-10$
 $C(4) = 6.90646490071546E-14$
 $C(5) = -7.96328820986401E-18$

800	11.843
900	12.724
1000	13.800
1100	15.100
1200	16.626
1300	18.363
1400	20.289
1500	22.376
1600	24.599
1700	26.951
1800	29.417
1900	31.988
2000	34.669
2100	37.441
2200	40.321
2300	43.331

Data is stored in 'C:\PYRO\P510895.TXT'

Reduction of file 'C:\PYRO\P510815.95 ' is complete

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP9: Ribbon Filament Lamp Data Reduction

17. Save T_{nom} and I_{corr} in data file "XXXXMMYY.TXT"
18. Use sample test data file below:

C:\PYRO\P510895.TXT	
800	11.843
900	12.724
1000	13.800
1100	15.100
1200	16.626
1300	18.363
1400	20.289
1500	22.376
1600	24.599
1700	26.951
1800	29.417
1900	31.988
2000	34.669
2100	37.441
2200	40.321
2300	43.331

19. Compare current calibration to previous calibration
20. Click on "Compare calibrations" from menu
21. Select previous calibration filename
22. Select current calibration filename
23. Select temperature scale ($^{\circ}\text{C}$, $^{\circ}\text{F}$, or K)
24. Select international temperature scale (ITS-48, ITS-68, ITS-90)
25. Enter wavelength (nm) for previous data file
26. Enter wavelength (nm) for current data file
27. Click on "Run" button
28. Enter data file name "XXXXmyy.MYY"
 - a. Define myy: month and year of previous calibration
 - b. Define MYy: month and year of current calibration
29. Store comparison data in log file "CMPMMDD.YY"
30. Use sample log file on next page:

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP9: Ribbon Filament Lamp Data Reduction

08-17-1995 11:50:32

Previous calibration data file = 'c:\files\p510393.txt'

Current calibration data file = 'c:\files\p510895.txt'

Comparing previous to current calibration

Temp [°C]	655.3nm		655.3nm		dI	DI/DT	WL Corr.	Change	Total	%Change Rad.	Scale Corr.
	Old [I] ITS-90	New [I] ITS-90	New [I] ITS-90	dI ITS-90				Temp. ITS-90	Change ITS-90		
	[A]	[A]	[A]	[A]	[A/°C]						[°C]
800	11.847	11.843	11.843	0.004	0.0078	0.00	0.51	0.51	0.98%	0.00	
900	12.724	12.724	12.724	0.000	0.0096	0.00	0.00	0.00	0.00%	0.00	
1000	13.799	13.800	13.800	-0.001	0.0119	0.00	-0.08	-0.08	-0.11%	0.00	
1100	15.098	15.100	15.100	-0.002	0.0142	0.00	-0.14	-0.14	-0.16%	0.00	
1200	16.623	16.626	16.626	-0.003	0.0164	0.00	-0.18	-0.18	-0.18%	0.00	
1300	18.361	18.363	18.363	-0.002	0.0185	0.00	-0.11	-0.11	-0.10%	0.00	
1400	20.287	20.289	20.289	-0.002	0.0201	0.00	-0.10	-0.10	-0.08%	0.00	
1500	22.372	22.376	22.376	-0.004	0.0216	0.00	-0.18	-0.18	-0.13%	0.00	
1600	24.600	24.599	24.599	0.001	0.0230	0.00	0.04	0.04	0.03%	0.00	
1700	26.956	26.951	26.951	0.005	0.0242	0.00	0.21	0.21	0.12%	0.00	
1800	29.423	29.417	29.417	0.006	0.0252	0.00	0.24	0.24	0.12%	0.00	
1900	31.992	31.988	31.988	0.004	0.0262	0.00	0.15	0.15	0.07%	0.00	
2000	34.666	34.669	34.669	-0.003	0.0272	0.00	-0.11	-0.11	-0.05%	0.00	
2100	37.442	37.441	37.441	0.001	0.0283	0.00	0.04	0.04	0.01%	0.00	
2200	40.329	40.321	40.321	0.008	0.0295	0.00	0.27	0.27	0.10%	0.00	
2300	43.340	43.331	43.331	0.009	0.0307	0.00	0.29	0.29	0.10%	0.00	

Comparison data stored in 'CMP0817.95'

Total change and %Change Radiance stored in 'P51_393.895'

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP9: Ribbon Filament Lamp Data Reduction

31. Use sample comparison data file below:

P51_393.895	08-17-1995	11:50:49 AM
Temp [°C]	dT [°C]	%dL
800	0.51	0.98%
900	0.00	0.00%
1000	-0.08	-0.11%
1100	-0.14	-0.16%
1200	-0.18	-0.18%
1300	-0.11	-0.10%
1400	-0.10	-0.08%
1500	-0.18	-0.13%
1600	0.04	0.03%
1700	0.21	0.12%
1800	0.24	0.12%
1900	0.15	0.07%
2000	-0.11	-0.05%
2100	0.04	0.01%
2200	0.27	0.10%
2300	0.29	0.10%

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP9: Ribbon Filament Lamp Data Reduction

32. Prepare calibration summary sheet
33. A sample calibration summary sheet is below:

**Radiance Temperature Lab
Calibration Summary Sheet**

Today's date: December 17, 1997

NIST Test No: 844/123456-95

Use sample report: 35050c.rfl

Company: Radiometric Systems Inc

Contact: John Doe

Phone: 301 123-4567

P.O.#: AB1234 dated January 1, 1995

Test number: 35050C Radiance Temperature Standard, Ribbon Filament Lamp

Test item: One General Electric Tungsten Ribbon Filament Lamp, Model # 30/T24/6,
Serial # P51

Temperatures measured: 800 °C to 2300 °C

Aligned at 1700 °C

Standard(s) used: SL20 (7.7788 A, 2.093 V), (Pcos(2θ)=0.001107, 12/92 TR2 p 123)

Temperature: 21 °C ± 1 °C

Relative humidity: 38% ± 5%

Calibration date: 5/24/95

Measurement date(s): 5/24/95 TR3 p42

Data taking program: C:\programs\datatak\pyro.mak\rtcal.frm 10/9/96

Calibration files: CAL0815.95, P510815.95

Data reduction program: C:\programs\lampredn\rtred.mak\lmpred.frm 5/25/95

Reduction summary: Red0525.95

Final data filename: P510895.TXT

Data analysis program: C:\programs\lampredn\rtred.mak\compare.frm 5/25/95

Reduction summary: Cmp0817.95

Previous calibration filename: P510393.TXT

Comparison filename: P51_393.895

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP9: Ribbon Filament Lamp Data Reduction

34. Enter data file name "XXXXmyy.MYY"
 - a. Prepare history file consisting of Microsoft EXCEL spreadsheet with plots of temperature difference from first calibration versus nominal temperature

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Variable Temperature Blackbody (BB) Alignment

Written by: Charles Gibson

Procedure: RTMLP10

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Variable Temperature Blackbody (BB) Alignment

CAUTION

COLLISION

1. Align PEP (see RTMLP2)
2. Remove field stop
3. Remove objective lens aperture, if any
4. Level aluminum plate ($<X$)
5. Use fiber optic light to illuminate BB
6. Translate BB along optical axis until center partition is in focus (Z)
7. Center BB on optical axis (X and Y)
8. Turn on BB (see RTMLP11)
9. Radiometrically center BB
 - a. Center $<Y$ about vertical axis
 - b. Center X along horizontal axis
 - c. Repeat i and ii until centered
 - d. Check centering
 - i. Monitor BB signal on graph logger
 - ii. Hold card in front of BB opening
 - iii. Check for equal drop-offs when blocking top, bottom, left, and right portions of BB opening

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Variable Temperature Blackbody (BB) Operation

Written by: Charles Gibson

Procedure: RTMLP11

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Variable Temperature Blackbody Operation


CAUTION
BRIGHT SOURCE
BURN
ELECTROCUTION

1. Turn On BB
 - a. Open tap water valve 1/8 turn
 - b. Check that blue box flow rate is about 2 gpm
 - c. Open chilled water valves
 - i. Open supply valve 1/4 turn
 - ii. Open return valve fully
 - d. Press reset button on blue water switch box to select chilled water
 - e. Check that blue box flow rate is about 5 gpm
 - f. Open argon cylinder valve fully
 - g. Turn on argon alarm and check that it trips at 10 cfh
 - h. Purge BB using argon flow rate of 70 cfh
 - i. Cover both BB openings for 5 s with your hands
 - ii. Uncover one end and wait 5 s
 - iii. Cover both BB openings for 5 s with your hands
 - iv. Uncover other end and wait 5 s
 - i. Repeat step h 5 times
 - j. Reduce argon flow rate to 10 cfh
 - k. Check that BB housing is cool and that argon is flowing
 - l. Insert correct detector filter for initial set point
 - i. Filter 0, 700 °C to 1300 °C, Low with insulation (LIN)
 - ii. Filter 1, 1300 °C to 1800 °C, Medium (MED)
 - iii. Filter 2, 1800 °C to 2500 °C, High (HI)
 - iv. Filter 3, 2500 °C to 3100 °C, Extra High(XH)
 - m. Switch control mode to [MANUAL]
 - n. Turn potentiometer to 0

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP11: Variable Temperature Blackbody (BB) Operation

- o. Press [START], [WATER ON], and [RUN] buttons
 - i. Note that BB does not control chilled water flow, because automatic BB shutoff switch has been disabled
- p. Turn potentiometer slowly until voltmeter reads between 150 VAC and 200 VAC (sets BB below 1300 °C)


2. Operate BB

- a. Run BB control program
 - i. Program is BBPROG.PRG in BetterBasic
 - ii. Enter BB at DOS prompt
 - iii. Enter measurement direction
 - (1) '0' to run high to low
 - (2) '1' to run low to high
 - iv. Enter set point unit
 - (1) 'C' for °C
 - (2) 'F' for °F
 - v.  Enter temperature range
 - (1) LOW
 - (2) LIN
 - (3) MED
 - (4) HI
 - (5) XH
 - vi. Insert correct filter and press [Enter]
- b. Press [Shift S] to enter initial set point
- c. Adjust potentiometer until temperature is 20 °C above set point
- d. Switch control mode to [PC] to maintain set point

WARNING! DO NOT INCREASE SET POINT BY STEPS LARGER THAN 20 °C TO REDUCE DAMAGE TO BB

- e. Select new set point (set point is in °C or °F depending on set point unit selected)
 - i. Press [Shift S] to enter set point
 - ii. Use numeric keypad to increment set point
 - (1) Press increment key
 - (a) Press 7 for 100 °C
 - (b) Press 4 for 10 °C
 - (c) Press 1 for 1 °C

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP11: Variable Temperature Blackbody (BB) Operation

- (d) Press 0 for 0.1 °C
 - (2) Press direction key to set new set point
 - (a) Press 8 to increase set point
 - (b) Press 2 to decrease set point
 - f. Pause program (suspending BB control) if new set point is not in current range
 - i. Use default set points depending on measurement direction
 - (1) Use high to low default set points
 - (a) Select XH for 2500 °C
 - (b) Select HI for 2300 °C
 - (c) Select MED for 1800 °C
 - (d) Select LOW for 1300 °C
 - (2) Use low to high default set points
 - (a) Select LOW for 800 °C
 - (b) Select MED for 1300 °C
 - (c) Select HI for 1800 °C
 - (d) Select XH for 2500 °C
 - ii. Switch control mode to [MANUAL] if new set point is not default value
 - iii. Loop program step 2.a.v. (see )
3. Use **EMERGENCY TURN OFF** procedure when necessary
 - a. Press [EMERGENCY STOP] button on power supply (chilled water and argon will continue flowing)
 - b. Continue with steps 4.b. to 4.j.
4. **Turn Off BB**
 - a. Switch control mode to [MANUAL]
 - b. Turn to potentiometer to 0
 - c. Press [EMERGENCY STOP] button on power supply (chilled water and argon will continue flowing)
 - d. Switch off main breaker on wall
 - e. Wait until furnace has stopped glowing
 - f. Close chilled water supply valve
 - g. Close chilled water return valve
 - h. Close cold water supply valve
 - i. Close argon cylinder valve
 - j. Close argon flow meter valve

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Disappearing Filament Optical Pyrometer Alignment

Written by: Charles Gibson

Procedure: RTMLP12

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Disappearing Filament Optical Pyrometer Alignment

1. Align BB (see RTMLP10)
2. Check that objective lens and eye piece are clean
3. Check batteries
4. Zero galvanometer
5. Use post with magnetic base, height adjustment, and angle bracket to mount pyrometer
6. Position post and level angle bracket mounting surface
7. Mount pyrometer to angle bracket

CAUTION

COLLISION

8. Use joystick to move BB to test pyrometer position
9. Use fiber optic light to illuminate BB
10. Set pyrometer at measuring distance from BB center partition
 - a. Use distance of 30.2 cm from center partition to end of graphite extension tube
 - b. Measure distance from front surface of objective lens holder to extension tube (for 1 m, distance is 69.8 cm)
11. Adjust pyrometer eye piece to focus on illuminated lamp filament
12. Focus pyrometer on center partition
13. Adjust pyrometer and/or move BB so that pyrometer views down center of BB cavity
14. Check distance
15. Adjust pyrometer and/or move BB so that pyrometer views down center of BB cavity
16. Record table position
17. Repeat for other test pyrometers

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Radiation Thermometer (RT) Alignment

Written by: Charles Gibson

Procedure: RTMLP13

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Radiation Thermometer (RT) Alignment

1. Align BB (see RTMLP10)
2. Check that objective lens and eye piece are clean
3. Use post with magnetic base, height adjustment, and angle bracket to mount RT
4. Position post and level angle bracket mounting surface
5. Mount RT to angle bracket

CAUTION

COLLISION

6. Use joystick to move BB to test RT position
7. Use fiber optic light to illuminate BB
8. Set RT at measuring distance from BB center partition
 - a. Use distance of 30.2 cm from center partition to end of graphite extension tube
 - b. Measure distance from front surface of objective lens holder to extension tube (for 1 m, distance is 69.8 cm)
 - c. Measure distance from alignment mark to extension tube for Land Cyclops RT (for 1 m, distance is 69.8 cm)
9. Adjust RT eye piece to operator's eye
10. Focus RT on center partition
11. Adjust RT and/or move BB so that RT views down center of BB cavity
12. Check distance
13. Adjust RT and/or move BB so that RT views down center of BB cavity
14. Record table position
15. Repeat steps 2 to 14 for each test RT
16. Turn on BB (see RTMLP11)
17. Check centering of each test RT
 - a. Monitor RT signal on graph logger or look at RT display panel
 - b. Hold card in front of BB opening
 - c. Check for equal drop-offs when blocking top, bottom, left, and right portions of BB opening

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Radiation Thermometer (Optical Pyrometer)
Calibration Procedure

Written by: Charles Gibson

Procedure: RTMLP14

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Radiation Thermometer (Optical Pyrometer) Calibration Procedure

1. Review calibration summary
 - a. Align standard lamp, BB, and test RTS
 - b. Measure from 2700 °C to 800 °C in 100 °C steps
 - c. Repeat 2300 °C, 1900 °C, 1500 °C, 1300 °C, 1100 °C, and 800 °C
2. Run general data taking program (see RTMLP1)
3. Align standard lamp (see RTMLP6)
4. Align BB (see RTMLP10)
5. Align RTS (see RTMLP13)
6. Run calibration program
 - a. Click on "Calibrations" from *Hardware.frm* menu to load *Rtcal.frm*
 - b. Select "Radiation thermometer" button
 - c. Enter data filename
 - d. Enter test filename suffix\$
 - i. Create test filenames from names on *Software.frm*
 - e. Select Yes and use suggested ranges (PMT high voltage and amplifier ranges are part of *Rtcal.frm* program)
 - f. Click on "Start" button
 - g. Enter data on *Software.frm*
 - i. Enter lamp and blackbody data
 - (1) Enter data for reference and BB
 - (2) Select "Use" button
 - (3) Enter lamp position
 - (4) Enter lamp name
 - (5) Select "Reference" button
 - (6) Enter lamp current (not necessary because lamp currents are automatically read from lamp monitoring computer)
 - (7) Enter power supply number

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP14: Radiation Thermometer (Optical Pyrometer)
Calibration Procedure

- ii. Enter radiation thermometer data
 - (1) Enter data for check standard and test RTS
 - (2) Select "Use" button
 - (3) Enter RT name
 - (4) Enter source table position
 - (5) Enter digital multimeter channel to read signal
 - (6) Enter RT type (O for optical pyrometers and R for radiation thermometers)
 - (7) Enter output type (P for panel or scale, D for digital, A for analog)
 - (8) Select instrument wavelength (nm)
- iii. Select temperature scale
- iv. Select measurement wavelength [nm]
- v. Enter number of repeats (number of test/standard measurements)
- vi. Enter number of readings (number of readings per DMM trigger)
- vii. Select lamp type
- viii. Select detector channel
- h. Click on "Save Configuration" button
- i. Select "Continue" button to return to *Rtcal.frm*
- j. Enter Standard temperature in °C in input box
- k. Enter Nominal radiance temperature in input box
- l. Read lamp currents automatically
- m. Enter BB set point
- n. Wait 5 min for BB to stabilize
- o. Move source table to reference source and measure full scale signal
- p. Move source table to dark signal position and measure dark signal
- q. Move source table to BB and measure full scale signal
- r. Move source table to dark signal position and measure dark signal on BB amplifier range if reference and BB amplifier ranges are different
- s. Repeat steps o to r "number of repeats" times
- t. Enter emissivity correction for test RT wavelength
 - i. Enter 0 for 655.3 nm
 - ii. Enter 1 for 900 nm
 - iii. Enter 2 for 1000 nm
 - iv. Enter 3 for true temperature

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP14: Radiation Thermometer (Optical Pyrometer)
Calibration Procedure

- u. Move BB in front of each RT and take reading
 - i. Read RTs with digital outputs 3 times
 - ii. Read RTs with analog outputs 3 times
 - iii. Enter display value for RTs with digital display panels
 - iv. Read OPs matching from bright to dark and then from dark to bright
 - (1) Enter scale value or press [Enter] to read output signal
 - (2) Read 2 sets
 - v. Select next action and then select "Continue" button
 - i. Repeat measurement with no changes
 - ii. Measure with new temperature
 - iii. Measure with new software changes
 - iv. Stop measurements
 - w. Repeat xxii for each measurement point
7. Store calibration data in log file "CALMMDD.YY"
8. Use sample log file on next page:

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP14: Radiation Thermometer (Optical Pyrometer)
Calibration Procedure

Date 12-13-1996 Time 9:48:10 AM

File Name = CAL1213.96

PMT Detector PMT HV = 449 V NRDS = 5 Integration Time = 8.35 s

Wavelength = 655.3 nm Nominal Rad Temp = 1800 °C

Name	SL20	TGBB
Rad. Temp	1254.96 °C	1799.01 °C
Current	7.7788 A	
Sp. Rad.	5.690E+02	2.472E+04
Ratio	1.000E+00	4.344E+01
Std. Ratio	0.06%	0.05%
Amp Range	9	8
Full Scale	1.016E+00 ±0.00%	4.411E+00 ±0.00%
	1.016E+00 ±0.00%	4.411E+00 ±0.01%
	1.015E+00 ±0.00%	4.410E+00 ±0.01%
Dark Sig.	1.510E-04	-1.166E-04
Mean	1.015E+00	4.411E+00
Std(Mean)	0.04%	0.01%
Range Cf.	9.9858E-02	1.0000E+00
SS Cf.	1.000E+00	9.987E-01
Corr. Mean	1.014E-01	4.405E+00
Position #	1	3
Table Pos.	37500	127500

Position CY52 TEST PYRO # 1

Reading # 1 1804 °C 0.0 °C

Mean = 1804.0 °C STD = 0.0 °C

Temperature Set for THERMOGAGE = 1800 °C

Measured Temperature at 655.3 nm = 1799.01 °C

Calculated Temperature at 655.3 nm = 1799.01 °C

Mean 1804.0 °C

Test Std 0.0 °C

Test Corr. -6.0 °C

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP14: Radiation Thermometer (Optical Pyrometer)
Calibration Procedure

9. Store RT data in test file "XXXXMMDD.YY"
10. Use sample test data file below:

```
c:\PYRO\CY521213.96
12-13-1996  9:59:06 AM
Standard used is SL20 at 1254.96 °C
Temperature in °C
Temperature measured at 655.3 nm Temperature calculated at 1000 °C
```

Tnom	Tbb	Tcorr	CY52	Std	Corr
1800	1799.01	1797.98	1804.0	0.0	-6.0
1700	1699.18	1698.26	1704.5	0.0	-6.2
1600	1599.47	1598.64	1603.0	0.0	-4.4
1500	1499.31	1498.56	1502.0	0.0	-3.4
1400	1399.24	1398.57	1401.0	0.0	-2.4
1300	1299.30	1298.71	1302.0	0.0	-3.3
1200	1199.36	1198.85	1201.0	0.0	-2.2
1100	1099.34	1098.89	1100.5	0.0	-1.6
1100	1099.25	1098.80	1100.5	0.0	-1.7
1000	999.15	998.76	1000.0	0.0	-1.2

11. Reduce data (see RTMLP15)

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Radiation Thermometer (RT) Data Reduction

Written by: Charles Gibson

Procedure: RTMLP15

Approved by: _____

Date: _____

Revised by: _____

Date: _____

Radiation Thermometer (RT) Data Reduction

1. Print out test file for RT check standard
 - a. Verify that RT check standard values are within uncertainty
2. Print out test file for BB RT check standard
 - a. Verify that BB RT check standard values are within uncertainty
3. Print out test files
4. Look in lab notebook and document any changes on print out and then edit file
5. Load EXCEL and open test file
6. Separate data by range or filter
7. Sort data
8. Calculate mean and standard deviation of repeats
9. Plot data
10. Save in data file "XXXXMMYY.XLS"
11. Use sample data file below:

TEST1096.XLS			11/7/96					
c:\PYRO\TEST1007.96								
Tnom	Tbb	Tcorr	TEST	Std	Corr	Corr Mean	Corr Std	Filter
°C	°C	°C	°C	°C	°C	°C	°C	
800	799.10	798.83	801	0	-2.2	-2.5	0.4	Clear
800	800.58	800.30	803	0	-2.7			
900	899.04	898.71	900.5	0	-1.8			
1000	999.13	998.74	1000	0	-1.3			
1100	1099.19	1098.74	1100	0	-1.3	-1.2	0.2	
1100	1100.48	1100.03	1101	0	-1.0			
1200	1199.32	1198.80	1200	0	-1.2			

12. Use, for repeats, values of blackbody temperature (T_{corr}), test RT reading (TEST), and test RT correction (Corr) that best represents mean of repeats so that actual reading is reported
13. Save data in text file (tab delimited) "XXXXMMYY.TXT"

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP15: Radiation Thermometer (RT) Data Reduction

14. Use sample file below:

Test1096.TXT			11/4/96	
Tnom	Tcorr	TEST	Corr	Filter
[°C]	[°C]	[°C]	[°C]	Clear
800	800.3	803	-2.7	
900	898.7	900.5	-1.8	
1000	998.7	1000	-1.3	
1100	1098.7	1100	-1.3	
1200	1198.8	1200	-1.2	

15. Prepare calibration summary sheet
16. Use sample calibration summary sheet on next page:

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP15: Radiation Thermometer (RT) Data Reduction

**Radiance Temperature Lab
Calibration Summary Sheet**

Today's date: December 17, 1997

NIST Test No: 844/123456-95

Use sample report: 35070sdc.1um

Company: Radiometric Systems Inc

Contact: Joe Calibration

Phone: 301 123-4567

P.O.#: AB1234 dated January 1, 1995

Test number: 35070S Special Test of Radiation Thermometer

Test item: One Minolta-Land Infrared Thermometer, Model # Cyclops 52, Serial # 0000000

Wavelength = 1000 nm

Distance: 100 cm from alignment mark

Target size: 5.6 mm

Focus: 1 meter Clear and dark lenses (sn L1)

Temperatures measured: 800 °C to 2500 °C

800 °C to 1200 °C with clear filter

1100 °C to 2500 °C with dark filter

Standard(s) used: SL20 (7.7788 A / 2.093 V), (Pcos(20)=0.001107, 12/92 TR2 p 123)

Land Cyclops 52 pyrometer, Serial no.: 202681

Thermogage 7000-1 pyrometer, Serial no.: 5127

Temperature: 21 °C ± 1 °C

Relative humidity: 38% ± 5%

Calibration date: 2/16/96

Measurement date(s): 2/15-16/96 TR3 p79

Data taking program: C:\programs\Lmpredn\pyro.mak\rtcal.frm 10/9/96

Calibration files: Cal1007.96, Test1007.96 Tgbb1007.96

Data reduction program: Excel 7.0

Data reduction file: Test1096.xls

Final data file: Test1096.txt

Optical Technology Division (844)
Radiance Temperature Measurements Laboratory Procedure
Procedure RTMLP15: Radiation Thermometer (RT) Data Reduction

17. Print out test file for RT check standard
 - a. Prepare history file consisting of Microsoft EXCEL spreadsheet with plots of temperature difference from first calibration versus nominal temperature

Appendix E: How to Request Radiance Temperature Measurement Services

1. Prepare a purchase order with the following:
 - a. Company name and address
 - b. Contact person's name and phone number.
 - c. Billing address
 - d. Service ID number requested (include range and points)
 - e. Test fee
 - f. Address to ship test item
 - g. Return shipping instructions (prepay and add, COD, charge to account with shipper)
 - i. the test fee does not include shipping costs
 - ii. customers are responsible for all shipping costs
 - iii. If no instructions are given, NIST will return the test item by common carrier, collect, and uninsured.
2. Send the purchase order to:

National Institute of Standards and Technology
Calibration Program
Building 820, Room 232
Gaithersburg, MD 20899-0001

Tel.: (301) 975-2002
FAX: (301) 869-3548
E-mail: calibrations@enh.nist.gov
3. Send the test item(s) to:

Charles Gibson
NIST
Building 221 / Room B208
Gaithersburg, MD 20899-0001

Tel.: (301) 975-2329
FAX: (301) 869-5700
E-mail: cgibson@nist.gov
WWW: <http://physics.nist.gov/Divisions/Div844/staff/Gp3/gibson.html>

Radiance temperature calibrations are performed in February, May, August, and November. Requests for calibration services are scheduled after the receipt of a purchase order.

NIST policy requires prepayment of all calibration services performed for non-U.S. organizations. Please contact the Calibration Program office to arrange payment.

THE SP 250 SERIES ON NIST MEASUREMENT SERVICES*

SP 250-1	Spectral Radiance Calibrations PB87179883	SP 250-21	Calibration of Beta-Particle Radiation Instrumentation PB88201579
SP 250-2	Far Ultraviolet Detector Standards PB87227609	SP 250-22	Platinum Resistance Thermometer Calibrations PB88138367
SP 250-3	Radiometric Standards in the Vacuum Ultraviolet PB87227625	SP 250-23	Liquid-in-Glass Thermometer Calibration Service PB89128888
SP 250-4	Fricke Dosimetry in High-Energy Electron Beams PB88110374	SP 250-24	Standard Cell Calibrations PB88123690
SP 250-5	Alpha-Particle Calibrations PB88168620	SP 250-25	Calibration Service for Inductive Voltage Dividers
SP 250-6	Regular Spectral Transmittance PB88108550	SP 250-26	NBS Phase Angle Calibration Services PB88225636
SP 250-7	Radiance Temperature Calibrations PB88123674	SP 250-27	AC-DC Difference Calibrations PB892222616
SP 250-8	Spectral Reflectance PB88109905	SP 250-28	Solid-State DC Voltage Standard Calibrations PB88168703
SP 250-9	Calibration of Beta-Particle-Emitting Ophthalmic Applicators PB88108535	SP 250-29	Traceable Frequency Calibrations PB88168364
SP 250-10	Radioactivity Calibrations with the "4 π " Gamma Ionization Chamber and Other Radioactivity Calibration Capabilities PB88123708	SP 250-30	GOES Satellite Time Code Dissemination: Description and Operation PB88168760
SP 250-11	Dosimetry for High Dose Applications PB88201587	SP 250-31	Mass Calibrations PB89153894
SP 250-12	Neutron Personnel Dosimetry PB87227617	SP 250-32	A Calibration Service for 30 MHz Attenuation and Phase Shift PB88238324
SP 250-13	Activation Foil Irradiation with Californium Fission Sources PB88217443	SP 250-33	A Calibration Service for Voltage Transformers and High-Voltage Capacitors PB882252903
SP 250-14	Activation Foil Irradiation by Reactor Cavity Fission Sources PB88217435	SP 250-34	High Vacuum Standard and Its Use PB89193841
SP 250-15	Photometric Calibrations PB88153747	SP 250-35	The Calibration of Thermocouples and Thermocouple Materials PB89209340
SP 250-16	Calibration of X-Ray and Gamma-Ray Measuring Instruments PB88211826	SP 250-36	A Calibration Service for Current Transformers PB91216770
SP 250-17	The NBS Photodetector Spectral Response Calibration Transfer Program PB88201595	SP 250-37	Photometric Calibrations PB97148472
SP 250-18	Neutron Source Strength Calibrations PB88211818	SP 250-38	NIST Leak Calibration Service PB92149772
SP 250-19	Calibration of Gamma-Ray-Emitting Brachytherapy Sources PB89193858	SP 250-39	NIST Pressure Calibration Service PB94164043
SP 250-20	Spectral Irradiance Calibrations PB88123781	SP 250-40	Absorbed-Dose Calibration of Ionization Chambers in a ⁶⁰ Co Gamma-Ray Beam SN003-003-03034-1 \$2.00

* Entries containing a stock number (SN003-003-) and price can be purchased from the Superintendent of Documents, U.S. Government Printing Office, Washington, DC 20402. GPO will accept checks, money orders, VISA, and MasterCard. For more information, or to place an order, call (202) 512-1800. Be sure to cite the stock number on all orders.

Entries containing PB numbers can be purchased from the National Technical Information Service, Springfield, VA 22161. NTIS will accept American Express in addition to the payment methods listed for GPO. For more information call (703)487-4650; to place an order call (800) 553-6487. Fax: (703) 321-8547. Be sure to cite the PB number on all orders.

Entries without stock or PB numbers are in preparation.

U.S. Department of Commerce
National Institute of Standards
and Technology
Gaithersburg, MD 20899-0001

Official Business
Penalty for Private Use \$300