







# NBS TECHNICAL NOTE 687

U.S. DEPARTMENT OF COMMERCE / National Bureau of Standards

## Measurement of RF Power-Absorption in Biological Specimens (10 to 100 MHz)

QC  
100  
U5753  
No. 687  
1976  
C.2

## NATIONAL BUREAU OF STANDARDS

The National Bureau of Standards<sup>1</sup> was established by an act of Congress March 3, 1901. The Bureau's overall goal is to strengthen and advance the Nation's science and technology and facilitate their effective application for public benefit. To this end, the Bureau conducts research and provides: (1) a basis for the Nation's physical measurement system, (2) scientific and technological services for industry and government, (3) a technical basis for equity in trade, and (4) technical services to promote public safety. The Bureau consists of the Institute for Basic Standards, the Institute for Materials Research, the Institute for Applied Technology, the Institute for Computer Sciences and Technology, and the Office for Information Programs.

**THE INSTITUTE FOR BASIC STANDARDS** provides the central basis within the United States of a complete and consistent system of physical measurement; coordinates that system with measurement systems of other nations; and furnishes essential services leading to accurate and uniform physical measurements throughout the Nation's scientific community, industry, and commerce. The Institute consists of the Office of Measurement Services, the Office of Radiation Measurement and the following Center and divisions:

Applied Mathematics — Electricity — Mechanics — Heat — Optical Physics — Center for Radiation Research: Nuclear Sciences; Applied Radiation — Laboratory Astrophysics<sup>2</sup> — Cryogenics<sup>2</sup> — Electromagnetics<sup>2</sup> — Time and Frequency<sup>2</sup>.

**THE INSTITUTE FOR MATERIALS RESEARCH** conducts materials research leading to improved methods of measurement, standards, and data on the properties of well-characterized materials needed by industry, commerce, educational institutions, and Government; provides advisory and research services to other Government agencies; and develops, produces, and distributes standard reference materials. The Institute consists of the Office of Standard Reference Materials, the Office of Air and Water Measurement, and the following divisions:

Analytical Chemistry — Polymers — Metallurgy — Inorganic Materials — Reactor Radiation — Physical Chemistry.

**THE INSTITUTE FOR APPLIED TECHNOLOGY** provides technical services to promote the use of available technology and to facilitate technological innovation in industry and Government; cooperates with public and private organizations leading to the development of technological standards (including mandatory safety standards), codes and methods of test; and provides technical advice and services to Government agencies upon request. The Institute consists of the following divisions and Centers:

Standards Application and Analysis — Electronic Technology — Center for Consumer Product Technology: Product Systems Analysis; Product Engineering — Center for Building Technology: Structures, Materials, and Life Safety; Building Environment; Technical Evaluation and Application — Center for Fire Research: Fire Science; Fire Safety Engineering.

**THE INSTITUTE FOR COMPUTER SCIENCES AND TECHNOLOGY** conducts research and provides technical services designed to aid Government agencies in improving cost effectiveness in the conduct of their programs through the selection, acquisition, and effective utilization of automatic data processing equipment; and serves as the principal focus within the executive branch for the development of Federal standards for automatic data processing equipment, techniques, and computer languages. The Institute consists of the following divisions:

Computer Services — Systems and Software — Computer Systems Engineering — Information Technology.

**THE OFFICE FOR INFORMATION PROGRAMS** promotes optimum dissemination and accessibility of scientific information generated within NBS and other agencies of the Federal Government; promotes the development of the National Standard Reference Data System and a system of information analysis centers dealing with the broader aspects of the National Measurement System; provides appropriate services to ensure that the NBS staff has optimum accessibility to the scientific information of the world. The Office consists of the following organizational units:

Office of Standard Reference Data — Office of Information Activities — Office of Technical Publications — Library — Office of International Relations — Office of International Standards.

<sup>1</sup> Headquarters and Laboratories at Gaithersburg, Maryland, unless otherwise noted; mailing address Washington, D.C. 20234.

<sup>2</sup> Located at Boulder, Colorado 80302.

3 1977  
acc. - ref.  
753  
687  
6  
5

# Measurement of RF Power - Absorption in Biological Specimens (10 to 100 MHz)

technical note, no. 687

Frank M. Greene

Electromagnetics Division  
Institute for Basic Standards  
National Bureau of Standards  
Boulder, Colorado 80302



---

U.S. DEPARTMENT OF COMMERCE, Elliot L. Richardson, Secretary

Edward O. Vetter, Under Secretary

Dr. Betsy Ancker-Johnson, Assistant Secretary for Science and Technology

U.S. NATIONAL BUREAU OF STANDARDS, Ernest Ambler, Acting Director

Issued November 1976

NATIONAL BUREAU OF STANDARDS TECHNICAL NOTE 687  
Nat. Bur. Stand. (U.S.), Tech Note 687, 29 pages (November 1976)

CODEN: NBTNAE

U.S. GOVERNMENT PRINTING OFFICE  
WASHINGTON: 1976

---

For sale by the Superintendent of Documents, U.S. Government Printing Office, Washington, D.C. 20402  
(Order by SD Catalog No. C13.46:687) Price 75 cents. (Add 25 percent additional for other than U.S. mailing).

## CONTENTS

	<u>Page</u>
1. INTRODUCTION-----	2
2. THEORY-----	3
2.1 Introduction-----	3
2.2 Forward and Reflected Power-----	4
2.3 Analysis of the Field Generator-----	5
3. MEASUREMENT APPLICATIONS-----	10
3.1 Key Equations and Curves-----	10
3.2 Measurement Example-----	10
4. OPERATING PROCEDURE-----	14
5. MEASUREMENT SENSITIVITY-----	15
5.1 Basic Sensitivity-----	15
5.2 Increasing the Sensitivity-----	15
6. MAGNETIC-INDUCTION HEATING IN A LOSSY DIELECTRIC-----	19
6.1 Introduction-----	19
6.2 Theoretical Background-----	19
6.3 RF Heating Density-----	21
6.4 Total RF Heating-----	21
7. THE CALORIMETER METHOD-----	23
8. CONCLUSIONS-----	24
9. REFERENCES-----	25







MEASUREMENT OF RF POWER-ABSORPTION IN BIOLOGICAL SPECIMENS

(10 to 100 MHz)

Frank M. Greene  
Electromagnetics Division  
National Bureau of Standards  
Boulder, Colorado 80302

ABSTRACT

A method is described for accurately determining the rf power being absorbed by a biological specimen during non-ionizing radiation-exposure testing using the NBS RF Near-Field Synthesizer in the frequency range 10 to 100 MHz. This method is based solely on measuring the forward and reflected power on the transmission line feeding the synthesizer. Commercially available rf wattmeters can be used, and an automatic data-acquisition system employed, if desired, to "read" the meters and rapidly calculate, display, and record the rf power flow.

The method has the advantage that the exact measuring point on the feed line is not critical, as it is with methods employing direct impedance measurements, and that the required measurements can be made without interfering with the exposure tests.

Key words: Electromagnetic field hazards; electromagnetic field synthesizer; electromagnetic radiation-exposure testing (non-ionizing); near-fields; rf biological hazards.

## 1. INTRODUCTION

A method will be described for determining the total amount of rf power being absorbed by a biological specimen during non-ionizing radiation-exposure testing using the NBS RF Near-Field Synthesizer in the frequency range 10 to 100 MHz [1]. This method is based solely on measuring the forward and reflected power on the transmission line feeding the field synthesizer. Commercially available rf wattmeters and an automatic data-acquisition system can be employed to further simplify the procedure.\*

The near-field synthesizer is installed and operated in a doubly-shielded copper room having a shielding efficiency of 60 dB or more. Any electromagnetic radiation from the synthesizer or test specimen is reduced by a factor of at least  $10^6$  by the shielding. Radiation losses are therefore considered to be completely negligible. Coupling between the synthesizer and the interior copper walls of the shielded room has been found to be essentially reactive. Very little resistance is coupled back into the synthesizer from the walls of the room.

The rationalized mks International System of Units (SI) is used throughout this report.

---

\*This report is the outgrowth of work originally sponsored by the USAF School of Aerospace Medicine at Brooks AF Base, and the National Institute for Occupational Safety and Health of the Department of HEW.

## 2. THEORY

2.1 Introduction. The NBS RF Near-Field Synthesizer consists of both an electric-field, and a magnetic-field generator. These two field generators consist of a balanced parallel-plate stripline and a balanced four-gap single-turn loop inductor for generating the "wanted" electric and magnetic fields, respectively [1]. Multiple feeds and balanced symmetry are used to minimize the "unwanted" magnetic and electric field components and to permit superposing the "wanted" fields without significant cross coupling. For the sake of producing hazard-level fields with a minimum of driving power, these field generators are resonated. The two generators can be operated simultaneously to produce a large range of values of E/H, or individually [1].

Each field generator has its own driving transmitter (source) and rf feed system, and the two generators can be considered to be operating essentially independently. The analysis that follows applies to each field generator individually. If both field generators are operating simultaneously, the total power absorbed by the biological specimen may be obtained by adding the power absorbed from each field generator. It has been found at NBS and elsewhere [2] that the magnetic field may represent a greater biological hazard than the electric field in the frequency range under consideration here. Therefore, this analysis will probably be used for the most part with the magnetic-field-generator portion of the near-field synthesizer.

The use of steady-state, sinusoidally time-varying rf driving sources is assumed throughout this report. It is not necessary that these sources be matched to their respective feed lines with this method since only the ratio of the reflected power to the forward power is involved. The sources comprise tuned Class-C rf power amplifiers which must be operated with as much as a 5:1 impedance mismatch in order to obtain the required rf power output with a plate efficiency of 75 to 80 percent [3].

2.2 Forward and Reflected Power. Either field generator is fed rf power by means of a reasonably low-loss 50-ohm transmission line. Either field generator can be tuned and adjusted independently to present a 50-ohm load to its feed line at the operating frequency before the biological specimen is put in place. After the specimen is inserted, the rf power that it absorbs can be determined directly in terms of the reflection that it produces.

The quantities to be discussed now are to be measured on the transmission line feeding either field generator. Since only measurements of the component forward and reflected waves are involved, the exact point at which they are measured is not critical on a low-loss line. This is one of the advantages of this method over methods based on direct impedance measurement.

The complex reflection coefficient,  $\Gamma$ , at the receiving end of a uniform lossless transmission line is given by [4,5]

$$\Gamma = \frac{Z_L - Z_O}{Z_L + Z_O} = \frac{Z'_L - 1}{Z'_L + 1}, \quad (1)$$

where  $Z'_L = Z_L/Z_O$ .  $Z_L$  is the load impedance presented by either field generator to its feed line, and  $Z_O$  is the characteristic impedance of the line. For this analysis  $Z_L$  and  $Z_O$  are both considered to be resistive as will be discussed later, making  $\Gamma$  real in the limit.

The forward power,  $P_f$ , flowing toward the load is given by

$$P_f = \frac{V_f^2}{Z_O}, \quad (2)$$

where  $V_f$  is the rms value of the magnitude of the forward-traveling voltage wave.

The reflected power,  $P_r$ , flowing back to the source is given by

$$P_r = \frac{V_r^2}{Z_O}, \quad (3)$$

where  $V_r$  is the rms value of the magnitude of the backward-traveling (reflected) voltage wave. Also

$$V_r = \Gamma V_f. \quad (4)$$

Therefore, the reflected power is also given by

$$P_r = \frac{\Gamma^2 V_f^2}{Z_0}. \quad (5)$$

The ratio of the reflected power to the forward power is therefore given by

$$\frac{P_r}{P_f} = \Gamma^2. \quad (6)$$

The load power,  $P_L$ , actually absorbed by the load,  $Z_L$ , is the arithmetic difference between the forward power,  $P_f$ , and the reflected power,  $P_r$ , if  $Z_0$  is resistive, i.e.,

$$P_L = P_f - P_r. \quad (7)$$

The ratio of the load power to the forward power is given by

$$\frac{P_L}{P_f} = 1 - \Gamma^2, \quad (8)$$

as obtained from eqs. (6) and (7).

2.3 Analysis of the Field Generator. Each field generator can be represented by an equivalent, tapped, parallel-resonant tank circuit, as shown in figure 1. It should be pointed out, however, that the basic results obtained from this model are also valid near resonance when used with transmission-line types, as well as waveguide types of resonators [6], as will be discussed in connection with eq. (11).



From Ohm's law it can be shown, with the aid of figure 1, that the power absorbed by the biological specimen,  $P_a$ , and the load power,  $P_L$ , are in the same ratio as the internal rf loss resistances involved. That is,

$$\frac{P_a}{P_L} = \frac{R_2}{R_1 + R_2} , \quad (9)$$

where  $R_1$  represents the loss resistance of the field generator itself prior to inserting the biological specimen, and  $R_2$  is the loss resistance coupled into the generator from the biological specimen after retuning.

It is assumed that the presence of the specimen and the subsequent retuning do not cause a sufficient change in the field distribution within the field generator to change its internal metallic losses. This is believed to be a reasonable assumption. The metallic conductivities involved are so high that the form of the fields is unchanged by the losses that result. The total loss can therefore be computed by adding directly the values of the metallic losses and the dielectric losses [7].

The internal resistance of the field generator is transformed at its input terminals at resonance in accordance with the following relationship, valid provided  $Q \geq 10$  [8]

$$Z_L \cong \frac{X_1^2}{R} , \quad (10)$$

where, from figure 1,  $X_1$  is the inductive reactance of the shunt matching inductor, and  $R$  is the total series rf resistance of the tank circuit. When eq. (10) is substituted into eq. (9) we obtain, after rearranging,

$$\frac{P_a}{P_L} = 1 - \frac{R_1}{R_1 + R_2} = 1 - \frac{Z_L}{Z_0} = 1 - Z'_L , \quad (11)$$

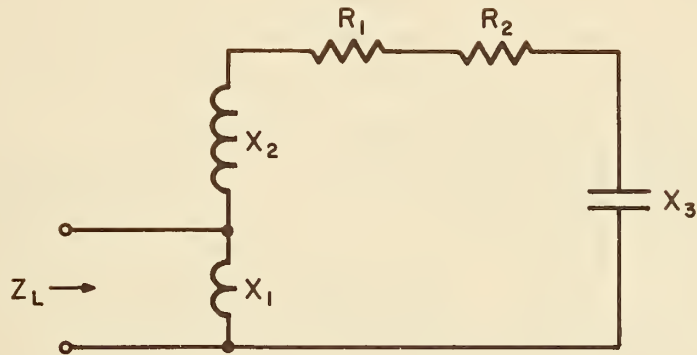


Figure 1. The equivalent lumped-constant circuit diagram of either the electric-field or the magnetic-field generator.

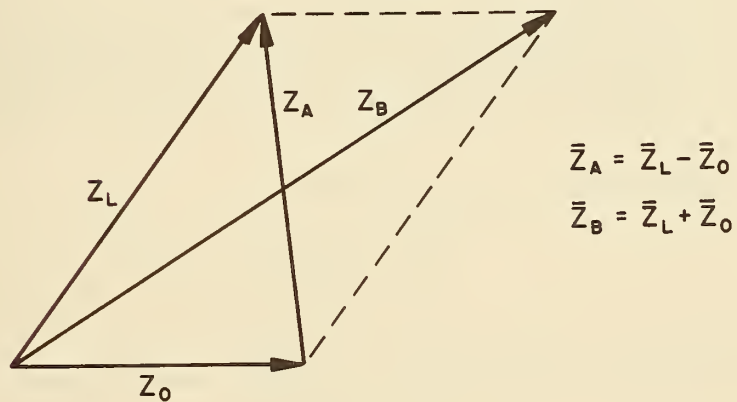


Figure 2. Vector diagram representing the numerator,  $\bar{Z}_L - \bar{Z}_O$ , and the denominator,  $\bar{Z}_L + \bar{Z}_O$ , of eq. (1a), where  $Z_O$  is considered to be resistive.



since  $Z_L = X_1^2 / (R_1 + R_2)$ , and  $Z_O = X_1^2 / R_1 = 50$  ohms, from eq. (10), with  $Z_L'$  limited to the range  $0 \leq Z_L' \leq 1$ . Equation (11) shows the basic relationship between the power absorbed by the biological specimen,  $P_a$ , and the input impedance of the field generator at resonance,  $Z_L$ , and can be used with any type of two-terminal resonator or cavity, as previously mentioned.

In practice, the field generator is first tuned and adjusted for zero reflected power (by adjusting both  $X_1$  and  $X_3$  of figure 1) to give an input impedance of  $Z_L = 50$  ohms before the biological specimen is put in place. Afterwards, the field generator is retuned for minimum reflected power (by adjusting  $X_3$  only). This cancels any reactance coupled in from the specimen, making  $Z_L$  resistive in both cases as shown in the following paragraph.

The complex receiving-end reflection coefficient is given by eq. (1), rewritten as

$$\bar{\Gamma} = \frac{\bar{Z}_L - \bar{Z}_O}{\bar{Z}_L + \bar{Z}_O} = \frac{\bar{Z}_A}{\bar{Z}_B} . \quad (1a)$$

The numerator is the vector difference and the denominator the vector sum of  $\bar{Z}_L$  and  $\bar{Z}_O$ , as shown in figure 2, where  $\bar{Z}_O$  is considered to be resistive. The numerator of eq. (1a) will be a minimum and the denominator a maximum simultaneously when  $\bar{Z}_L$  becomes resistive, thus minimizing  $\bar{\Gamma}$  as can be readily seen from figure 2.

Solving eq. (1) for  $Z_L'$  gives

$$Z_L' = \frac{1 + \Gamma}{1 - \Gamma} . \quad (12)$$

Substituting eq. (12) into (11) gives

$$\frac{P_a}{P_L} = \frac{2|\Gamma|}{1 + |\Gamma|} , \quad (13)$$

with  $\Gamma$  limited to the range  $-1 \leq \Gamma \leq 0$ . Using eq. (6), eq. (13) can also be written

$$\frac{P_a}{P_L} = \frac{2\sqrt{P_r/P_f}}{1 + \sqrt{P_r/P_f}}. \quad (14)$$

To obtain the ratio of the power absorbed by the specimen,  $P_a$ , to the forward power,  $P_f$ , eq. (13) must be multiplied by eq. (8), giving

$$\frac{P_a}{P_f} = 2|\Gamma| [1 - |\Gamma|]. \quad (15)$$

Using eq. (6) again, eq. (15) can also be written

$$\frac{P_a}{P_f} = 2\sqrt{P_r/P_f} [1 - \sqrt{P_r/P_f}]. \quad (16)$$

Finally, eq. (8) can be included in the form

$$\frac{P_L}{P_f} = 1 - \frac{P_r}{P_f}. \quad (17)$$

### 3. MEASUREMENT APPLICATIONS

3.1 Key Equations and Curves. Equations (14), (16), and (17) are the key equations to be used in this technique. As can be seen, the power absorbed by the biological specimen,  $P_a$ , can be obtained relative to either the load power,  $P_L$ , or the forward power,  $P_f$ , knowing only the ratio of the reflected power to the forward power,  $P_r/P_f$ , as summarized in table I.

Equations (14) and (17) are shown plotted in figure 3 as curves (a) and (b), respectively. Equations (16) and (17) are shown plotted in figure 4 as curves (a) and (b), respectively. It is interesting to note from curve (a) of figure 4 that the power absorbed by the specimen,  $P_a$ , can never exceed half of the forward power,  $P_f$ . This maximum occurs when the reflected power is 25 percent of the forward power.

3.2 Measurement Example. An example will be given to illustrate the use of the equations and curves. It is assumed that the field generator is first tuned and adjusted for zero reflected power to give an input impedance of 50 ohms before the biological specimen is put in place. Afterwards, the generator is retuned for minimum reflected power to cancel any reactance coupled in from the specimen, as previously discussed.

#### Example

Assume that the forward power,  $P_f = 100$  watts, and the reflected power,  $P_r = 11.1$  watts. Then  $P_r/P_f = 0.111$ . The load power,  $P_L$ , can be obtained from eq. (17) or curve (b) of figure 3, if desired, as  $P_L = 100 - 11.1 = 88.9$  watts. From eq. (14), or curve (a) of figure 3 (corresponding to  $P_r/P_f = 0.111$ ) the ratio  $P_a/P_L = 0.500$ . That is, half of the load power, or 44.4 watts, is absorbed by the biological specimen, and the other half goes into heating the field generator. This same value of absorbed power,  $P_a$ , can also be obtained directly from eq. (16), or curve (a) of figure 4, giving  $P_a/P_f = 0.444$ , or  $P_a = 44.4$  watts as before.

Table I

RF Power Absorbed by a Biological Specimen

<u>Reflected Power</u>	<u>Load Power</u>	<u>Power Absorbed by Specimen</u>	<u>Power Absorbed by Specimen</u>
Eq. (6)	Eq. (8)	Eq. (14)	Eq. (16)
$\frac{P_r}{P_f} \times 100$	$\frac{P_L}{P_f} \times 100$	$\frac{P_a}{P_L} \times 100$	$\frac{P_a}{P_f} \times 100$
(%)	(%)	(%)	(%)
0.00	100.00	0.00	0.00
0.10	99.90	6.13	6.12
0.20	99.80	8.56	8.54
0.50	99.50	13.21	13.14
1.00	99.00	18.18	18.00
2.00	98.00	24.78	24.28
4.00	96.00	33.33	32.00
6.00	94.00	39.35	36.99
8.00	92.00	44.10	40.57
10.00	90.00	48.05	43.25
11.11	88.89	50.00	44.44
15.00	85.00	55.83	47.46
20.00	80.00	61.80	49.44
25.00	75.00	66.67	50.00
30.00	70.00	70.78	49.54
40.00	60.00	77.49	46.49
50.00	50.00	82.84	41.42
60.00	40.00	87.30	34.92
70.00	30.00	91.11	27.33
80.00	20.00	94.43	18.89
90.00	10.00	97.37	9.74
100.00	0.00	100.00 (a)	0.00

(a) Theoretical limiting value.



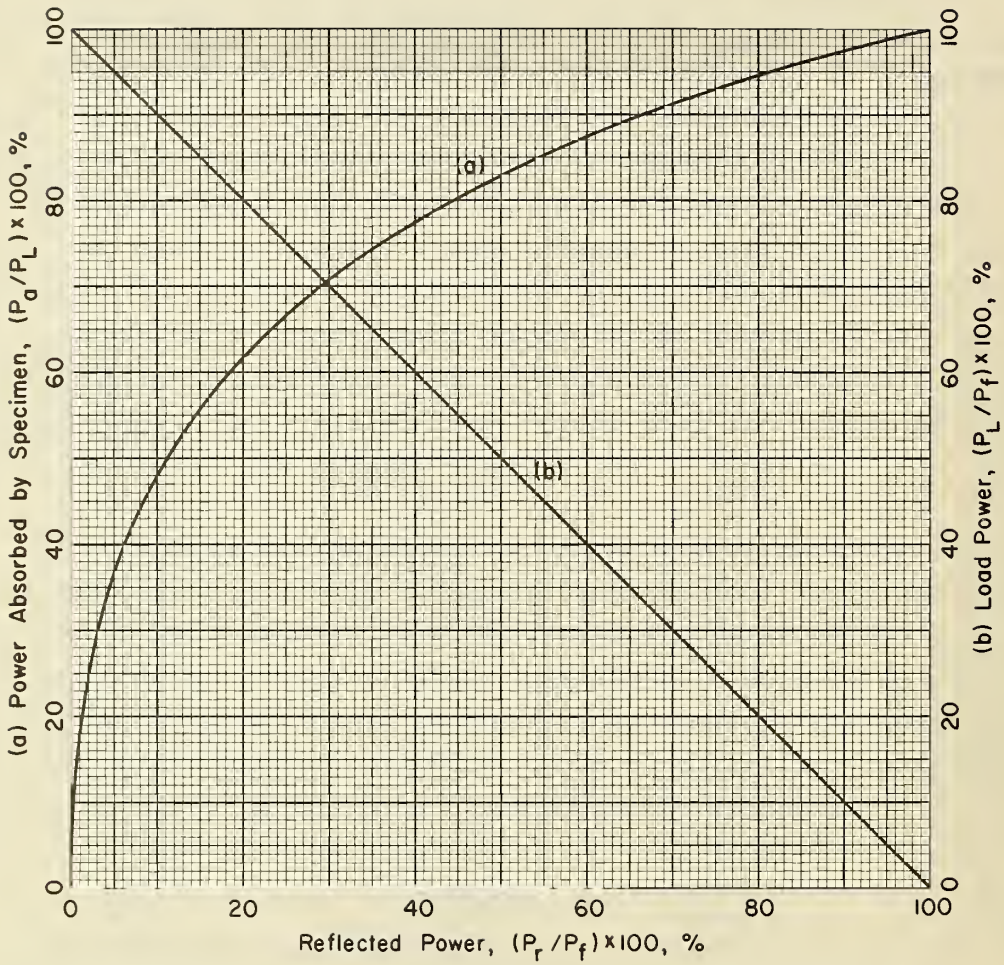


Figure 3. (a) RF power absorbed by the biological specimen,  $(P_a/P_L) \times 100, \%$ ; and (b) load power,  $(P_L/P_f) \times 100, \%$ ; both vs. the reflected power,  $(P_r/P_f) \times 100, \%$ ; where  $P_f$  is the forward power in watts.

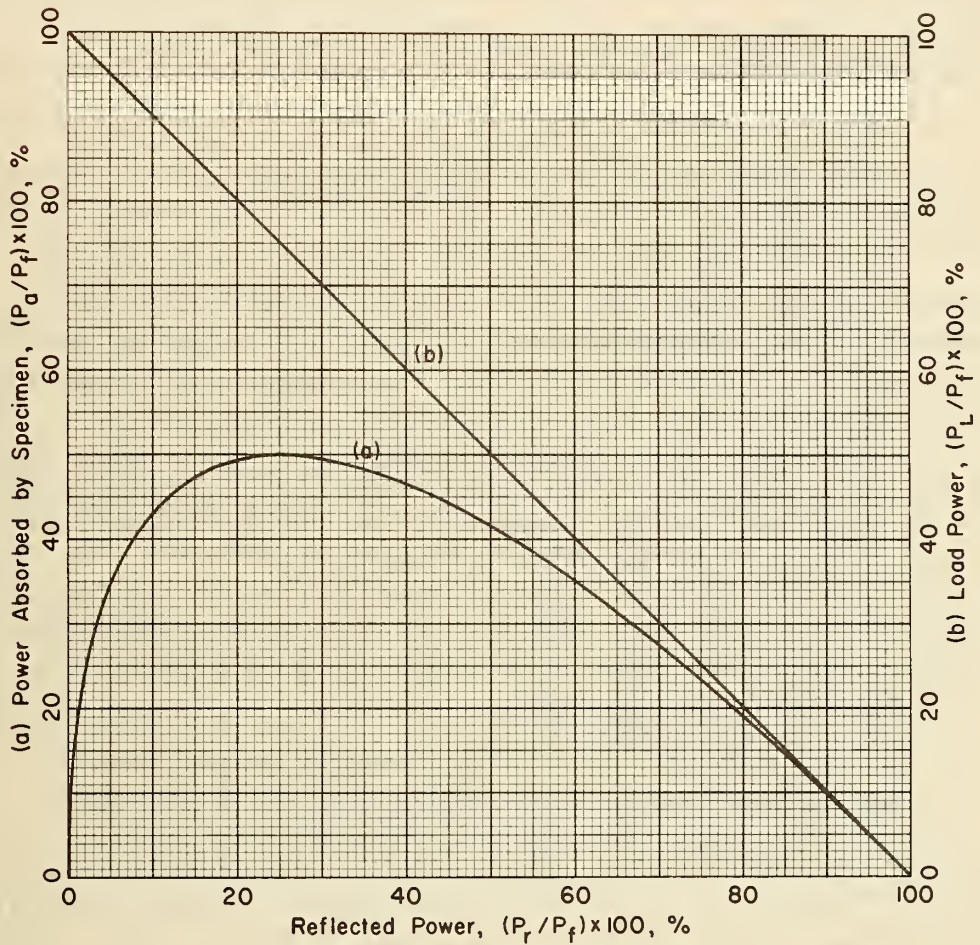


Figure 4. (a) RF power absorbed by the biological specimen,  $(P_a/P_f) \times 100, \%$ ; and (b) load power,  $(P_L/P_f) \times 100, \%$ ; both vs. the reflected power,  $(P_r/P_f) \times 100, \%$ ; where  $P_f$  is the forward power in watts.

#### 4. OPERATING PROCEDURE

- (a) The rf directional couplers of the wattmeters being used should be inserted in the feed-line to each field generator, usually within one or two meters of the generator input terminals. The exact location is, however, not critical with this method, as discussed in section 2.2.
- (b) The input channels of the automatic data-acquisition system can be connected directly to the respective dc terminals of the microammeters used to indicate forward and reflected power. Store the calibration data of these wattmeters in the associated computer as part of the data acquisition program (i.e., rf watts input to the coupler vs. dc millivolts output).
- (c) The field generator should be adjusted at the desired operating frequency (in the absence of the biological specimen) in accordance with the usual procedure [1].
- (d) Adjust the generator matching networks for zero reflected power, so that the generator presents a load impedance  $Z_L = 50$  ohms to its feed line.
- (e) The rf power input to the field generator should then be reduced and the biological specimen placed at its center.
- (f) The rf power should then be gradually increased and the tuning capacitors of the generator matching units carefully readjusted for minimum reflected power. This cancels the reactance coupled into the field generator from the specimen as discussed in section 2.3
- (g) Read the forward power,  $P_f$ , and the reflected power,  $P_r$ . Calculate the ratio,  $P_r/P_f$ , and the load power,  $P_L = P_f - P_r$ . The remainder of the calculations should be made as outlined in section 3.2.



## 5. MEASUREMENT SENSITIVITY

5.1 Basic Sensitivity. The basic measurement sensitivity [9] or resolution,  $\delta$ , of the method being described is given by the ratio of the reflected power,  $P_r/P_f$ , to the power absorbed by the biological specimen,  $P_a/P_f$ , at the final operating point on the characteristic curve (eq. (16)), i.e.,

$$\text{Measurement sensitivity} = \delta = \frac{P_r/P_f}{P_a/P_f} . \quad (18)$$

From eq. (16) or curve (a) of figure 4, it can be seen that the measurement sensitivity is zero when  $P_a/P_f = 0$ , and increases with increasing values of  $P_a/P_f$  to a sensitivity of 0.500 when  $P_a/P_f = 0.500$ , probably the upper limit of the useful measurement range. While the sensitivity is deemed to be adequate over most of the measuring range, it can be increased for low values of  $P_a/P_f$  by: (a) using a more sensitive wattmeter element for measuring the reflected power (readily available); (b) reducing the power and increasing the coupling between the specimen and the field generator to increase  $P_r/P_f$ ; or (c) choosing the initial operating point higher on the characteristic curve.

5.2 Increasing the Sensitivity. The latter method of increasing the measurement sensitivity involves adjusting the input impedance of the field generator,  $Z_L$ , to an initial value greater than 50 ohms prior to placing the biological specimen in the generator. The resulting relationship between  $P_r/P_f$  and  $P_a/P_f$  can be determined by rewriting eq. (11) as

$$\frac{P_a}{P_L} = 1 - \frac{Z_L'}{N} , \quad (19)$$

where  $N = Z_L/50$ .  $Z_L'$  is limited to real values in the range  $0 \leq Z_L' \leq N$ . Equation (12) is then substituted into eq. (19) and

the result multiplied by eq. (8), giving

$$\frac{P_a}{P_f} = \frac{(N-1) - 2\Gamma - (N+1)\Gamma^2}{N}, \quad (20)$$

where  $\Gamma$  is limited to real values in the range  $-1 \leq \Gamma \leq \frac{N-1}{N+1}$ . Using eq. (6), eq. (20) can be written

$$\frac{P_a}{P_f} = \frac{(N-1) \pm 2\sqrt{P_r/P_f} - (N+1)(P_r/P_f)}{N}, \quad (21)$$

where the (+) sign is used when  $\Gamma$  is negative ( $Z_L < 50 \Omega$ ), and the (-) sign when  $\Gamma$  is positive ( $Z_L > 50 \Omega$ ).

The rf power meters used in this technique, of course, cannot distinguish between positive and negative values of  $\Gamma$ , so that some other means must be employed. One such simple technique would be to slowly increase the value of  $Z_L$ . If this causes  $P_r$  to increase, then  $Z_L > 50 \Omega$ . If an increase in  $Z_L$  causes a decrease in  $P_r$ , then  $Z_L < 50 \Omega$ .

Equation (21) is shown plotted in figure 5 for values of  $N = 1, 2, 3, 5, 10$ , and  $\infty$ , for the permissible range of values of  $\Gamma$  given for eq. (20). The upper part of each curve results from negative values of  $\Gamma$ , and the lower part from positive values. As can be seen, the measurement sensitivity (when  $P_a/P_f = 0$ ) has a value of  $\delta = 0$  when  $N = 1$ , increasing to a limiting value of  $\delta = 1$  as  $N$  is increased. This represents an improvement in sensitivity for values of  $P_a/P_f < 0.5$ .

As can be seen from figure 5, the reflected power (caused by the initial mismatch) is given by the intercept of each curve with the  $P_a/P_f = 0$  axis. These values of  $P_r/P_f$  can be determined from the expression,

$$\frac{P_r}{P_f} = \left[ \frac{N-1}{N+1} \right]^2, \quad (22)$$

and are 0.000, 0.111, 0.250, 0.444, 0.669, and 1.000, respectively, for the values of  $N = 1, 2, 3, 5, 10$ , and  $\infty$ . The measurement

sensitivity can be determined approximately from the inverse slope of each curve at its intercept point. As the biological specimen is introduced into the field generator,  $Z_L$  will decrease in value, improving the match. This technique should probably be investigated further and may require modification of the present tuning procedure used with the near-field synthesizer.

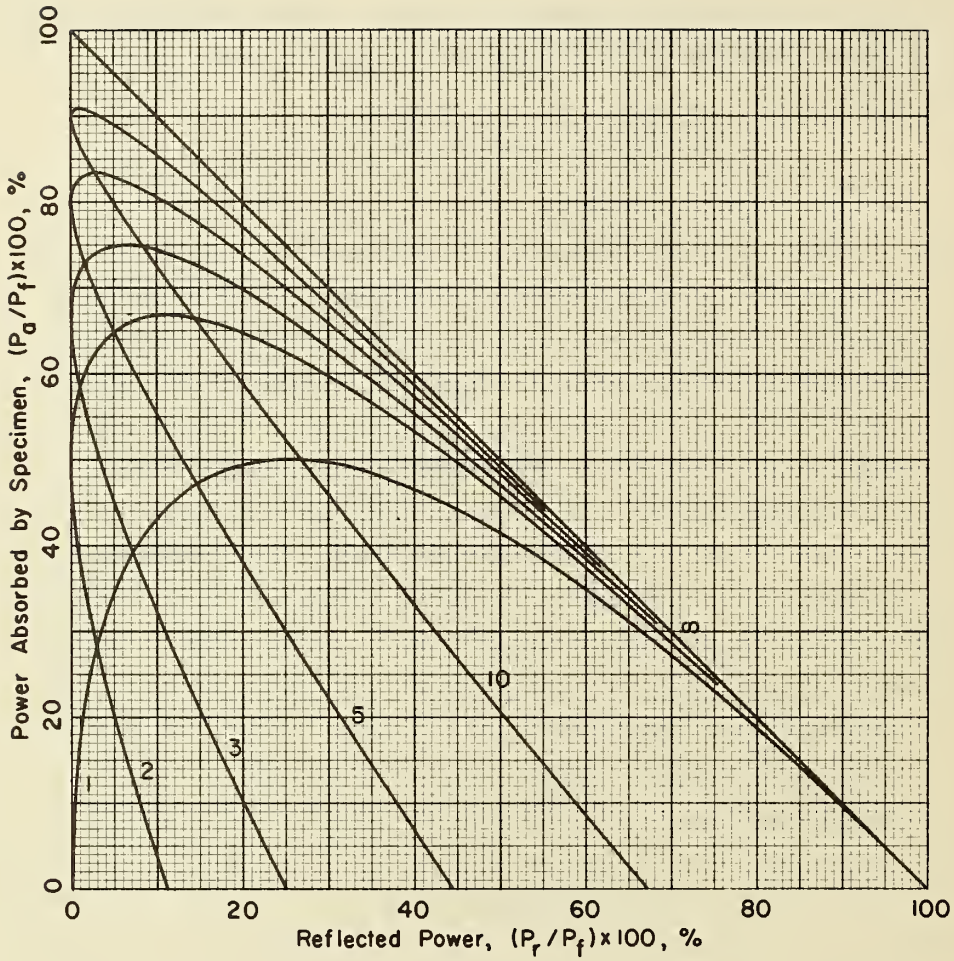


Figure 5. RF power absorbed by the biological specimen,  $(P_a/P_f) \times 100, \%$ , vs. the reflected power,  $(P_r/P_f) \times 100, \%$ , for values of  $N = 1, 2, 3, 5, 10$ , and  $\infty$ , where  $N = Z_L/50$  from eq. (19).



## 6. MAGNETIC-INDUCTION HEATING IN A LOSSY DIELECTRIC

6.1 Introduction. The quasi-static "magnetic-induction" heating caused by an essentially uniform rf magnetic field,  $H$ , in a homogeneous, lossy dielectric of permittivity,  $\epsilon$ , and conductivity,  $\sigma$ , will be determined in this section. The dielectric specimen is in the form of a right circular cylinder of radius,  $r$ , and height,  $h$ , as shown in figure 6, with  $H$  parallel to the axis of the cylinder. The dimensions involved in the dielectric are assumed to be small compared with the operating wavelength. As pointed out previously, the heating caused by the magnetic component of the field may often be greater than that caused by the electric component in the frequency range under consideration here [2].

6.2 Theoretical Background. The electric field,  $E$ , and the magnetic field,  $H$ , within the dielectric are related at any point by Maxwell's equation

$$\nabla \times E = - \frac{\partial B}{\partial t} . \quad (23)$$

Integrating eq. (23) over the cross-sectional area of the cylinder,  $A = \pi r^2$ , and applying Stokes's theorem gives, if  $H$  is varying sinusoidally in time,

$$\oint_0^{2\pi r} E \cdot d\ell = -j\omega\mu HA, \quad (24)$$

or

$$2\pi r E_\phi = -j\omega\mu H_z A. \quad (25)$$

It should be noted that eq. (24) and eq. (25) are independent of the dielectric properties of the cylinder.

Solving eq. (25) for the magnitude of the electric field gives

$$|E| = \frac{\omega\mu H r}{2} . \quad (26)$$

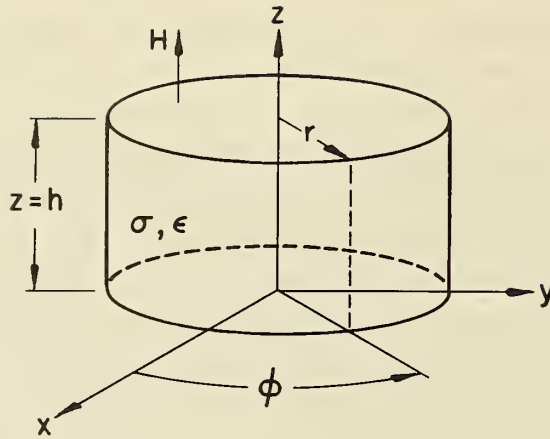


Figure 6. The geometry of the lossy dielectric specimen is in the form of a right circular cylinder of radius,  $r$ , and height,  $h$ , with  $H$  parallel to the axis of the cylinder.

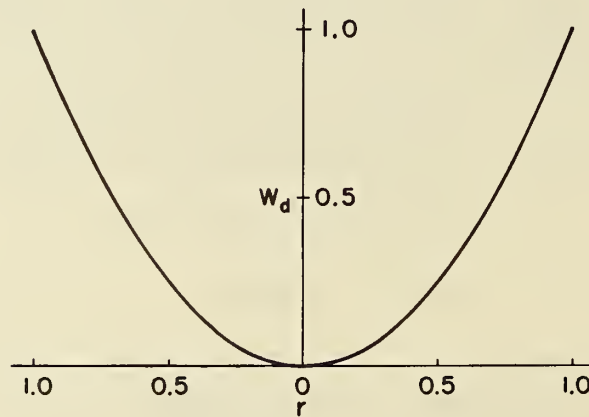


Figure 7. Normalized radial variation of the density of the rf heating,  $W_d$ , in a cylindrical lossy dielectric specimen from eq. (28), with  $H$  parallel to the axis of the cylinder.

The electric-field lines are in the form of closed circles concentric with the axis of the cylinder and lie in planes normal to H. As can be seen, the magnitude of E is proportional to the radius of the circular E-line in question.

6.3 RF-Heating Density. The density of the rf heating in the dielectric cylinder is given by

$$W_d = E^2 \sigma, \quad W/m^3. \quad (27)$$

The component of the rf heating due to the magnetic field within the cylinder can be determined by substituting eq. (26) in eq. (27), giving

$$W_d = \frac{\sigma \omega^2 \mu^2 H^2 r^2}{4}, \quad W/m^3. \quad (28)$$

Note that the resulting density of the rf heating, and hence the hazard, is proportional to  $\omega^2$ , and  $H^2$ , and is independent of the dielectric constant,  $\epsilon$ , of the specimen. As can be seen, this heating is not uniformly distributed within the cylindrical specimen, but is zero along the axis and increases as the square of the radius,  $r$ , as indicated by eq. (28) and shown in figure 7. There is no variation in  $W_d$  in either the  $\phi$  or  $z$  directions.

Equation (28) neglects internal thermal-conduction processes within the specimen, including blood flow in the case of living tissue, which tend to equalize differences in the heat distribution. Thermal radiation must also be taken into account as indicated later in eq. (33) of section 7. In a practical sense, the total rf heating within the specimen may be more meaningful than the rf heating density, since the former averages the rf heating density and may correlate better with the actual temperature rise within the specimen.

6.4 Total RF Heating. The total rf heating in the dielectric cylinder will now be determined. The heating,  $dW$ , in a differential volume,  $dV$ , is given by

$$dW = E^2 \sigma dV = E^2 \sigma r dr d\phi dz. \quad (29)$$



Substituting eq. (26) in eq. (29) gives

$$dW = \frac{\sigma \omega^2 \mu^2 H^2 r^3}{4} dr d\phi dz. \quad (30)$$

Integrating eq. (30) over the volume of the dielectric cylinder, gives for the total rf heating,

$$W_t = \frac{\sigma \omega^2 \mu^2 H^2 \pi r^4 h}{8}, \text{ watts.} \quad (31)$$

Note that the total heating is proportional to  $r^4$ , whereas the heating density given by eq. (28) is proportional to  $r^2$ .

#### Example

Assume the following values for the various parameters in eq. (31):

$$\begin{aligned} \sigma &= 1 \text{ mho/m (1\% saline solution);} \\ \epsilon_r &= 80; \\ f &= 26 \times 10^6, \text{ Hz;} \\ \mu &= 4\pi \times 10^{-7}, \text{ H/m;} \\ H &= 50 \text{ A/m (rms);} \\ r &= 0.1 \text{ m; and} \\ h &= 0.1 \text{ m.} \end{aligned}$$

The total heating,  $W_t$ , generated in the dielectric cylinder due to the value of magnetic field assumed above is 414 watts. It is interesting to note that over 90 percent of this total heat is developed in the outer half of the cylinder radius, as determined from eq. (31), or that half of the heat is developed in the outer 16 percent of the radius.

## 7. THE CALORIMETER METHOD

The value of the rf power absorbed by a specimen can be verified using a calorimeter technique if a one-percent saline (Ringer's) solution [10] is used as the specimen. This test solution (made up with distilled water) should be used at room temperature in a low-dielectric-constant polyfoam-insulated container. The walls should have a thickness of 2.5 cm or more to minimize heat loss due to thermal convection. The rf power being absorbed by the solution (in the absence of thermal radiation and convection) is given by

$$W = 69.77 w(T_2 - T_1)/(t_2 - t_1), \quad (32)$$

where  $W$  = net power absorbed by the water, watts.

$w$  = weight of the water, kg.

$(T_2 - T_1)$  = temperature increase of the water, °C.

$(t_2 - t_1)$  = elapsed time of the experiment, minutes.

If the weight of the water is 3.14 kg, from the example in section 6;  $W = 414$  watts; and  $t_2 - t_1 = 1$  min; the resulting temperature rise will be  $1.89^\circ\text{C}$ . Due to the small temperature rise involved, the effect of thermal radiation can be shown to be negligible from the Stefan-Boltzmann equation for black-body radiation [11]

$$W_{bb} = \sigma T^4, \quad \text{W/m}^2, \quad (33)$$

where  $\sigma = 5.672 \times 10^{-8}$ , and  $T$  = absolute temperature, K.

This will provide a worst-case test since a black-body radiates more than any other at or near room temperature. The body is in thermal equilibrium when its surroundings are at room temperature,  $T_1$ . The net radiation will therefore be proportional to the difference of the fourth powers of the two absolute temperatures involved,  $T_1$  and  $T_2$ , where  $T_2 = T_1 + 1.89^\circ$ , from the above example. This gives a value of  $W_{bb} = 1.11$  watts (for a surface area of  $0.126 \text{ m}^2$ ) which is negligible compared to the 414 watts used in the above example.

## 8. CONCLUSIONS

A relatively simple method is described in this report for accurately determining the rf power being absorbed by a biological specimen during radiation exposure testing in the frequency range 10 to 100 MHz. The method will permit researchers to monitor rf power flow into the specimen measuring solely forward and reflected power on the feed line to the field generator with conventional rf wattmeters. The measurement accuracy is limited only by the quality of the directional couplers used. A commercially available automatic data-acquisition system can be used to "read" the meters and rapidly calculate, display, and record the rf power flow.

The method has the advantage that the exact measuring point on the feed line is not critical, as it is with methods employing direct impedance measurements. The required forward-and reflected-power measurements can be made at a full operating rf power level of several hundred watts without interfering with the exposure tests. Additional measurements to determine the effect of changes in the field distribution within the field generator, caused by the presence of the biological specimen, would be helpful.

## 9. REFERENCES

- [1] Greene, F.M., Development and Construction of an Electromagnetic Near-Field Synthesizer, NBS Technical Note 652 (May 1974). For sale by the Superintendent of Documents, U.S. Government Printing Office, Washington, D.C. 20402 (Order by SD Catalog No. C13.46:652) \$0.50.
- [2] Guy, A.W., Electromagnetic Power Deposition in Man Exposed to High-Frequency Fields and the Associated Thermal and Physiologic Consequences (University of Washington, Seattle, Washington 98105). Report SAM-TR-73-13 to USAF School of Aerospace Medicine, Brooks Air Force Base, Texas 78235 (December 1973).
- [3] Terman, F.E., Electronic and Radio Engineering, pp. 466-473 (McGraw-Hill Book Co., Inc., New York, N.Y., 1955).
- [4] Guillemin, E.A., Communication Networks, Vol. II, pp. 33-79 (John Wiley & Sons, Inc., New York, N.Y., 1935).
- [5] Ramo, S., Whinnery, J.R., and VanDuzer, T., Fields and Waves in Communication Electronics, pp. 23-53 (John Wiley & Sons, Inc., New York, N.Y., 1965).
- [6] Ibid., pp. 565-568.
- [7] Montgomery, C.G., Dicke, R.H., and Purcell, E.M., Principles of Microwave Circuits, M.I.T. Radiation Laboratory Series, Vol. 8, p. 367 (McGraw-Hill Book Co., Inc., New York, N.Y., 1948).
- [8] Terman, F.E., op. cit., pp. 44-57.
- [9] IEEE Standard Dictionary of Electrical and Electronics Terms, p. 520. The Institute of Electrical and Electronics Engineers, Inc. (John Wiley & Sons, Inc., New York, N.Y., 1972).
- [10] Kall, A.R., Radiation Hazards Caused by High-Power, High-Frequency Fields, p. III-27 (Ark Electronics Corporation, Willow Grove, Pennsylvania) report to United States Information Agency, Washington, D.C. (June 1968).
- [11] Machol, R.E., System Engineering Handbook, p. 16-3 (McGraw-Hill Book Co., Inc., New York, N.Y., 1965).

U.S. DEPT. OF COMM. BIBLIOGRAPHIC DATA SHEET	1. PUBLICATION OR REPORT NO. NBS TN-687	2. Gov't Accession No.	3. Recipient's Accession No.
4. TITLE AND SUBTITLE Measurement of RF Power-Absorption in Biological Specimens (10 to 100 MHz)		5. Publication Date November 1976	6. Performing Organization Code 276.03
7. AUTHOR(S) Frank M. Greene		8. Performing Organ. Report No. NBS TN-687	
9. PERFORMING ORGANIZATION NAME AND ADDRESS NATIONAL BUREAU OF STANDARDS DEPARTMENT OF COMMERCE WASHINGTON, D.C. 20234		10. Project/Task/Work Unit No. 2763281	11. Contract/Grant No.
12. Sponsoring Organization Name and Complete Address (Street, City, State, ZIP) Same as item 9		13. Type of Report & Period Covered Final 1/1/76 to 10/1/76	14. Sponsoring Agency Code
15. SUPPLEMENTARY NOTES This report is the outgrowth of work originally sponsored by the USAF School of Aerospace Medicine at Brooks Air Force Base, and the National Institute for Occupational Safety and Health of the Department of HEW.			
16. ABSTRACT (A 200-word or less factual summary of most significant information. If document includes a significant bibliography or literature survey, mention it here.)  A method is described for accurately determining the rf power being absorbed by a biological specimen during non-ionizing radiation-exposure testing using the NBS RF Near-Field Synthesizer in the frequency range 10 to 100 MHz. This method is based solely on measuring the forward and reflected power on the transmission line feeding the synthesizer. Commercially available rf wattmeters can be used, and an automatic data-acquisition system employed, if desired, to "read" the meters and rapidly calculate, display, and record the rf power flow.  The method has the advantage that the exact measuring point on the feed line is not critical, as it is with methods employing direct impedance measurements, and that the required measurements can be made without interfering with the exposure tests.			
17. KEY WORDS (six to twelve entries; alphabetical order; capitalize only the first letter of the first key word unless a proper name; separated by semicolons) Electromagnetic field hazards; electromagnetic field synthesizer; electromagnetic radiation-exposure testing (non-ionizing); near-fields; rf biological hazards.			
18. AVAILABILITY <input checked="" type="checkbox"/> Unlimited  <input type="checkbox"/> For Official Distribution. Do Not Release to NTIS  <input checked="" type="checkbox"/> Order From Sup. of Doc., U.S. Government Printing Office Washington, D.C. 20402, SD Cat. No. C13.46:687  <input type="checkbox"/> Order From National Technical Information Service (NTIS) Springfield, Virginia 22151		19. SECURITY CLASS (THIS REPORT)  UNCLASSIFIED	21. NO. OF PAGES  29
		20. SECURITY CLASS (THIS PAGE)  UNCLASSIFIED	22. Price  \$0.75



# NBS TECHNICAL PUBLICATIONS

## PERIODICALS

**JOURNAL OF RESEARCH** reports National Bureau of Standards research and development in physics, mathematics, and chemistry. It is published in two sections, available separately:

• **Physics and Chemistry (Section A)**

Papers of interest primarily to scientists working in these fields. This section covers a broad range of physical and chemical research, with major emphasis on standards of physical measurement, fundamental constants, and properties of matter. Issued six times a year. Annual subscription: Domestic, \$17.00; Foreign, \$21.25.

• **Mathematical Sciences (Section B)**

Studies and compilations designed mainly for the mathematician and theoretical physicist. Topics in mathematical statistics, theory of experiment design, numerical analysis, theoretical physics and chemistry, logical design and programming of computers and computer systems. Short numerical tables. Issued quarterly. Annual subscription: Domestic, \$9.00; Foreign, \$11.25.

**DIMENSIONS/NBS (formerly Technical News Bulletin)**—This monthly magazine is published to inform scientists, engineers, businessmen, industry, teachers, students, and consumers of the latest advances in science and technology, with primary emphasis on the work at NBS. The magazine highlights and reviews such issues as energy research, fire protection, building technology, metric conversion, pollution abatement, health and safety, and consumer product performance. In addition, it reports the results of Bureau programs in measurement standards and techniques, properties of matter and materials, engineering standards and services, instrumentation, and automatic data processing.

Annual subscription: Domestic, \$9.45; Foreign, \$11.85.

## NONPERIODICALS

**Monographs**—Major contributions to the technical literature on various subjects related to the Bureau's scientific and technical activities.

**Handbooks**—Recommended codes of engineering and industrial practice (including safety codes) developed in cooperation with interested industries, professional organizations, and regulatory bodies.

**Special Publications**—Include proceedings of conferences sponsored by NBS, NBS annual reports, and other special publications appropriate to this grouping such as wall charts, pocket cards, and bibliographies.

**Applied Mathematics Series**—Mathematical tables, manuals, and studies of special interest to physicists, engineers, chemists, biologists, mathematicians, computer programmers, and others engaged in scientific and technical work.

**National Standard Reference Data Series**—Provides quantitative data on the physical and chemical properties of materials, compiled from the world's literature and critically evaluated. Developed under a world-wide program coordinated by NBS. Program under authority of National Standard Data Act (Public Law 90-396).

## BIBLIOGRAPHIC SUBSCRIPTION SERVICES

The following current-awareness and literature-survey bibliographies are issued periodically by the Bureau: **Cryogenic Data Center Current Awareness Service.** A literature survey issued biweekly. Annual subscription: Domestic, \$20.00; Foreign, \$25.00.

**Liquified Natural Gas.** A literature survey issued quarterly. Annual subscription: \$20.00.

**NOTE:** At present the principal publication outlet for these data is the Journal of Physical and Chemical Reference Data (JPCRD) published quarterly for NBS by the American Chemical Society (ACS) and the American Institute of Physics (AIP). Subscriptions, reprints, and supplements available from ACS, 1155 Sixteenth St. N.W., Wash. D. C. 20056.

**Building Science Series**—Disseminates technical information developed at the Bureau on building materials, components, systems, and whole structures. The series presents research results, test methods, and performance criteria related to the structural and environmental functions and the durability and safety characteristics of building elements and systems.

**Technical Notes**—Studies or reports which are complete in themselves but restrictive in their treatment of a subject. Analogous to monographs but not so comprehensive in scope or definitive in treatment of the subject area. Often serve as a vehicle for final reports of work performed at NBS under the sponsorship of other government agencies.

**Voluntary Product Standards**—Developed under procedures published by the Department of Commerce in Part 10, Title 15, of the Code of Federal Regulations. The purpose of the standards is to establish nationally recognized requirements for products, and to provide all concerned interests with a basis for common understanding of the characteristics of the products. NBS administers this program as a supplement to the activities of the private sector standardizing organizations.

**Consumer Information Series**—Practical information, based on NBS research and experience, covering areas of interest to the consumer. Easily understandable language and illustrations provide useful background knowledge for shopping in today's technological marketplace.

*Order above NBS publications from: Superintendent of Documents, Government Printing Office, Washington, D.C. 20402.*

*Order following NBS publications—NBSIR's and FIPS from the National Technical Information Services, Springfield, Va. 22161.*

**Federal Information Processing Standards Publications (FIPS PUBS)**—Publications in this series collectively constitute the Federal Information Processing Standards Register. Register serves as the official source of information in the Federal Government regarding standards issued by NBS pursuant to the Federal Property and Administrative Services Act of 1949 as amended, Public Law 89-306 (79 Stat. 1127), and as implemented by Executive Order 11717 (38 FR 12315, dated May 11, 1973) and Part 6 of Title 15 CFR (Code of Federal Regulations).

**NBS Interagency Reports (NBSIR)**—A special series of interim or final reports on work performed by NBS for outside sponsors (both government and non-government). In general, initial distribution is handled by the sponsor; public distribution is by the National Technical Information Services (Springfield, Va. 22161) in paper copy or microfiche form.

**Superconducting Devices and Materials.** A literature survey issued quarterly. Annual subscription: \$20.00. Send subscription orders and remittances to the preceding bibliographic services to National Bureau of Standards, Cryogenic Data Center (275.02) Boulder, Colorado 80302.

**U.S. DEPARTMENT OF COMMERCE**  
**National Bureau of Standards**  
Washington, D.C. 20234

OFFICIAL BUSINESS

Penalty for Private Use, \$300

POSTAGE AND FEES PAID  
U.S. DEPARTMENT OF COMMERCE  
COM-215



Fourth Class Mail



**75 YEARS**  
**NBS**  
1901-1976